

2025 PRODUCT CATALOGUE.

Contents

Product



Examination / Treatment
Couches & Chairs
Pg 2-17



Healthcare Certified
Chairs & Stools
Pg 18-33



Medical Trolleys &
Storage Solutions
Pg 34-59



Clinical Bins & Sack
Holders
Pg 60-62



Mobile Patient Privacy
Screens
Pg 63-64



Medicine Cabinets
Pg 65-72



Ward Furniture
Pg 73-79



Stainless Steel
Healthcare Furniture
Pg 80-85



HTM Compliant Cabinets
& Worksurfaces
Pg 86-92

Other



Upholstery Colours
Pg 93-94



About Us
Pg 95-96



Material Cleaning
Guidance
Pg 97

clean
simple
examination

Examination/Treatment Couches & Chairs

Our examination and treatment couches and chairs are durable, easy to maintain, and designed for patient comfort. Suitable for varied healthcare settings, they provide reliable, multi-purpose solutions that enhance patient care and support clinical efficiency.

Variable Height Couches



Fusion Eco Variable Height Couches
Pg 3-4



Fusion Drop End Multi-Discipline Couches
Pg 5

Treatment Chairs



Fusion Treatment Chair
Pg 6



Fusion Gynae Variable Height
Pg 7-8



Fusion Podiatry Chairs
Pg 9



Fusion Phlebotomy Chairs
Pg 10



Phlebotomy Treatment Chairs
Pg 11

Fixed Height Couches



Treatment Couch Systems
Pg 12-15



Examination Couches
Pg 16-17

Accessories



Couch Steps
Pg 17

Fusion Eco 3 Section Variable Height Couches

product information



Electric or Hydraulic height adjustment 45-96cm



Maximum user weight (40st.) 250kg



Gas strut assisted multi-position head and foot sections



Back rest and leg rest adjustment +70°/-35°



Easy to use handset (Electric models only)



Retractable positioning castors



Adjustable height levelling foot



Fully upholstered head & foot sections



Double layer density foam for support & comfort



Luxury 65mm deep contoured upholstery



Anti-bacterial/anti-microbial vinyl



Paper roll holder (supplied as standard)



Easy clean design and materials



Fire (Flammability) Testing Flammability BS7176 Medium Hazard BS5852 : 2006 Clause 11 (composites)

accessories * these accessories can only be fitted when ordered with the couch



Sun-FECO/LEG1 Leg Supports



Sun-FECO/PHB1 Phlebotomy Arm (each)



Sun-FECO/FSW1* Foot Switch



Sun-FECO/SUP1 Gynae Stirrups



Sun-FECO/CBH1* Breathing Hole



Sun-FECO/CSR* Folding Side Rails (pair)



Sun-CHP1 Head Pillow



Sun-FECO/HEO* Electrically Operated Head Section Upgrade



Sun-FECO/FEO* Electrically Operated Foot Section Upgrade



Sun-INST1* Installation Service & Packaging Recycling




Sun-DISPOSAL1** Removal & Disposal of Old Couch


installation

** Removal from ground floor only - weight & removal restrictions apply, please enquire


dimensions

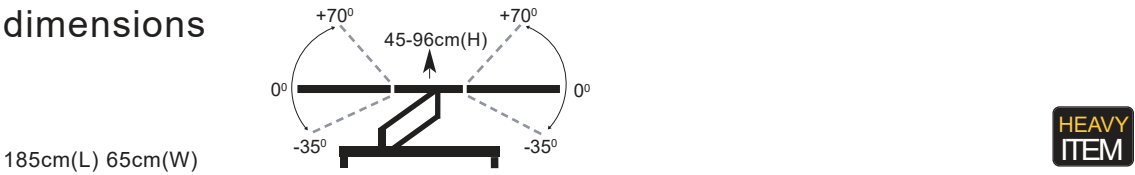


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

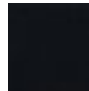


BLACK




BLUEBELL




BLUEBERRY




FEATHER



LAPIS



SNOW-WHITE



SUNFLOWER

+23 MORE AVAILABLE


please see inside back cover for full range and specification of our vinyl range

Fusion Eco 3 Section Variable Height Couches

The new Fusion Eco couch range provides most of the functionality you would expect for our standard Fusion range at a more affordable price to suit individual budgets. Fusion Eco is simply, strong, stable, aesthetically pleasing and easy to clean.



Fully upholstered head & foot sections with gas assisted adjustment

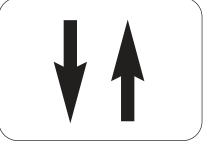

Electric or Hydraulic height adjustment



Increased maximum user weight to 250kg



Optional Electric Head & Foot Sections (electric models only)



Warranty extended to 15 years on the frame. For full warranty details, please see the T&Cs.



Code: Sun-FECO3E
Fusion Eco 3 Section Electric Couch (Shown with optional accessories)



Code: Sun-FECO3H
Fusion Eco 3 Section Hydraulic Couch



Hydraulic height adjustment

Fusion Eco 2 Section Variable Height Couches

product information



Electric or Hydraulic height adjustment 45-96cm



250kg
Maximum user weight (40st.)



Gas strut assisted multi-position head and foot sections



Back rest and leg rest adjustment +70°/-35°



Easy to use handset (Electric models only)



Retractable positioning castors



Adjustable height levelling foot



Fully upholstered head & foot sections



Double layer density foam for support & comfort



Luxury 65mm deep contoured upholstery



Anti-bacterial/anti-microbial vinyl



Paper roll holder (supplied as standard)




Easy clean design and materials



Anti-bacterial Vinyl
Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)

accessories * these accessories can only be fitted when ordered with the couch




Sun-FECO/LEG1 Leg Supports




Sun-FECO/PHB1 Phlebotomy Arm (each)




Sun-FECO/FSW1* Foot Switch

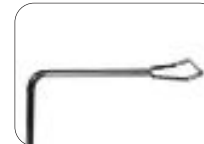


Sun-FECO/SUP1 Gynae Stirrups




Sun-CHP1 Head Pillow




Sun-FECO/CBH1* Breathing Hole



Sun-FECO/CSR* Folding Side Rails (pair)



Sun-FECO/HEO* Electrically Operated Head Section Upgrade




Sun-INST1* Installation Service & Packaging Recycling



Sun-DISPOSAL1** Removal & Disposal of Old Couch

installation





** Removal from ground floor only - weight & removal Restrictions apply, please enquire


dimensions




185cm(L) 65cm(W)




anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

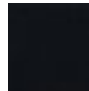


BLACK




BLUEBELL




BLUEBERRY




FEATHER




LAPIS



SNOW-WHITE



SUNFLOWER



+23 MORE AVAILABLE


please see inside back cover for full range and specification of our vinyl range

Fusion Eco 2 Section Variable Height Couches

The new Fusion Eco couch range provides most of the functionality you would expect for our standard Fusion range at a more affordable price to suit individual budgets. Fusion Eco is simply, strong, stable, aesthetically pleasing and easy to clean.



Fully upholstered head section with gas assisted adjustment

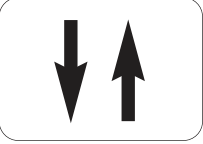

Electric or Hydraulic height adjustment



Increased maximum user weight to 250kg



Optional Electric Head Section

(electric models only)



Warranty extended to 15 years on the frame. For full warranty details, please see the T&Cs.



Code: Sun-FECO2E
Fusion Eco 2 Section Electric Couch (Shown with optional accessories)




Code: Sun-FECO2H
Fusion Eco 2 Section Hydraulic Couch




Hydraulic height adjustment

Fusion Drop End Multi-Discipline Couch


product information

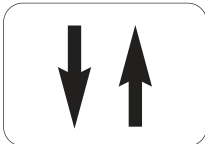


Electric or Hydraulic height adjustment




Maximum user weight (40st.)




Optional breathing hole & plug in foot section for physio use




50mm Retractable positioning castors




Adjustable height levelling foot




Gas assisted head & foot sections as standard




Electric head/foot section adjustment (Sun-FDECE2)




Tilting Seat for Trendelenburg (Sun-FDECE2)




Hoist clearance 13cm




Range of additional accessories for flexibility



Fully upholstered head & foot sections




Luxury 65mm deep contoured upholstery




Anti-bacterial/anti-microbial vinyl




Paper roll holder (supplied as standard)



Easy clean design and materials



Anti-bacterial Vinyl

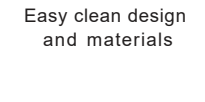

Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)

accessories * these accessories can only be fitted when ordered with the couch



Sun-CHP1 Head Pillow



Sun-FDEC/CBH2* Breathing Hole - in foot section



Sun-FDEC/LEG1 Leg Supports(pair)



Sun-FDEC/PHB1 Phlebotomy Arm (pair)



Sun-FDEC/PAS1 Padded Arm Supports(pair)



Sun-FDEC/SUP1 Stirrups (pair)



Sun-FDEC/CSR* Folding Side Rails (pair)




Sun-FDEC/FSW1* Foot Switch




Sun-DISPOSAL1** Removal & Disposal of Old Couch

** Removal from ground floor only - weight and removal restrictions apply


dimensions

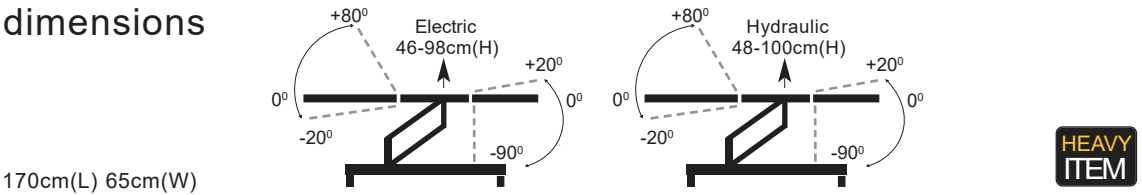


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

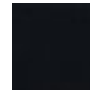


BLACK




BLUEBELL




BLUEBERRY




FEATHER




LAPIS



SNOW-WHITE



SUNFLOWER



+23 MORE AVAILABLE

please see inside back cover for full range and specification of our vinyl range

Fusion Drop End Multi-Discipline Couch

The new Fusion Drop End Couch range has been designed using our extensive knowledge to provide a trusted, ergonomic design providing comfort for the patient, ease of access via the drop end design and comfortable patient positioning for the clinician. The couch offers amazing flexibility and can be configured into a wide range of specialist treatment uses including physio, gynae, phlebotomy or treatment chair via the addition of relevant accessories.



Code: Sun-FDECE2
Fusion Drop End Multi-Discipline Couch
Electric Height Adjustment
Electric back and foot sections



Warranty extended to 15 years on the frame. For full warranty details, please see the T&Cs.



Code: Sun-FDECE1
Fusion Drop End Multi-Discipline Couch
Electric Height Adjustment
Gas assisted head and foot sections



Hydraulic height adjustment

Code: Sun-FDECH1
Fusion Drop End Multi-Discipline Couch
Hydraulic Height Adjustment
Gas assisted head and foot sections

Fusion Treatment Chairs

product information



Powered tilting seat
Sun-FTRTE2/3/5/6 only)



-20°
TRENDELENBURG
Tilting Seat for Trendelenburg
Sun-FTRTE2/3/5/6 only)



Head adjustment
+80° / -20°



Adjustable foot section
+20° / -90°



2 Adjustable/removable
padded arm supports



250kg
Maximum user
weight (40st.)



Easy to use handset



50mm Retractable
positioning castors



Adjustable height
levelling foot



Dual sided foot pedals
(Hydraulic models only)



Hoist clearance
13cm



Fully upholstered head
& foot sections



Luxury 65mm deep
contoured upholstery



Anti-bacterial/
anti-microbial vinyl



Paper roll holder
(supplied as standard)



Easy clean design
and materials



West Yorkshire
Materials
Testing Service
Anti-bacterial Vinyl

Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11
(composites)

accessories * these accessories can only be fitted when ordered with the couch



Sun-FTRT/PHB1
Phlebotomy Arm (pair)



Sun-CHP1
Head Pillow



Sun-FTRT/CBH1*
Breathing Hole - in foot section
(Sun-FTRTE1,2,3 & Sun-FTRTH1 only)



Sun-FTRT/FSW1*
Foot Switch
(Sun-FTRTE1 & 4)



Sun-FTRT/FSW2*
Foot Switch 2 Section
(Sun-FTRTE2 & 5)



Sun-FTRT/FSW3*
Foot Switch 3 Section
(Sun-FTRTE3 & 6)




installation
Sun-INST1*
Installation Service &
Packaging Recycling




Sun-DISPOSAL1*/**
Removal & Disposal
of Old Couch

** Removal from
ground floor only -
weight & removal
restrictions apply,
please enquire


dimensions

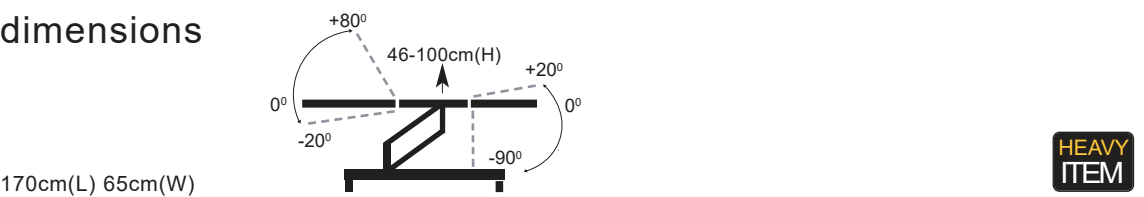


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

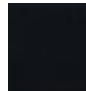


BLACK




BLUEBELL




BLUEBERRY




FEATHER




LAPIS



SNOW-WHITE



SUNFLOWER



**+23 MORE
AVAILABLE**

please see inside back cover for full range and specification of our vinyl range

Fusion Treatment Chairs

The new Fusion Treatment Chair range has been designed using our extensive knowledge to provide a trusted, ergonomic design providing comfort for the patient, ease of access via the split leg design and comfortable patient positioning for the clinician.



Powered tilting seat



Adjustable foot section



Recovery position



Code: Sun-FTRTE3
Fusion Treatment Chair - Electric Height Single Foot Section
Powered head section and powered tilting seat

Electric Height Single Foot Section

Code: Sun-FTRTE1 - Gas assisted head section and fixed seat
Code: Sun-FTRTE2 - Gas assisted head section and powered tilting seat
Code: Sun-FTRTE3 - Powered head section and powered tilting seat

Electric Height Split Foot Section

Code: Sun-FTRTE4 - Gas assisted head section and fixed seat
Code: Sun-FTRTE5 - Gas assisted head section and powered tilting seat
Code: Sun-FTRTE6 - Powered head section and powered tilting seat

Hydraulic Height Adjustment with Fixed Seat

Code: Sun-FTRTH1 - Gas assisted head section and single foot section
Code: Sun-FTRTH4 - Gas assisted head section and split foot section



Warranty extended
to 15 years on the
frame. For full
warranty details,
please see the
T&Cs.



Code: Sun-FTRTH4
Fusion Treatment Chair - Hydraulic Height Adjustment Fixed Seat
Gas assisted head section and split foot section

product information



Electric or Hydraulic height adjustment
45-98cm



225kg
Maximum user weight (35st.)



Gas assisted back
(-20 - +80°)



Easy to use handset
(Electric models only)



Rear retractable positioning castors



100mm Front braked castors



Cut away lithotomy section for easier access



100mm Deep stainless steel operators tray



Fully adjustable leg supports



Fold down arms for patient comfort & security



Double layer density foam for support & comfort



Anti-bacterial/anti-microbial vinyl



Easy clean design and materials




Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)


accessories * these accessories can only be fitted when ordered with the couch




Sun-FGYN/FSW1*
Foot Switch
(For Gynae 1)




Sun-FGYN/CHP1
Head Pillow




Sun-FGYN/CBH1*
Breathing Hole




Sun-FGYN/SUP1
Stirrups (pair)



Sun-FGYN/WB
Wall Bracket for Foot Section



Sun-FGYN/FS
Detachable Foot Section Extension




Sun-INST1*
Installation Service & Packaging Recycling



Sun-DISPOSAL1***
Removal & Disposal of Old Couch

** Removal from ground floor only - weight and removal restrictions apply, please enquire


dimensions

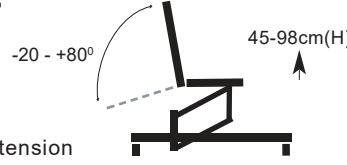


137cm(L) 65cm(W)
187cm with foot extension

45-98cm(H)

-20 - +80°

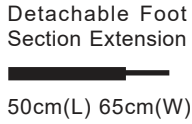


Detachable Foot Section Extension
50cm(L) 65cm(W)




HEAVY ITEM


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

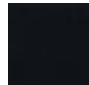


BLACK




BLUEBELL




BLUEBERRY




FEATHER




LAPIS



SNOW-WHITE



SUNFLOWER



+23 MORE AVAILABLE

please see inside back cover for full range and specification of our vinyl range

Fusion Gynae 1 - Electric & Hydraulic Couches

The new Gynae couch range features an ergonomic design providing comfort for the patient and ease of access and patient positioning for the clinician. Fusion Gynae is simply, strong, stable, aesthetically pleasing and easy to clean.



Code: Sun-FGYNE1
Fusion Gynae1 - 2 section couch
Electric powered height, gas assisted back,
no seat tilt with arm supports& leg supports



Warranty extended to 15 years on the frame. For full warranty details, please see the T&Cs.



Code: Sun-FGYNH1
Fusion Gynae1 - 2 section couch
Hydraulic powered height, gas assisted back,
no seat tilt with arm supports& leg supports



Hydraulic height adjustment

product information



Electric height adjustment
45-98cm



225kg
Maximum user weight (35st.)



Electric Powered Tilting Seat 0 to -15°
(Trendelenburg)



Fully powered head section
(Sun-FGYNE3)



Easy to use handset



Rear retractable positioning castors



100mm Front braked castors



Cut away lithotomy section for easier access



100mm Deep stainless steel operators tray



Fully adjustable leg supports



Fold down arms for patient comfort & security



West Yorkshire Materials Testing Service
Anti-bacterial Vinyl



eze
Easy clean design and materials

Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)

accessories * these accessories can only be fitted when ordered with the couch



Sun-FGYN/FSW2*
Foot Switch (Gynae 2)



Sun-FGYN/FSW3*
Foot Switch (Gynae 3)



Sun-FGYN/CHP1
Head Pillow



Sun-FGYN/CBH1*
Breathing Hole



Sun-INST1*
Installation Service & Packaging Recycling



Sun-FGYN/SUP1
Stirrups (pair)



Sun-FGYN/FS
Detachable Foot Section Extension



Sun-FGYN/WB
Wall Bracket for Foot Section



Sun-DISPOSAL1**
Removal & Disposal of Old Couch

** Removal from ground floor only - weight & removal restrictions apply, please enquire

dimensions




45-98cm(H)
-20 - +80°
Seat Tilt 0 to -15°
137cm(L) 65cm(W)
187cm with foot extension


Detachable Foot Section Extension
50cm(L) 65cm(W)




anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

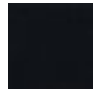

BLACK




BLUEBELL




BLUEBERRY




FEATHER



LAPIS



SNOW-WHITE



SUNFLOWER

+23 MORE AVAILABLE

please see inside back cover for full range and specification of our vinyl range

Fusion Gynae 2 & 3 - Electric Couches

The new Gynae couch range features an ergonomic design providing comfort for the patient and ease of access and patient positioning for the clinician. Fusion Gynae is simply, strong, stable, aesthetically pleasing and easy to clean.



Code: Sun-FGYNE2
Fusion Gynae2 - 2 section couch
Electric powered height, gas assisted back, powered tilting seat with arm supports& leg supports



Warranty extended to 15 years on the frame. For full warranty details, please see the T&Cs.



Code: Sun-FGYNE3
Fusion Gynae3 - 2 section couch
Electric powered height, back and tilting seat with arm supports& leg supports

Fusion Podiatry Chair

product information



Electric or Hydraulic height adjustment



Maximum user weight (40st.)



50mm Retractable positioning castors




Gas assisted head and split foot sections



Electric head section adjustment (Sun-FPODE2)



Adjustable leg length extension up to 15cm



Tilting Seat for Trendelenburg (Sun-FPODE1 & Sun-FPODE2)



Stay level swing back padded arm rests for easy patient access



Adjustable height levelling foot



Hoist clearance 13cm



Fully upholstered head & foot sections



Luxury 65mm deep contoured upholstery



Anti-bacterial/anti-microbial vinyl



Paper roll holder (supplied as standard)



Easy clean design and materials



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)



Anti-bacterial Vinyl

accessories * these accessories can only be fitted when ordered with the couch



Sun-CHP1 Head Pillow



Sun-FPOD/FSW1* Foot Switch (FPODE1)



Sun-FPOD/FSW2* Foot Switch (FPODE2)




Sun-FPOD/FSW2* Foot Switch (FPODE2)




Sun-DISPOSAL1/** Removal & Disposal of Old Couch

** Removal from ground floor only - weight & removal restrictions apply, please enquire

dimensions

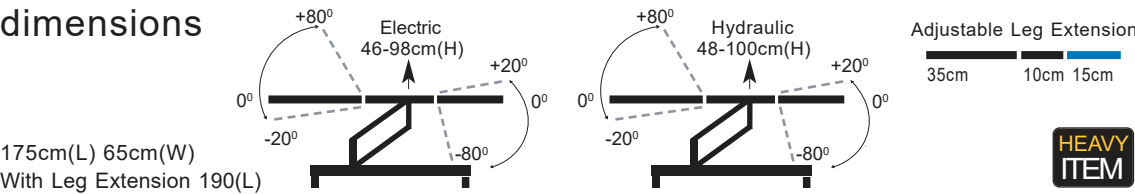

anti-bacterial vinyl upholstery colours - please specify when ordering



APPLE BLACK BLUEBELL BLUEBERRY FEATHER LAPIS SNOW-WHITE SUNFLOWER

+23 MORE AVAILABLE

please see inside back cover for full range and specification of our vinyl range

Fusion Podiatry Chair

The new Fusion Podiatry Chair range has been designed using our extensive knowledge to provide a trusted, ergonomic design providing comfort for the patient, ease of access via the split leg design, Trendelenburg movement as standard on electric models and comfortable patient positioning for the clinician.



Code: Sun-FPODE1
Fusion Podiatry Chair
Electric Height Adjustment
Gas assisted head and split adjustable length foot sections
(Powered Trendelenburg -20° Movement)

Code: Sun-FPODE2 (shown above)
Fusion Podiatry Chair
Electric Height Adjustment
Electric head and split adjustable length foot sections
(Powered Trendelenburg -20° Movement)



Warranty extended to 15 years on the frame. For full warranty details, please see the T&Cs.



Code: Sun-FPODH1
Fusion Podiatry Chair
Hydraulic Height Adjustment
Gas assisted head and split adjustable length foot sections
(No Trendelenburg Movement)



Hydraulic height adjustment

Fusion Phlebotomy Chairs


product information



Electric or Hydraulic height adjustment




250kg
Maximum user weight (40st.)




Gas assisted head & foot sections




Electric head section & foot sections adjustment option (Sun-FPHBE2 & 3)



-20°
TRENDLENBURG
Tilting Seat for Trendelenburg (Sun-FPHBE3)




Two fully adjustable phlebotomy arms as standard



50mm Retractable positioning castors



Adjustable height levelling foot



Dual sided foot pedals (Hydraulic models only)



Hoist clearance 13cm




Fully upholstered head & foot sections



Luxury 65mm deep contoured upholstery



Anti-bacterial/anti-microbial vinyl




Paper roll holder (supplied as standard)



Easy clean design and materials



Anti-bacterial Vinyl

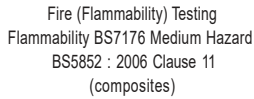

Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)

accessories * these accessories can only be fitted when ordered with the couch



Sun-CHP1 Head Pillow



Sun-FPHB/FSW1* Foot Switch (Sun-FPHBE1)



Sun-FPHB/FSW3* Foot Switch 3 Section (Sun-FPHBE2)



Sun-FPHB/FSW4* Foot Switch 4 Section (Sun-FPHBE3)




Sun-INST1* Installation Service & Packaging Recycling




Sun-DISPOSAL1** Removal & Disposal of Old Couch

** Removal from ground floor only - weight & removal restrictions apply, please enquire


dimensions

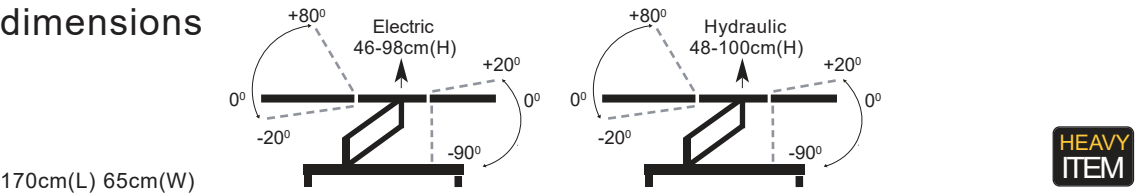


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

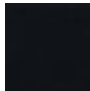


BLACK



BLUEBELL




BLUEBERRY




FEATHER



LAPIS



SNOW-WHITE



SUNFLOWER

+23 MORE AVAILABLE

please see inside back cover for full range and specification of our vinyl range

Fusion Phlebotomy Chairs

The new Fusion Phlebotomy Chair range has been designed using our extensive knowledge to provide a trusted, ergonomic design providing comfort for the patient and comfortable patient positioning for the clinician.



Easy removal of the phlebotomy arms allows the chair to be used as a standard drop end couch



In emergency situations the chair goes head down and feet up for patient recovery



Code: Sun-FPHBE3
Fusion Phlebotomy Chair - Electric Height Adjustment
Electric back and foot sections and tilting seat

Electric Height Adjustment

Code: Sun-FPHBE1 - Electric Height Adjustment gas assisted head and foot sections
Code: Sun-FPHBE2 - Electric Height Adjustment electric back and foot sections
Code: Sun-FPHBE3 - Electric back and foot sections and tilting seat



Warranty extended to 15 years on the frame. For full warranty details, please see the T&Cs.



Code: Sun-FPHBH1
Fusion Phlebotomy Chair - Hydraulic Height Adjustment
Gas assisted head and foot sections



Phlebotomy Chairs

product information



Maximum user weight (17st.)



Lightweight design only 8kgs



Multi-position arm rest, made from a hygienic PU moulding



Arm rest can be positioned either side (2 brackets supplied)



Ergonomically designed moulded seat & back (Sun-PCHA only)



Anti-bacterial/anti-microbial vinyl (Sun-PCHA/VYL only)



40mm Deep foam for support and comfort (Sun-PCHA/VYL only)



Easy clean design and materials



Supplied fully assembled



CE marked



Independently tested



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)




British made

accessories




Sun-PCHA/Arm
Additional Arm Rest

dimensions

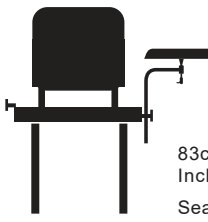

83cm(H) 70cm(W) 56cm(D)
Including arm brackets(excluding arm)
Seat Height 43.5cm(H)

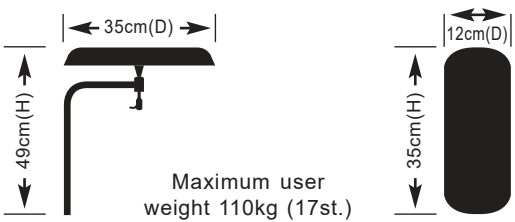

Arm Rest
35cm(D)
49cm(H)
12cm(D)
35cm(H)
Maximum user weight 110kg (17st.)

moulded seat & back colours - please specify when ordering



Black




Blue




Grey


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

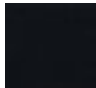


BLACK




BLUEBELL




BLUEBERRY




FEATHER




LAPIS



SNOW-WHITE



SUNFLOWER



please see inside back cover for full range and specification of our vinyl range

Phlebotomy Chairs

Durable, stable and ergonomically designed, phlebotomy chairs are easy to use and adjust and can be used left, right or dual sided with an additional arm, whilst also being easy to clean.




Code: Sun-PCHA
Phlebotomy chair
Moulded seat and back, one arm rest



ODE: Sun-PCHA/VYL
hlebotomy chair
inyl upholstery, one arm rest

Treatment Couch Systems

product information

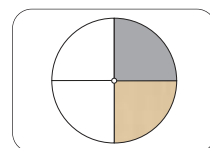

3 Cabinet finishes



Maximum user weight (26st.)



Multi-position self locking backrest



Luxury 65mm deep contoured upholstery



Anti-bacterial/anti-microbial vinyl



Paper roll holder (supplied as standard)



2 Pull-out work surfaces



Cupboards include pull-out work surface & basket, shelf and lock



Locking doors supplied with 2 keys



Durable silver epoxy coated frame



Cabinets supplied fully assembled



Easy clean design and materials



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)



FIRA Standard 6250 2005
Specification:
Adhesion Performance
General/Severe

accessories * these accessories can only be fitted when ordered with the couch



Sun-CBH1*
Breathing Hole



Sun-CHP1
Head Pillow



Sun-LEG1*
Leg Supports



installation

Sun-INST1*
Installation Service & Packaging Recycling



Sun-DISPOSAL1**
Removal & Disposal of Old Couch



Sun-STP1
Single Tier Step (oval)




Sun-STP2
Double Tier Step (oval)



Sun-STP3
Single Tier Step (rectangular)

** Removal from ground floor only - weight and removal restrictions apply, please enquire

dimensions

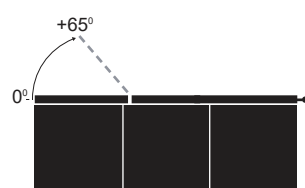

83cm(H) 193cm(L) 60cm(W)
with paper roll holder

HEAVY
ITEM

anti-bacterial vinyl upholstery colours - please specify when ordering



please see inside back cover for full range and specification of our vinyl range

Specialist 2 Section Treatment Couch Systems

Combining functionality with storage and comfort, Specialist Couches assist in making examination and treatment procedures easier and more efficient.



Code: Sun-CS1
Specialist treatment couch system with 2 drawerline units, 1 drawer pack, 2 section upholstery with optional head pillow and leg supports



White (high gloss) cabinet finish



Grey cabinet finish



Beech effect cabinet finish

product information



3 Cabinet finishes



Maximum user weight (26st.)



Multi-position self locking backrest



Luxury 65mm deep contoured upholstery



Anti-bacterial/anti-microbial vinyl



Paper roll holder (supplied as standard)



2 Pull-out work surfaces



Cupboards include pull-out work surface & basket, shelf and lock



Locking doors supplied with 2 keys



Durable silver epoxy coated frame



Cabinets supplied fully assembled



Easy clean design and materials



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)



FIRA Standard 6250 2005
Specification:
Adhesion Performance
General/Severe

accessories * these accessories can only be fitted when ordered with the couch



Sun-CBH1* Breathing Hole



Sun-CHP1 Head Pillow



Sun-LEG1* Leg Supports



Sun-INST1* Installation Service & Packaging Recycling



Sun-DISPOSAL1** Removal & Disposal of Old Couch



Sun-STP1 Single Tier Step (oval)




Sun-STP2 Double Tier Step (oval)



Sun-STP3 Single Tier Step (rectangular)

** Removal from ground floor only - weight and removal restrictions apply, please enquire


dimensions

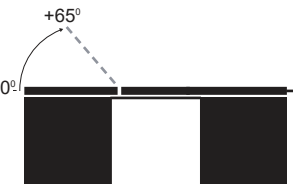


83cm(H) 193cm(L) 60cm(W)
with paper roll holder

HEAVY ITEM


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

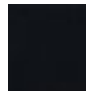


BLACK




BLUEBELL




BLUEBERRY




FEATHER




LAPIS



SNOW-WHITE



SUNFLOWER



please see inside back cover for full range and specification of our vinyl range

Specialist 2 Section Treatment Couch Systems

Combining functionality with storage and comfort, Specialist Couches assist in making examination and treatment procedures easier and more efficient.



Code: Sun-CS3
Specialist treatment couch system with 1 drawerline unit, 1 drawer pack, 2 section upholstery with optional head pillow and leg supports



White (high gloss) cabinet finish




Grey cabinet finish



Beech effect cabinet finish

product information

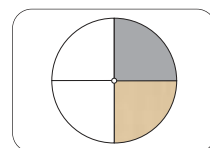

3 Cabinet finishes



Maximum user weight (26st.)



Multi-position self locking backrest



Luxury 65mm deep contoured upholstery



Anti-bacterial/anti-microbial vinyl



Paper roll holder (supplied as standard)



2 Pull-out work surfaces



Cupboards include pull-out work surface & basket, shelf and lock



Locking doors supplied with 2 keys



Durable silver epoxy coated frame



Cabinets supplied fully assembled



Easy clean design and materials



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)



FIRA Standard 6250 2005
Specification:
Adhesion Performance
General/Severe

accessories * these accessories can only be fitted when ordered with the couch



Sun-CBH1* Breathing Hole



Sun-CHP1 Head Pillow



Sun-LEG1* Leg Supports



Sun-INST1* Installation Service & Packaging Recycling



Sun-DISPOSAL1** Removal & Disposal of Old Couch



Sun-STP1 Single Tier Step (oval)




Sun-STP2 Double Tier Step (oval)



Sun-STP3 Single Tier Step (rectangular)

** Removal from ground floor only - weight and removal restrictions apply, please enquire


dimensions

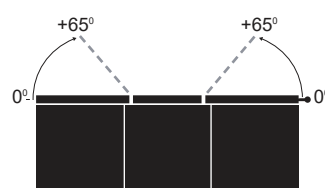

83cm(H) 193cm(L) 60cm(W)
with paper roll holder



anti-bacterial vinyl upholstery colours - please specify when ordering

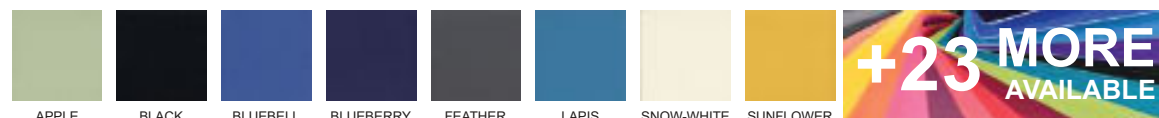

please see inside back cover for full range and specification of our vinyl range

Specialist 3 Section Treatment Couch Systems

Combining functionality with storage and comfort, Specialist Couches assist in making examination and treatment procedures easier and more efficient.



Code: Sun-CS1/3S
Specialist treatment couch system with 2 drawerline units, 1 drawer pack, 3 section upholstery with optional head pillow



White (high gloss) cabinet finish



Grey cabinet finish



Beech effect cabinet finish

product information



3 Cabinet finishes



Maximum user weight (26st.)



Multi-position self locking backrest



Luxury 65mm deep contoured upholstery



Anti-bacterial/anti-microbial vinyl



Paper roll holder (supplied as standard)



2 Pull-out work surfaces



Cupboards include pull-out work surface & basket, shelf and lock



Locking doors supplied with 2 keys



Durable silver epoxy coated frame



Cabinets supplied fully assembled



Easy clean design and materials



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)



FIRA Standard 6250 2005
Specification:
Adhesion Performance
General/Severe

accessories * these accessories can only be fitted when ordered with the couch



Sun-CBH1* Breathing Hole



Sun-CHP1 Head Pillow



Sun-LEG1* Leg Supports



Sun-INST1* Installation Service & Packaging Recycling



Sun-DISPOSAL1** Removal & Disposal of Old Couch



Sun-STP1 Single Tier Step (oval)




Sun-STP2 Double Tier Step (oval)



Sun-STP3 Single Tier Step (rectangular)

** Removal from ground floor only - weight and removal restrictions apply, please enquire


dimensions

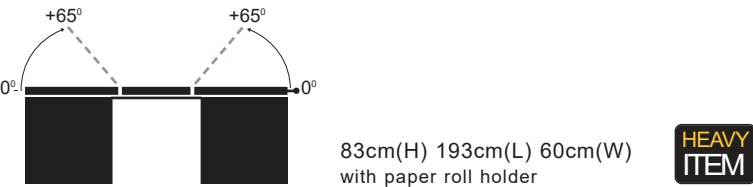


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

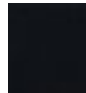


BLACK




BLUEBELL




BLUEBERRY




FEATHER




LAPIS



SNOW-WHITE



SUNFLOWER



please see inside back cover for full range and specification of our vinyl range

Specialist 3 Section Treatment Couch Systems

Combining functionality with storage and comfort, Specialist Couches assist in making examination and treatment procedures easier and more efficient.



Code: Sun-CS3/3S
Specialist treatment couch system with 1 drawerline unit, 1 drawer pack, 3 section upholstery with optional head pillow



White (high gloss) cabinet finish

Grey cabinet finish



Beech effect cabinet finish

product information



3 Cabinet finishes



Maximum user weight (26st.)



Multi-position self locking backrest



Luxury 65mm deep contoured upholstery



Anti-bacterial/anti-microbial vinyl



Paper roll holder (supplied as standard)



Durable silver epoxy coated frame



Cupboards include pull-out work surface & basket, shelf and lock



60cm diameter legs for stability with adjustable feet for uneven floors



Choice of integrated storage solutions



Optional couch steps (see pages 42-43)



Easy clean design and materials



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)



FIRA
Independently tested by FIRA
Specification:
Adhesion Performance
General/Severe

accessories * these accessories can only be fitted when ordered with the couch



Sun-CBH1* Breathing Hole



Sun-CHP1 Head Pillow



Sun-MCD1 Single Drawer Unit



Sun-CPS2 Pull-out Work Surface



Sun-MPT51 Vista 15 Under Couch Trolley



Sun-CMU1 Mobile Storage Unit (with locking door)



Sun-PB1* Pull-out Basket for Storage Unit



Sun-LEG1* Leg Supports (pair)



Sun-INST1* Installation Service & Packaging Recycling




Couch Steps (see nextpage)

** Removal from ground floor only - weight and removal restrictions apply, please enquire


dimensions



79cm(H) 193cm(L) 60cm(W)
with paper roll holder




HEAVY ITEM




Mobile Storage Unit
58cm(H) 50cm(W) 60cm(D)


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

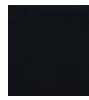


BLACK




BLUEBELL




BLUEBERRY




FEATHER




LAPIS



SNOW-WHITE



SUNFLOWER



+23 MORE AVAILABLE

please see inside back cover for full range and specification of our vinyl range

Practitioner Deluxe Examination Couches

Stylish, practical and comfortable, Practitioner Couches make examination and treatment procedures easier, combined with contemporary aesthetics.



Code: Sun-MEC1
Practitioner deluxe examination couch
with optional pull-out work surface, drawer units,
mobile storage unit, leg supports, step and head pillow



White (high gloss) cabinet finish



Grey cabinet finish



Beech effect cabinet finish



Code: Sun-MPT51
Vista 15 under couch trolley
4 single depth trays

product information



Low height
only 50cm high



Maximum user
weight (26st.)



Multi-position
self locking backrest



Single layer density foam
for support & comfort



Luxury 65mm deep
contoured upholstery



Anti-bacterial/
anti-microbial vinyl



Paper roll holder
(supplied as standard)



Durable silver epoxy
coated frame



60cm diameter
legs for stability



4 adjustable feet
for uneven floors



Easy clean design
and materials




Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998




Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11
(composites)


accessories * these accessories can only be fitted when ordered with the couch




Sun-CBH1*
Breathing Hole




Sun-CHP1
Head Pillow



Sun-INST1*
Installation Service &
Packaging Recycling



Sun-DISPOSAL1**
Removal & Disposal
of Old Couch



Couch Step

** Removal from ground floor only - weight and removal restrictions apply, please enquire

dimensions




0°
+65°




HEAVY
ITEM

50cm(H) 193cm(L) 60cm(W)
with paper roll holder


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

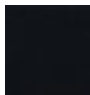


BLACK




BLUEBELL




BLUEBERRY




FEATHER




LAPIS



SNOW-WHITE



SUNFLOWER



+23 MORE
AVAILABLE

please see inside back cover for full range and specification of our vinyl range

Alberti Rest Couches

Designed for resting, recovery or where a reduced height application is required, Alberti Couches are strong, durable, easy to clean while being exceptionally comfortable.



Code: Sun-RC1
Alberti rest couch with optional head pillow

Couch Steps


Designed for using with fixed height couches, couch steps are the perfect solution for assisting patients to raise or lower themselves from or to the floor.




Code: Sun-STP1
Single tier couch step (oval)

Code: Sun-STP2
Double tier couch step (oval)


Code: Sun-STP3
Single tier couch step (rectangular)



Sun-STP1
16cm(H) 54cm(W) 32cm(D)

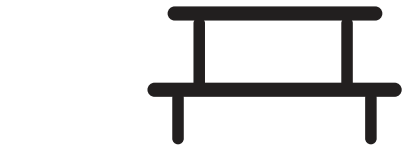

Sun-STP2
40cm(H) 57cm(W) 48cm(D)



Sun-STP3
22cm(H) 53cm(W) 31cm(D)



Maximum user
weight (26st.)



Anti-slip surfaces



Anti-tipping designs



Rubber feet
prevent movement



Unique oval design
(STP1 & STP2 only)

* Please contact us regarding availability for quantities of 10 or more.

*Stable
contemporary
support*

Healthcare Certified Chairs & Stools

Our healthcare-certified chairs and stools deliver durable, hygienic, and ergonomic comfort for patients, visitors, and professionals, featuring easy-to-clean materials, adjustable support, and long-lasting performance across diverse healthcare settings.

Individual Seating



Neptune Chairs
Pg 19



Venus Chairs
Pg 20



Galaxy Chairs
Pg 20



Aurora Chairs
Pg 21



Aurora Bariatric Chairs
Pg 21



Atlas Chairs
Pg 22-23



Atlas Bariatric Chairs
Pg 24

Modular/Bench Seating



Neptune Visitor 3 Seat Modules
Pg 25



Neptune Visitor 4 Seat Modules
Pg 26



Neptune Visitor 5 Seat Modules
Pg 27



Venus Visitor 3 Seat Modules
Pg 28



Neptune Visitor 4 Seat Module
Pg 29



Neptune Visitor 5 Seat Module
Pg 30

Operator Chairs & Stools



Gas-lift Chairs
Pg 31



Gas-lift Stools
Pg 32-33

product information



Maximum user weight
(20st.) per seat




One piece
durable moulded
polypropylene seat




Comfortable foam
upholstered
seat pad



Comfortable and durable
moulded arm rests
(Sun-SEAT71 only)

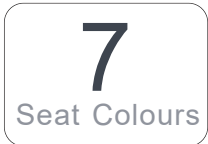

Available in 7
moulded seat colours

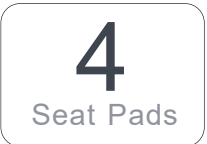

Choice of 4
upholstered seat pads



Anti-bacterial
vinyl and inter/vene
upholstery



Supplied fully
assembled



British made



Easy clean design
and materials



Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998




Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11
(composites)


moulded seat and upholstery colours - please specify when ordering

please see inside back cover for full specification of our upholstery range

moulded seat colours

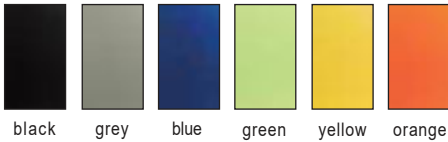

anti-bacterial vinyl upholstery colours - please specify when ordering

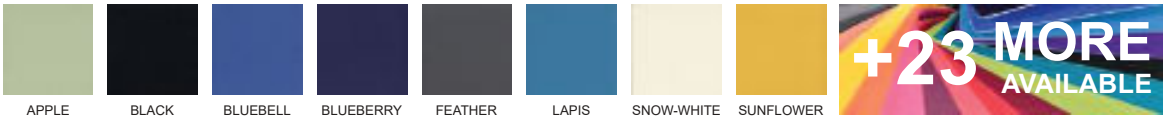

please see inside back cover for full range and specification of our vinyl range

Neptune Chairs

Neptune offers a contemporary range of chairs and seating modules, with great design and stability whilst being durable, hygienic and easy to clean. Available in 7 modern and vibrant colours with the option of upholstered seat pads in 2 anti-bacterial materials, vinyl and inter/vene™.



Code: Sun-SEAT70/ORANGE/VYL/BLK
Neptune visitor chair in orange with
black vinyl upholstered seat pad



Code: Sun-SEAT71/GREEN/VYL/GRY
Neptune visitor chair with arms in green
and grey vinyl upholstered seat pad



Code: Sun-SEAT70/YELLOW/BLK
Neptune visitor chair in yellow with black



Code: Sun-SEAT71/BLUE/GRY Neptune
visitor chair with arms in blue and grey

Sun-SEAT70



78cm(H)
47cm(W)
47cm(D)
Seat Height 47cm

Sun-SEAT71



78cm(H)
61cm(W)
47cm(D)
Seat Height 47cm

product information - venus visitor chairs



Maximum user weight (17st.)



Durable moulded polypropylene seats



Plastic feet prevent movement/protect floors



Safe and easy to stack (max. 10 high)



Supplied fully assembled



Chair linking kits available (please see below)



Easy clean design and materials




British made




Delivery normally within 3-5 working days *

* Please contact us regarding availability for quantities of 10 or more.


product information - galaxy visitor chairs




Maximum user weight (17st.)




Comfortable deep foam upholstery




Rounded/upholstered seat corners




Safe and easy to stack (max. 5 high)




Plastic feet prevent movement/protect floors




Vinyl or inter/vene Anti-bacterial upholstery




Supplied fully assembled




Supplied fully assembled




BSEN 16139: 2016 Furniture Strength, Durability and Safety Requirements




Anti-bacterial Vinyl Strength Test BSEN ISO 13934-1: 2013 BSEN ISO 13973-3: 2000 BSEN ISO 12947-2: 1998



Fire (Flammability) Testing Flammability BS7176 Medium Hazard BS5852 : 2006 Clause 11 (composites)




British made




Delivery normally within 3-5 working days *

* Please contact us regarding availability for quantities of 10 or more.


anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE

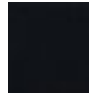


BLACK




BLUEBELL




BLUEBERRY




FEATHER




LAPIS



SNOW-WHITE



SUNFLOWER



please see inside back cover for full range and specification of our vinyl range

Venus Chairs


Multi-purpose economy chairs, ergonomically designed moulded seats and backs utilise modern materials, combining smooth surfaces, rounded corners with comfort, durability and ease of cleaning.



Code: Sun-SEAT1 - Venus visitor chair with moulded polypropylene seat


Code: Sun-SEAT2 - Venus visitor chair and arms with moulded polypropylene seat

Sun-SEAT1




79cm(H) 54cm(W) 56cm(D)
Seat Height 44cm(H)

moulded polypropylene

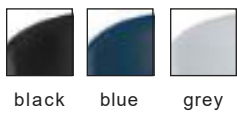


black blue grey

Sun-SEAT2



79cm(H) 58cm(W) 56cm(D)
Seat Height 44cm(H)

moulded polypropylene

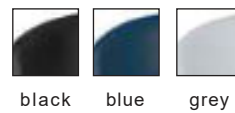

black blue grey

Galaxy Chairs

Multi-purpose upholstered economy chairs with comfortable deep foam ergonomically designed seat & back.



Sun-SEAT32



82cm(H)
54cm(W) 52cm(D)
Seat Height 44cm(H)


Sun-SEAT33




82cm(H)
58cm(W) 52cm(D)
Seat Height 44cm(H)

Code: Sun-SEAT32VYL - Galaxy visitor chair in anti-bacterial vinyl upholstery

Code: Sun-SEAT33VYL - Galaxy visitor chair with arms in anti-bacterial vinyl upholstery



Sun-SEAT/LINK1
Chair Linking Kit
for Chairs with No Arms



Sun-SEAT/LINK2
Chair Linking Kit
for Chairs with Arms

product information - aurora visitor chairs



Maximum user weight (17st.)



Comfortable deep foam upholstery



Vinyl or inter/vene Anti-bacterial upholstery




Solid lacquered beech wood frame




3 stackable

Safe and easy to stack (max. 3 high)

product information - aurora bariatric chair



Maximum user weight (34st.)



Comfortable deep foam upholstery



Solid beech arms



Solid beech arms



Ergonomic design

additional product information



Supplied fully assembled



Easy clean design And materials



BSEN 16139: 2016 Furniture Strength, Durability and Safety Requirements (Sun-SEAT35 & Sun-SEAT36 only)



Anti-bacterial Vinyl Strength Test BSEN ISO 13934-1: 2013 BSEN ISO 13973-3: 2000 BSEN ISO 12947-2: 1998



Fire (Flammability) Testing Flammability BS7176 Medium Hazard BS5852 : 2006 Clause 11 (composites)



British made



Delivery normally within 3-5 working days */**

* Please contact us regarding availability for quantities of 10 or more.

** Sun-SEAT35, Sun-SEAT36 & Sun-SEAT37 only.

anti-bacterial vinyl upholstery colours - please specify when ordering




APPLE




BLACK




BLUEBELL



BLUEBERRY




FEATHER



LAPIS



SNOW-WHITE



SUNFLOWER



please see inside back cover for full range and specification of our vinyl range

Aurora Chairs


Deluxe beech frame upholstered chairs with comfortable deep foam padded seat & back, with or without arms.



Code: Sun-SEAT35VYL - Aurora visitor chair with anti-bacterial vinyl upholstery


Code: Sun-SEAT36VYL - Aurora visitor chair and arms with anti-bacterial vinyl upholstery

Sun-SEAT35



84cm(H)
51cm(W)
57cm(D)
Seat Height 45cm

Sun-SEAT36



84cm(H)
53cm(W)
57cm(D)
Seat Height 45cm

Aurora Bariatric Chair

Our Aurora bariatric chair combines strength, stability and comfort with ease of access and is easy to clean.



Sun-CHA75



88cm(H)
84cm(W)
58cm(D)
Seat Width 72cm
Seat Height 45cm

Code: Sun-CHA75VYL - Aurora Bariatric 222kg (34st.) arm chair with anti-bacterial vinyl upholstery

product information

Maximum user weight (27.5st.)

Comfortable deep foam upholstery

Solid beech frame

Solid beech arms

Removable seat pad for cleaning

Optional pressure relieving foam upgrade

Optional drop-arms to aid lateral transfer (mid & high back chairs)

Wings for head support (Sun-CHA52 only)

Optional house keeping wheels (fitted to rear legs only)

Vinyl or inter/vene Anti-bacterial upholstery

British made

Supplied fully assembled

Easy clean design and materials

BSEN 16139: 2016 Furniture Strength, Durability and Safety Requirements

Anti-bacterial Vinyl Strength Test BSEN ISO 13934-1: 2013 BSEN ISO 13973-3: 2000 BSEN ISO 12947-2: 1998

Fire (Flammability) Testing Flammability BS7176 Medium Hazard BS5852 : 2006 Clause 11 (composites)

Anti-bacterial Vinyl

accessories * these accessories can only be fitted when ordered with the chair

Sun-CHA48/PRF Pressure Relieving Foam Upgrade for Seat & Back

Sun-CHA49/PRF Pressure Relieving Foam Upgrade for Seat & Back

Sun-CHA50/PRF Pressure Relieving Foam Upgrade for Seat & Back

Sun-CHA50/HKW* House Keeping Wheels (fitted to rear legs only)

anti-bacterial vinyl upholstery colours - please specify when ordering

APPLE

BLACK

BLUEBELL

BLUEBERRY

FEATHER

LAPIS

SNOW-WHITE

SUNFLOWER

please see inside back cover for full range and specification of our vinyl range

Atlas Chairs

Deluxe solid beech frame upholstered seating range with deep foam padded seat and back for optimum comfort.



Code: Sun-CHA48VYL - Low back side chair with anti-bacterial vinyl upholstery



Code: Sun-CHA49VYL - Low back arm chair with anti-bacterial vinyl upholstery



Code: Sun-CHA50VYL - Mid back arm chair with anti-bacterial vinyl upholstery



Code: Sun-CHA50DAVYL - Mid back arm chair with drop arms in anti-bacterial vinyl upholstery

Sun-CHA48

94cm(H)
54cm(W)
63cm(D)

Sun-CHA49

94cm(H)
60cm(W)
60cm(D)

Sun-CHA50

104cm(H)
60cm(W)
60cm(D)

Sun-CHA50DA

104cm(H)
60cm(W)
60cm(D)

Seat Heights 48cm - Internal Seat Widths 49cm

Atlas Chairs

Maximum user weight (27.5st.)

Solid beech frame

Solid beech arms

Removable seat pad for cleaning

Optional pressure relieving foam upgrade

Optional drop-arms to aid lateral transfer (mid & high back chairs)

Wings for head support (Sun-CHA52 only)

Optional house keeping wheels (fitted to rear legs only)

Vinyl or inter/vene Anti-bacterial upholstery

Supplied fully assembled

Easy clean design and materials

BSEN 16139: 2016 Furniture Strength, Durability and Safety Requirements

Anti-bacterial Vinyl Strength Test BSEN ISO 13934-1: 2013 BSEN ISO 13973-3: 2000 BSEN ISO 12947-2: 1998

Fire (Flammability) Testing Flammability BS7176 Medium Hazard BS5852 : 2006 Clause 11 (composites)

British made

accessories * these accessories can only be fitted when ordered with the chair

Sun-CHA51/PRF Pressure Relieving Foam Upgrade for Seat & Back

Sun-CHA52/PRF Pressure Relieving Foam Upgrade for Seat & Back

Sun-CHA50/HKW* House Keeping Wheels (fitted to rear legs only)

anti-bacterial vinyl upholstery colours - please specify when ordering

APPLE

BLACK

BLUEBELL

BLUEBERRY

FEATHER

LAPIS

SNOW-WHITE

SUNFLOWER

please see inside back cover for full range and specification of our vinyl range

Atlas Chairs

Deluxe solid beech frame upholstered seating range with deep foam padded seat and back for optimum comfort.

Code: Sun-CHA51VYL - High back arm chair with anti-bacterial vinyl upholstery

Code: Sun-CHA51DAVYL - High back arm chair with drop arms in anti-bacterial vinyl upholstery

Code: Sun-CHA52VYL - High back arm chair with wings in anti-bacterial vinyl upholstery

Code: Sun-CHA52DAVYL - High back arm chair with wings and drop arms in anti-bacterial vinyl upholstery

Sun-CHA51

117cm(H)
60cm(W)
60cm(D)

Sun-CHA51DA

117cm(H)
60cm(W)
60cm(D)

Sun-CHA52

118cm(H)
60cm(W)
60cm(D)

Sun-CHA52DA

118cm(H)
60cm(W)
60cm(D)

Seat Heights 48cm - Internal Seat Widths 49cm

product information

Maximum user weight (40st.)

Comfortable deep foam upholstery

Solid beech frame

Solid beech arms

Removable seat pad for cleaning

Optional pressure relieving foam upgrade

Optional drop-arms to aid lateral transfer (mid & high back chairs)

Wings for head support (Sun-CHA56 only)

Optional house keeping wheels (fitted to rear legs only)

Vinyl or inter/vene Anti-bacterial upholstery

Supplied fully assembled

Easy clean design and materials

Independently tested by FIRA

BSEN 16139: 2016 Furniture Strength, Durability and Safety Requirements

Independently tested by FIRA

Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998

Anti-bacterial Vinyl

Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)

British made

accessories * these accessories can only be fitted when ordered with the chair

Sun-CHA55/PRF Pressure Relieving Foam Upgrade for Seat & Back

Sun-CHA56/PRF Pressure Relieving Foam Upgrade for Seat & Back

Sun-CHA50/HKW* House Keeping Wheels (fitted to rear legs only)

anti-bacterial vinyl upholstery colours - please specify when ordering

APPLE

BLACK

BLUEBELL

BLUEBERRY

FEATHER

LAPIS

SNOW-WHITE

SUNFLOWER

please see inside back cover for full range and specification of our vinyl range

Atlas Bariatric Chairs

Deluxe solid beech frame upholstered seating range with deep foam padded seat and back for optimum comfort.



Code: Sun-CHA55VYL - High back bariatric arm chair with anti-bacterial vinyl upholstery



Code: Sun-CHA55DAVYL - High back bariatric arm chair with drop arms in anti-bacterial vinyl upholstery



Code: Sun-CHA56VYL - High back bariatric arm chair with wings in anti-bacterial vinyl upholstery



Code: Sun-CHA56DAVYL - High back bariatric arm chair with wings and drop arms in anti-bacterial vinyl upholstery

Sun-CHA55

110cm(H)
93cm(W)
73cm(D)

Sun-CHA55DA

110cm(H)
93cm(W)
73cm(D)

Sun-CHA56

126cm(H)
93cm(W)
73cm(D)

Sun-CHA56DA

126cm(H)
93cm(W)
73cm(D)

Seat Heights 51cm - Internal Seat Widths 81cm

Neptune Visitor 3 Seat Modular Seating

product information

Maximum user weight
(20st.) per seat

One piece
durable moulded
polypropylene seats

Comfortable foam
upholstered
seat pads

Strong ergonomic
chrome arms

Stable ergonomic
chrome legs

Adjustable feet for
uneven floors

7

Seat Colours

2

Seat Pads

Vinyl or inter/vene
Anti-bacterial upholstery

easy assembly

Minimal assembly
required
(tools supplied)

installation

Optional installation
service available

West Yorkshire
Material
Testing Service

Anti-bacterial Vinyl

British Made
Quality Assured

eze


Easy clean design
and materials

FIRA

Independently
tested by FIRA

Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998


dimensions

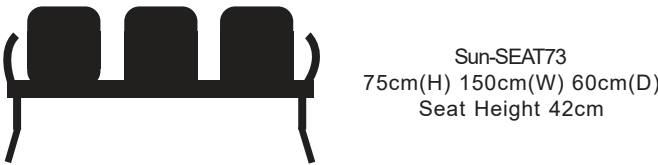

moulded seat and upholstery colours - please specify when ordering

please see inside back cover for full specification of our upholstery range

moulded seat colours

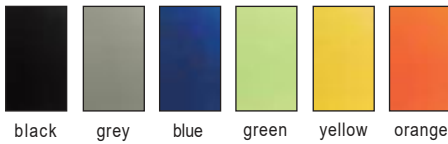

anti-bacterial vinyl upholstery colours - please specify when ordering



please see inside back cover for full range and specification of our vinyl range

Neptune 3 Seat Module

Neptune offers a contemporary range of seating modules and chairs, with great design and stability whilst being durable, hygienic and easy to clean. Available in 7 modern and vibrant colours with the option of upholstered seat pads in 2 anti-bacterial materials, vinyl and inter/vene™.



Code: Sun-SEAT73/3/BLACK
Neptune visitor 3 seat module in black - 2 chrome arms



Code: Sun-SEAT73/3/GREEN/VYL/BLK
Neptune visitor 3 seat module in green with 3 black
seat pads in anti-bacterial vinyl - 2 chrome arms



Code: Sun-SEAT73/3/YELLOW/GRY
Neptune visitor 3 seat module in yellow with 3 grey
seat pads in anti-bacterial inter/vene - 2 chrome arms

Neptune Visitor 3 Seat Modular Seating

product information

127kg

Maximum user weight (20st.) per seat

One piece durable moulded polypropylene seats

Comfortable foam upholstered seat pads

Strong ergonomic chrome arms

Stable ergonomic chrome legs

Adjustable feet for uneven floors

7
Seat Colours

2
Seat Pads

Vinyl or inter/vene
Anti-bacterial upholstery

easy assembly

Minimal assembly required (tools supplied)

installation

West Yorkshire
Materials
Testing Service
Anti-bacterial Vinyl

Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)

British Made
Quality Assured

eze

Easy clean design
and materials

FIRA
Independently
tested by FIRA

Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998


dimensions




moulded seat and upholstery colours - please specify when ordering

please see inside back cover for full specification of our upholstery range

moulded seat colours

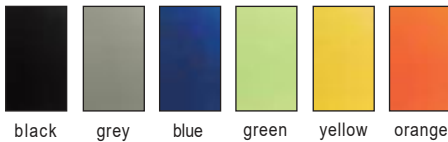

anti-bacterial vinyl upholstery colours - please specify when ordering

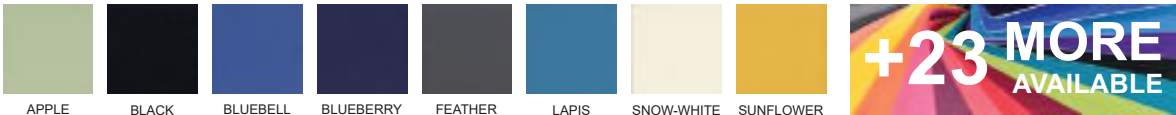

please see inside back cover for full range and specification of our vinyl range

Neptune 4 Seat Module

Neptune offers a contemporary range of seating modules and chairs, with great design and stability whilst being durable, hygienic and easy to clean. Available in 7 modern and vibrant colours with the option of upholstered seat pads in 2 anti-bacterial materials, vinyl and inter/vene™.



Code: Sun-SEAT74/4/ORANGE- Neptune visitor 4 seat module in orange - 2 chrome arms



Code: Sun-SEAT74/4/IVORY/VYL/GRY
Neptune visitor 4 seat module in ivory with 4 grey seat pads in anti-bacterial vinyl - 2 chrome arms



Code: Sun-SEAT74/4/BLUE/BLK
Neptune visitor 4 seat module in blue with 4 black seat pads - 2 chrome arms

Neptune Visitor 3 Seat Modular Seating

product information

Maximum user weight (20st.) per seat

One piece durable moulded polypropylene seats

Comfortable foam upholstered seat pads

Strong ergonomic chrome arms

Stable ergonomic chrome legs

Adjustable feet for uneven floors

7

Seat Colours

2

Seat Pads

Vinyl or inter/vene Anti-bacterial upholstery

easy assembly

Minimal assembly required (tools supplied)

installation

Optional installation service available

West Yorkshire Fire Service Testing Service

Anti-bacterial Vinyl

British Made Quality Assured

eze


Easy clean design and materials

FIRA

Independently tested by FIRA

Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998


dimensions

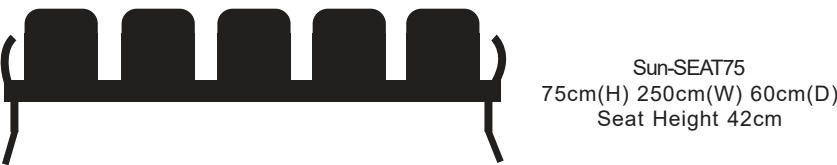

moulded seat and upholstery colours - please specify when ordering

please see inside back cover for full specification of our upholstery range

moulded seat colours

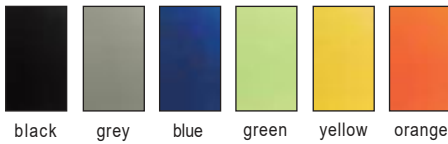

anti-bacterial vinyl upholstery colours - please specify when ordering



please see inside back cover for full range and specification of our vinyl range

Neptune 5 Seat Modules

Neptune offers a contemporary range of seating modules and chairs, with great design and stability whilst being durable, hygienic and easy to clean. Available in 7 modern and vibrant colours with the option of upholstered seat pads in 2 anti-bacterial materials, vinyl and inter/vene™.



Code: Sun-SEAT75/5/YELLOW - Neptune visitor 5 seat module in yellow - 2 chrome arms



Code: Sun-SEAT75/5/GREY/VYL/BLK
Neptune visitor 5 seat module in grey with 5 black seat pads in anti-bacterial vinyl - 2 chrome arms



Code: Sun-SEAT75/5/ORANGE/GRY
Neptune visitor 5 seat module in orange with 5 grey seat pads - 2 chrome arms

product information



Maximum user weight
(17st.) per seat



Durable moulded
polypropylene seats



Comfortable deep
foam upholstered seats



Rounded/upholstered
seat corners



Vinyl or inter/vene
Anti-bacterial upholstery




Adjustable feet for
uneven floors



Durable silver epoxy
coated frame



British made

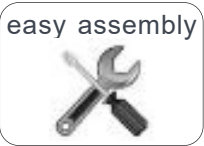

easy assembly
Minimal assembly
required (tools supplied)



Optional installation
service available



Easy clean design
and materials




Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11
(composites)

dimensions

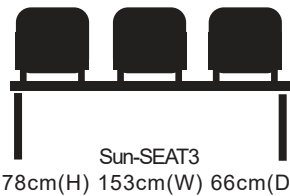

upholstery colours - please specify when ordering


please see inside back cover for full specification of our upholstery range

polypropylene seats



black blue grey

anti-bacterial vinyl upholstery colours - please specify when ordering

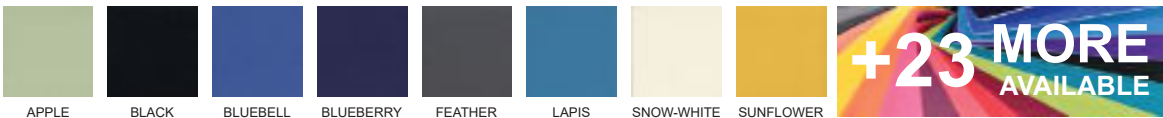

please see inside back cover for full range and specification of our vinyl range

Venus 3 Seat Modules

Contemporary range of 3 seat seating modules, which are durable, hygienic and easy to clean. Available in a wide range of colours and anti-bacterial fabrics.



Code: Sun-SEAT3 - Venus visitor 3 seat module with moulded polypropylene seats



Code: Sun-SEAT3VYL - Venus visitor 3 seat module with anti-bacterial vinyl upholstery
Code: Sun-SEAT3TVYL - Venus visitor 2 seat & 1 table module with anti-bacterial vinyl upholstery

Venus Visitor 4 Seat Modular Seating

product information



Maximum user weight (17st.) per seat



Durable moulded polypropylene seats



Comfortable deep foam upholstered seats



Rounded/upholstered seat corners



Vinyl or inter/vene Anti-bacterial upholstery



Adjustable feet for uneven floors



Durable silver epoxy coated frame



Magazine table option



easy assembly
Minimal assembly required (tools supplied)



installation
Optional installation service available



British made



Easy clean design and materials




Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)
Anti-bacterial Vinyl


dimensions

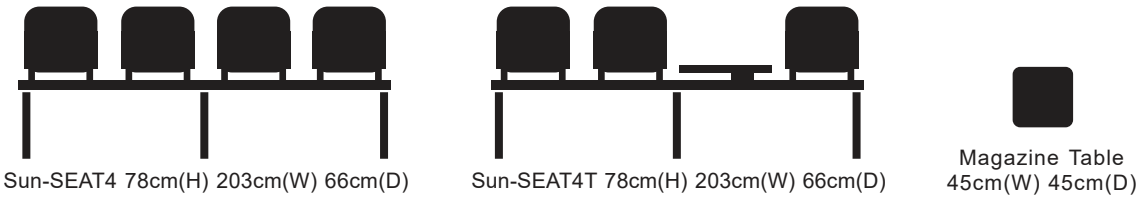


upholstery colours - please specify when ordering

please see inside back cover for full specification of our upholstery range

polypropylene seats

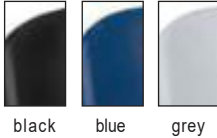

anti-bacterial vinyl upholstery colours - please specify when ordering

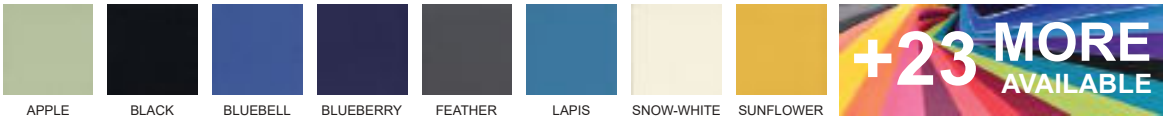

please see inside back cover for full range and specification of our vinyl range

Venus Visitor 4 Seat Modular Seating

Contemporary range of 4 seat seating modules, which are durable, hygienic and easy to clean. Available in a wide range of colours and anti-bacterial fabrics.



Code: Sun-SEAT4 - Venus visitor 4 seat module with moulded polypropylene seats
Code: Sun-SEAT4T- Venus visitor 3 seat & 1 table module with moulded polypropylene seats



Code: Sun-SEAT4VYL - Venus visitor 4 seat module with anti-bacterial vinyl upholstery
Code: Sun-SEAT4TVYL - Venus visitor 3 seat & 1 table module with anti-bacterial vinyl upholstery

Venus Visitor 5 Seat Modular Seating

product information

Maximum user weight (17st.) per seat

Durable moulded polypropylene seats

Comfortable deep foam upholstered seats

Rounded/upholstered seat corners

Vinyl or inter/vene Anti-bacterial upholstery

Adjustable feet for uneven floors

Durable silver epoxy coated frame

Magazine table option

easy assembly

installation

British made

Easy clean design and materials

Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998

Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)


Venus Visitor 5 Seat Modular Seating

Contemporary range of 5 seat seating modules, which are durable, hygienic and easy to clean. Available in a wide range of colours and anti-bacterial fabrics.




Code: Sun-SEAT5 - Venus visitor 5 seat module with moulded polypropylene seats
Code: Sun-SEAT5T- Venus visitor 4 seat & 1 table module with moulded polypropylene seats


dimensions

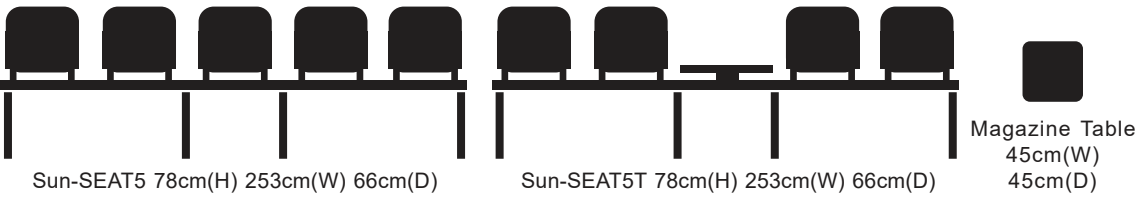

upholstery colours - please specify when ordering
please see inside back cover for full specification of our upholstery range
polypropylene seats

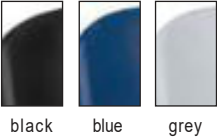

anti-bacterial vinyl upholstery colours - please specify when ordering

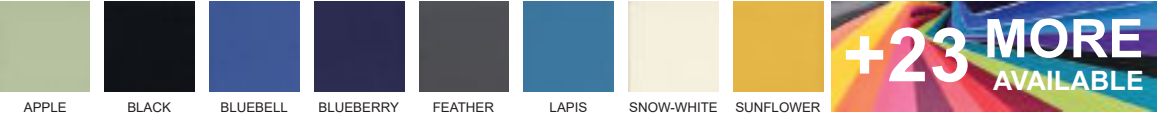


please see inside back cover for full range and specification of our vinyl range




Code: Sun-SEAT5VYL - Venus visitor 5 seat module with anti-bacterial vinyl upholstery
Code: Sun-SEAT5TVYL - Venus visitor 4 seat & 1 table module with anti-bacterial vinyl upholstery

Gas-lift Chairs


product information




Maximum user weight (17st.)




Anti-bacterial/ anti-microbial vinyl




Comfortable padded seat and back




40mm Deep foam for support and comfort




Hygienic rounded corners




Strong and durable ergonomic arms (CHA2 & CHA4)




Adjustable height chrome foot ring provides leg support (CHA3, CHA4 & CHA5)



Gas-assisted adjustment



Five castor base for stability and manoeuvrability



Glides reduce lateral movement (CHA5 only)

HEIGHT RANGE

42-57cm

Sun-CHA1

HEIGHT RANGE

42-57cm

Sun-CHA2

HEIGHT RANGE

48-68cm

Sun-CHA3

HEIGHT RANGE


48-68cm

Sun-CHA4


HEIGHT RANGE

55-81cm

Sun-CHA5




Easy clean design and materials



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)




British made




Delivery normally within 3-5 working days *

* Please contact us regarding availability for quantities of 10 or more.

dimensions

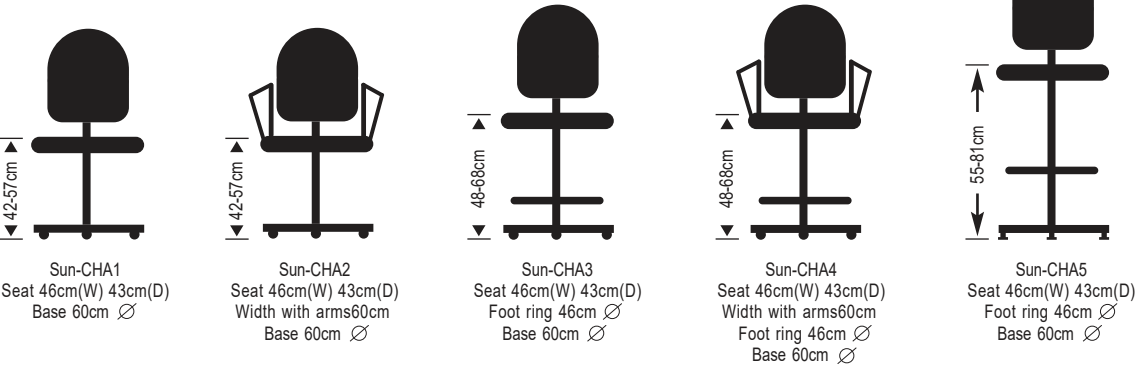

anti-bacterial vinyl upholstery colours - please specify when ordering

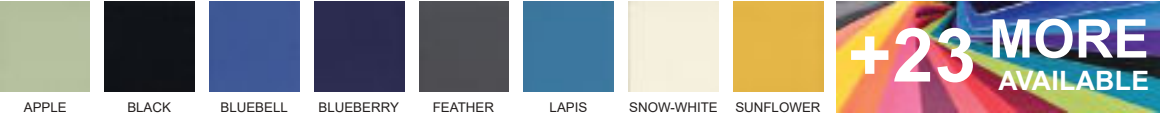

please see inside back cover for full range and specification of our vinyl range

Gas-Lift Chairs

Adjustable for support and comfort, Operator Chairs are designed for everyday general healthcare use, therefore are comfortable, strong, durable and easy to clean.



Code: Sun-CHA1 - Gas-lift chair



Code: Sun-CHA2 - Gas-lift chair with arms



Code: Sun-CHA3 - Gas-lift chair & foot ring



Code: Sun-CHA4 - Gas-lift chair with arms & foot ring



Code: Sun-CHA5 - High level gas-lift chair & foot ring

Gas-lift Stools

product information

110kg

Maximum user weight (17st.)

Anti-bacterial/ anti-microbial vinyl

Luxury 65mm deep upholstery

2 layer foam for support and comfort

Adjustable height chrome foot ring (STO3 & STO3/GLIDES)

Gas-assisted adjustment

HEIGHT RANGE

42-58cm

Sun-STO2

HEIGHT RANGE

39-55cm

Sun-STO2/GLIDES

HEIGHT RANGE

51-69cm

Sun-STO3

HEIGHT RANGE

48-66cm

Sun-STO3/GLIDES

Five castor base for stability and manoeuvrability

Glides reduce lateral movement (STO2/GLIDES & STO3/GLIDES only)

Easy clean design and materials

Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998


Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)

British made


Delivery normally within 3-5 working days *

* Please contact us regarding availability for quantities of 10 or more.

dimensions

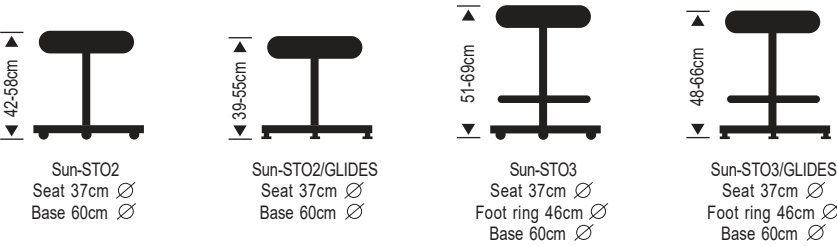

anti-bacterial vinyl upholstery colours - please specify when ordering

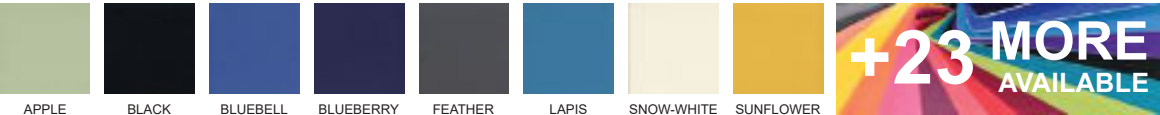

please see inside back cover for full range and specification of our vinyl range

Gas-lift Stools

Designed for everyday general healthcare use, Examination Stools can be static or mobile, are simple to adjust and comfortable, while being strong, durable and easy to clean.



Code: Sun-STO2
Gas-lift stool with 5 castor base



Code: Sun-STO2/GLIDES
Gas-lift stool with 5 glides



Code: Sun-STO3
Gas-lift stool & foot ring with 5 castor base



Code: Sun-STO3/GLIDES
Gas-lift stool & foot ring with 5 glides

Gas-lift Stools with Back Rests

product information

110kg

Maximum user weight (17st.)

Anti-bacterial/ anti-microbial vinyl

Comfortable padded back rest

Luxury 65mm deep seat with upholstered back rest

2 layer foam for support and comfort

Adjustable height chrome foot ring (Sun-STO5 & Sun-STO6)

Gas-assisted adjustment

HEIGHT RANGE 42-58cm

Sun-STO4

HEIGHT RANGE 53-79cm

Sun-STO5

HEIGHT RANGE 42-58cm

Sun-STO6

Five castor base for stability & manoeuvrability (Sun-STO4 & Sun-STO6)

Glides reduce lateral movement (Sun-STO5)

Easy clean design and materials

FIRA Independently tested by FIRA

Anti-bacterial Vinyl Strength Test
BSEN ISO 13934-1: 2013
BSEN ISO 13973-3: 2000
BSEN ISO 12947-2: 1998

West Yorkshire Materials Testing Service

Anti-bacterial Vinyl


Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11 (composites)

British made


Delivery normally within 3-5 working days *

* Please contact us regarding availability for quantities of 10 or more.

dimensions

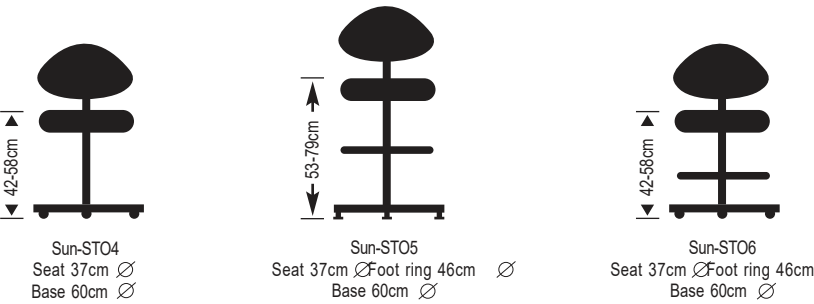

anti-bacterial vinyl upholstery colours - please specify when ordering

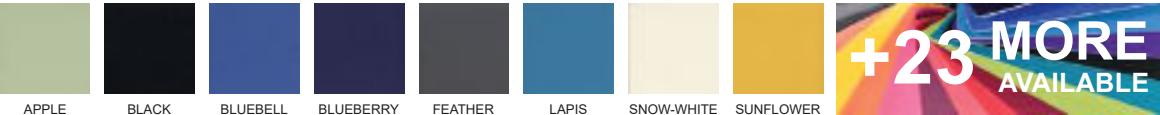

please see inside back cover for full range and specification of our vinyl range

Gas-lift Stools with Back Rests

Designed for everyday general healthcare use, Examination Stools with back rests can be static or mobile, are simple to adjust and comfortable, while being strong, durable and easy to clean.



Code: Sun-STO4
Gas -lift stool with back rest



Code: Sun-STO6
Gas -lift stool with back rest and foot ring



Code: Sun-STO5
High level gas -lift stool with back rest, foot ring and glides

*effective
ergonomic
design*

Medical Trolleys & Storage Solutions

Our British-made medical trolleys and storage solutions improve healthcare efficiency, offering versatile, secure, and durable options for medication, records, surgical, and dressing management, supporting organised and accessible workflows.

Vista Trolleys & Storage Modules



Vista Trolleys
Pg 35-41



Vista Clinical Trolleys
Pg 42-45



Vista Storage Modules
Pg 46-47

Surgical & Dressing Trolleys



Surgical Trolleys
Pg 48-49



Surgical Trolleys
(Mirror Polished)
Pg 50-51



Dressing Trolleys
Pg 52 - 53

Ward Drug & Medicine Dispensing Trolley



EPMA Ward Drug & Medicine
Dispenser Trolley
Pg 54



Ward Drug & Medicine
Dispenser Trolley
Pg 55

Dosage System Trolleys



Unit Dosage System (UDS)
Trolley
Pg 56



Monitor Dosage System (MDS)
Trolley
Pg 57

Medical Notes & Records Trolleys



Medical Notes Trolley
Pg 58



Records Transfer Trolley
Pg 59

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



100mm castors with front brakes (Vista 15 only)



Optional tamper resistant locking door (Vista 15 only)



Extra low level (68cm), ideal for under couches (Vista 15 only)



Stainless steel handle (Vista 15 only)



Removable top shelf with integrated upstand (Vista 15 only)



Easy clean design and materials



British Made

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors (front 2 braked)



Sun-BUF2 *
4 Protective Buffers (Vista 10 Trolleys)



Sun-BUF3 *
4 Protective Buffers (Vista 15 Trolleys)



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 */**
Sharps Box Bracket (Daniels - All Sizes)



Sun-MPT/BIN *
2.5 Litre Bin (Requires Medi-rail)



Sun-MPT/Rail1 *
Medi-rail



Sun-MPT/AGH1 */**
Alcohol Gel Holder



Sun-GBH1 */**
Glove Box Holder



Sun-MPT51/LD */***
Tamper Resistant Locking Door (Vista 15 Trolleys)



Sun-DIV1
Single Depth Tray Dividers (per set)
Vista 10 & 15



Sun-DIV4
Double Depth Tray Dividers (per set)
Vista 15





Sun-MPT51/CV
Trolley Cover (Vista 15 Trolleys)




Sun-MPT1/CV
Trolley Cover (Vista 10 Trolleys)

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions

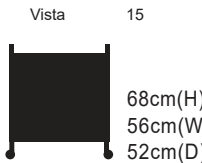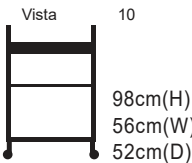

Single Depth Trays
8cm(H) 47cm(W)
42cm(D)
Double Depth Trays
16cm(H) 47cm(W)
42cm(D)

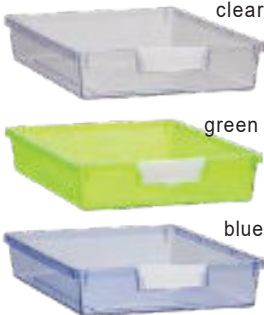

tray colours

please specify when ordering
(Tray labels and covers included)

Vista 10 & 15 Trolleys

Vista 10 is an economy priced trolley without any compromise on quality, while Vista 15 is designed to fit under a couch, both trolleys are functional, durable and easy to clean.



Code: Sun-TRO1
Vista 10 Trolley - 1 tray & 1 shelf



Code: Sun-MPT51
Vista 15 Trolley - 4 single depth trays



Code: Sun-MPT51/LD
Locking door for Vista 15 trolley



Code: Sun-MPT52
Vista 15 Trolley - 2 single/1 double depth trays



Code: Sun-MPT53
Vista 15 Trolley - 2 double depth trays

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



Narrow width (39cm) compact design



100mm castors with front brakes



Optional tamper resistant locking door



Removable top shelf with integrated upstand



Easy clean design and materials



Supplied fully assembled



British made

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors
(front 2 braked)



Sun-BUF2 *
4 Protective Buffers



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 */**
Sharps Box Bracket
(Daniels - All Sizes)



Sun-MPT/AGH1 */**
Alcohol Gel Holder



Sun-GBH1 */**
Glove Box Holder



Sun-MPT/Rail1 *
Medi-rail



Sun-MPT/BIN *
2.5 Litre Bin
(Requires Medi-rail)



Sun-DIV2
Double Depth
Tray Dividers (per set)



Sun-DIV3
Single Depth
Tray Dividers (per set)




Sun-MPT5/LD */***
Tamper Resistant
Locking Door




Sun-MPT5/CV
Trolley Cover

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions

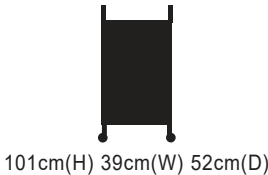

101cm(H) 39cm(W) 52cm(D)

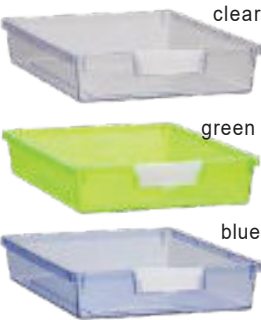

tray colours
please specify when
ordering
(Tray labels and covers
included)

Single Depth Trays
8cm(H) 31cm(W) 42cm(D)
Double Depth Trays
16cm(H) 31cm(W) 42cm(D)

Vista 30 Trolleys

Designed for environments with limited space, Vista 30 provides standard working height with compact design, optimum storage, and is functional, durable and easy to clean.



Code: Sun-MPT10
Vista 30 Trolley - 6 single depth trays



Code: Sun-MPT5
Vista 30 Trolley - 2 single/2 double depth trays



Code: Sun-MPT6
Vista 30 Trolley - 3 double depth trays



Code: Sun-MPT5/LD
Locking door for Vista 30 trolley

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



Low level (80cm), ideal for under worktops



100mm castors with front brakes



Optional tamper resistant locking door



Removable top shelf with integrated upstand



Stainless steel handle



Easy clean design and materials



British made

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors
(front 2 braked)



Sun-BUF3 *
4 Protective Buffers



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 */**
Sharps Box Bracket
(Daniels - All Sizes)



Sun-MPT/AGH1 */**
Alcohol Gel Holder



Sun-DIV1
Single Depth
Tray Dividers (per set)



Sun-DIV4
Double Depth
Tray Dividers (per set)



Sun-GBH1 */**
Glove Box Holder



Sun-MPT/Rail1 *
Medi-rail



Sun-MPT/BIN *
2.5 Litre Bin
(Requires Medi-rail)



Sun-MPT40/LD */***
Tamper Resistant
Locking Door



Sun-MPT40/CV
Trolley Cover

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions



80cm(H) 56cm(W) 52cm(D)



clear



green



blue

tray colours

please specify when ordering
(Tray labels and covers included)

Single Depth Trays
8cm(H) 47cm(W) 42cm(D)
Double Depth Trays
16cm(H) 47cm(W) 42cm(D)

Vista 40 Trolleys

Designed to fit under a worktop, Vista 40 can be stored neatly away while providing optimum number of trays for storage, and is functional, durable and easy to clean.



Code: Sun-MPT40 Vista 40 Trolley
5 single depth trays



Code: Sun-MPT41 Vista 40 Trolley
3 single/1 double depth trays



Code: Sun-MPT42 Vista 40 Trolley
1 single/2 double depth trays



Code: Sun-MPT40/LD
Locking door for Vista 40 trolley

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



100mm castors with front brakes



Optional tamper resistant locking door



Removable top shelf with integrated upstand



Easy clean design and materials



Supplied fully assembled



British made

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors (front 2 braked)



Sun-BUF2 *
4 Protective Buffers



Sun-MPT/OXY1 */**
Oxygen Bottle Holder (Size D)



Sun-MPT/AGH1 */**
Alcohol Gel Holder



Sun-GBH1 */**
Glove Box Holder



Sun-MPT/DRIP1 */**
Adjustable Height Drip Pole



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 */**
Sharps Box Bracket (Daniels - All Sizes)



Sun-MPT/BIN *
2.5 Litre Bin (Requires Medi-rail)



Sun-MPT1/LD */***
Tamper Resistant Locking Door



Sun-DIV1
Single Depth Tray Dividers (per set)



Sun-DIV4
Double Depth Tray Dividers (per set)



Sun-MPT/Rail1 *
Medi-rail




Sun-MPT1/CV
Trolley Cover

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions



101cm(H) 56cm(W) 52cm(D)

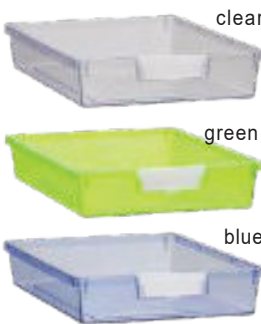

tray colours
please specify when ordering
(Tray labels and covers included)

Single Depth Trays
8cm(H) 47cm(W) 42cm(D)
Double Depth Trays
16cm(H) 47cm(W) 42cm(D)

Vista 50 Trolleys

Suitable for everyday general healthcare use, Vista 50 provides standard working height with generous working space, optimum tray storage, and is functional, durable and easy to clean.



Code: Sun-MPT1
Vista 50 Trolley - 6 single depth trays



Code: Sun-MPT2
Vista 50 Trolley - 4 single/1 double depth trays



Code: Sun-MPT13
Vista 50 Trolley - 3 double depth trays



Code: Sun-MPT1/LD
Locking door for Vista 50 trolley

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



100mm castors with front brakes



Optional tamper resistant locking door



Removable top shelf with integrated upstand



Easy clean design and materials



Supplied fully assembled



British made

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors (front 2 braked)



Sun-BUF2 *
4 Protective Buffers



Sun-MPT/OXY1 */**
Oxygen Bottle Holder (Size D)



Sun-MPT/AGH1 */**
Alcohol Gel Holder



Sun-GBH1 */**
Glove Box Holder



Sun-MPT/DRIP1 */**
Adjustable Height Drip Pole



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 */**
Sharps Box Bracket (Daniels - All Sizes)



Sun-MPT/BIN *
2.5 Litre Bin (Requires Medi-rail)



Sun-MPT55/LD */**
Tamper Resistant Locking Door



Sun-DIV1
Single Depth Tray Dividers (per set)



Sun-DIV4
Double Depth Tray Dividers (per set)



Sun-MPT/Rail1 *
Medi-rail



Sun-MPT55/CV
Trolley Cover

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions



110cm(H) 56cm(W) 52cm(D)



clear



green



blue

tray colours

please specify when ordering (Tray labels and covers included)

Single Depth Trays
8cm(H) 47cm(W) 42cm(D)
Double Depth Trays
16cm(H) 47cm(W) 42cm(D)

Vista 55 Trolleys

Suitable for everyday general healthcare use, Vista 55 provides standard working height with generous working space, smart design optimises storage with upto seven trays, and is functional, durable and easy to clean.



Code: Sun-MPT55
Vista 55 Trolley - 7 single depth trays



Code: Sun-MPT56
Vista 55 Trolley - 3 single/2 double depth trays



Code: Sun-MPT57
Vista 55 Trolley - 1 single/3 double depth trays



Code: Sun-MPT55/LD
Locking door for Vista 55 trolley

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



100mm castors with front brakes



Optional tamper resistant locking doors



Removable top shelf with integrated upstand



Easy clean design and materials



Supplied fully assembled



British made

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors
(front 2 braked)



Sun-BUF2 *
4 Protective Buffers



Sun-MPT/OXY1 */**
Oxygen Bottle Holder
(Size D)



Sun-MPT/BIN *
2.5 Litre Bin
(Requires Medi-rail)



Sun-MPT60/LD */**
2 Tamper Resistant
Locking Doors



Sun-MPT/DRIP1 */**
Adjustable Height
Drip Pole



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 */**
Sharps Box Bracket
(Daniels - All Sizes)



Sun-MPT/AGH1 */**
Alcohol Gel Holder



Sun-GBH1 */**
Glove Box Holder



Sun-MPT/Rail1 *
Medi-rail



Sun-DIV2
Double Depth
Tray Dividers (per set)



Sun-DIV3
Single Depth
Tray Dividers (per set)



Sun-MPT60/CV
Trolley Cover

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions



101cm(H) 70cm(W) 52cm(D)



clear



green



blue

tray colours

please specify when ordering
(Tray labels and covers included)

Single Depth Trays
8cm(H) 31cm(W) 42cm(D)
Double Depth Trays
16cm(H) 31cm(W) 42cm(D)

Vista 60 Trolleys

Designed to optimise storage where space is limited, Vista 60 provide additional working space and compact storage with up to 12 trays, and is functional, durable and easy to clean.



Code: Sun-MPT60
Vista 60 Trolley - 12 single depth trays



Code: Sun-MPT61
Vista 60 Trolley - 6 single/3 double depth trays



Code: Sun-MPT62
Vista 60 Trolley - 6 double depth trays



Code: Sun-MPT60/LD
2 Locking doors for Vista 60 trolley

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



100mm castors with front brakes



Optional tamper resistant locking doors



Removable top shelf with integrated upstand



Easy clean design and materials



Supplied fully assembled



British made

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors
(front 2 braked)



Sun-BUF2 *
4 Protective Buffers



Sun-MPT/OXY1 */**
Oxygen Bottle Holder
(Size D)



Sun-MPT/BIN *
2.5 Litre Bin
(Requires Medi-rail)



Sun-MPT4/LD */***
2 Tamper Resistant
Locking Doors



Sun-MPT/DRIP1 */**
Adjustable Height
Drip Pole



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 */**
Sharps Box Bracket
(Daniels - All Sizes)



Sun-MPT/AGH1 */**
Alcohol Gel Holder



Sun-GBH1 */**
Glove Box Holder



Sun-MPT/Rail1 *
Medi-rail



Sun-DIV1
Single Depth
Tray Dividers (per set)



Sun-DIV4
Double Depth
Tray Dividers (per set)



Sun-MPT4/CV
Trolley Cover

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions



101cm(H) 103cm(W) 52cm(D)



Single Depth Trays
8cm(H) 47cm(W) 42cm(D)
Double Depth Trays
16cm(H) 47cm(W) 42cm(D)

tray colours

please specify when ordering
(Tray labels and covers included)

Vista 100 Trolleys

Designed for the most demanding healthcare environments, Vista 100 combines maximum working space with optimum tray storage, and functional, durable and easy to clean.



Code: Sun-MPT16
Vista 100 Trolley - 12 single depth trays



Code: Sun-MPT4
Vista 100 Trolley - 8 single/2 double depth trays



Code: Sun-MPT17
Vista 100 Trolley - 6 single/3 double depth trays



Code: Sun-MPT4/LD
2 Locking doors for Vista 100 trolley

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



Removable mirror polished recessed stainless steel tray (304 grade)



Narrow width (39cm) compact design



100mm castors with front brakes



4 Protective buffers as standard



Optional tamper resistant locking door



Easy clean design and materials



British made



Supplied fully assembled

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors (front 2 braked)



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 */**
Sharps Box Bracket (Daniels - All Sizes)



Sun-MPT/BIN *
2.5 Litre Bin (Requires Medi-rail)



Sun-CVT30/LD */***
Tamper Resistant Locking Door



Sun-MPT/AGH1 */**
Alcohol Gel Holder



Sun-GBH1 */**
Glove Box Holder



Sun-MPT/Rail1 *
Medi-rail



Sun-DIV2
Double Depth
Tray Dividers (per set)




Sun-DIV3
Single Depth
Tray Dividers (per set)




Sun-MPT5/CV
Trolley Cover

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions

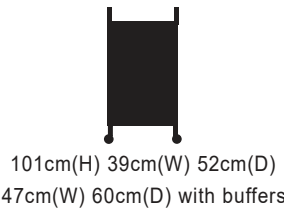

Single Depth Trays
8cm(H) 31cm(W) 42cm(D)
Double Depth Trays
16cm(H) 31cm(W) 42cm(D)

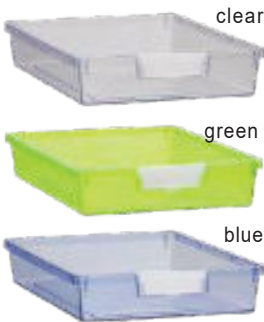

tray colours
please specify when
ordering
(Tray labels and covers
included)

Vista 30 Clinical Trolleys

Designed for environments with limited space, Vista 30 Clinical provides standard working height with compact design, optimum storage, a removable recessed stainless steel tray and is functional, durable and easy to clean.



Code: Sun-CVT30
Vista 30 Clinical trolley - 6 single depth trays



Code: Sun-CVT31
Vista 30 Clinical trolley -
2 single/2 double depth trays



Code: Sun-CVT32
Vista 30 Clinical trolley - 3 double depth trays



Code: Sun-CVT30/LD
Locking door for Vista 30 clinical trolley

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



Removable mirror polished recessed stainless steel tray (304 grade)



Low level (80cm), ideal for under worktops



100mm castors with front brakes



Optional tamper resistant locking door



Stainless steel handle



4 Protective buffers as standard



Easy clean design and materials



Supplied fully assembled



British made

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors
(front 2 braked)



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 **
Sharps Box Bracket
(Daniels - All Sizes)



Sun-MPT/BIN *
2.5 Litre Bin
(Requires Medi-rail)



Sun-CVT40/LD ***
Tamper Resistant
Locking Door



Sun-MPT/AGH1 ***
Alcohol Gel Holder



Sun-GBH1 **
Glove Box Holder



Sun-MPT/Rail1 *
Medi-rail



Sun-DIV1
Single Depth
Tray Dividers (per set)




Sun-DIV4
Double Depth
Tray Dividers (per set)




Sun-MPT40/CV
Trolley Cover

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions

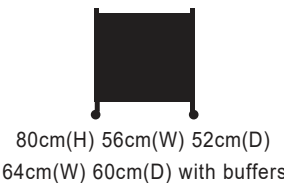

80cm(H) 56cm(W) 52cm(D)
64cm(W) 60cm(D) with buffers

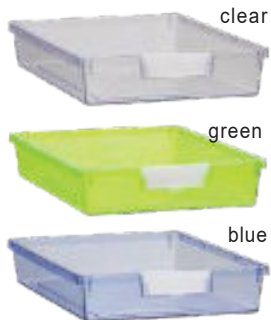

tray colours
please specify when ordering
(Tray labels and covers included)

Single Depth Trays
8cm(H) 47cm(W) 42cm(D)
Double Depth Trays
16cm(H) 47cm(W) 42cm(D)

Vista 40 Clinical Trolleys

Designed to fit under a worktop, Vista 40 Clinical can be stored neatly away, while providing optimum number of trays for storage, removable recessed stainless steel tray and is functional, durable and easy to clean.



Code: Sun-CVT40
Vista 40 Clinical trolley
5 single depth trays



Code: Sun-CVT41
Vista 40 Clinical trolley
3 single/1 double depth tray



Code: Sun-CVT42
Vista 40 Clinical trolley
1 single/2 double depth trays



Code: Sun-CVT40/LD
Locking door for Vista 40 clinical trolley

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



Removable mirror polished recessed stainless steel tray (304 grade)



100mm castors with front brakes



Optional tamper resistant locking door



4 Protective buffers as standard



Easy clean design and materials



British made



Supplied fully assembled

accessories

* These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors
(front 2 braked)



Sun-MPT/OXY1 **
Oxygen Bottle Holder
(Size D)



Sun-MPT/DRIP1 **
Adjustable Height
Drip Pole



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 **
Sharps Box Bracket
(Daniels - All Sizes)



Sun-MPT/AGH1 **
Alcohol Gel Holder



Sun-GBH1 **
Glove Box Holder



Sun-MPT/Rail1 *
Medi-rail



Sun-MPT/BIN *
2.5 Litre Bin
(Requires Medi-rail)



Sun-CVT50/LD ***
Tamper Resistant
Locking Door



Sun-DIV4
Double Depth
Tray Dividers (per set)




Sun-DIV1
Single Depth
Tray Dividers (per set)




Sun-MPT1/CV
Trolley Cover

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions

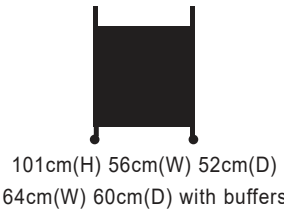

101cm(H) 56cm(W) 52cm(D)
64cm(W) 60cm(D) with buffers

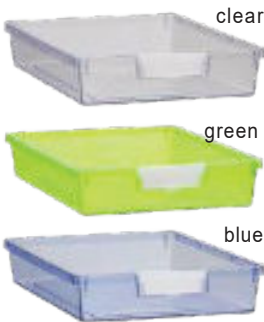

tray colours
please specify when
ordering
(Tray labels and covers
included)

Single Depth Trays
8cm(H) 47cm(W) 42cm(D)
Double Depth Trays
16cm(H) 47cm(W) 42cm(D)

Vista 50 Clinical Trolleys

Suitable for everyday general healthcare use, Vista 50 Clinical provides standard working height with generous working space, optimum tray storage, removable recessed stainless steel tray and is functional, durable and easy to clean.



Code: Sun-CVT50
Vista 50 Clinical trolley - 6 single depth trays



Code: Sun-CVT51
Vista 50 Clinical trolley
4 single/1 double depth tray



Code: Sun-CVT52
Vista 50 Clinical Trolley - 3 double depth trays



Code: Sun-CVT50/LD
Locking door for Vista 50 clinical trolley

product information



Tilt & locking tray runner system



Trays integrate with all other Vista products



Removable mirror polished recessed stainless steel tray (304 grade)



100mm castors with front brakes



Optional tamper resistant locking doors



4 Protective buffers as standard



Easy clean design and materials



British made



Supplied fully assembled

accessories

* These accessories can only be fitted when ordered with the trolley



Sun-LC3A *
4 Antistatic Castors (front 2 braked)



Sun-MPT/OXY1 **
Oxygen Bottle Holder (Size D)



Sun-MPT/DRIP1 ***
Adjustable Height Drip Pole



Sun-MPTS1 *
Folding Shelf



Sun-SBOX1 **
Sharps Box Bracket (Daniels - All Sizes)



Sun-MPT/AGH1 ***
Alcohol Gel Holder



Sun-GBH1 ***
Glove Box Holder



Sun-MPT/Rail1 *
Medi-rail



Sun-MPT/BIN *
2.5 Litre Bin (Requires Medi-rail)



Sun-CVT60/LD ***
2 Tamper Resistant Locking Doors



Sun-DIV3
Single Depth Tray Dividers (per set)



Sun-DIV2
Double Depth Tray Dividers (per set)



Sun-MPT60/CV
Trolley Cover

** Includes Medi-rail *** Locking door option includes a secure top and back panel

dimensions



101cm(H) 70cm(W) 52cm(D)
78cm(W) 60cm(D) with buffers



tray colours

please specify when ordering
(Tray labels and covers included)

Single Depth Trays
8cm(H) 31cm(W) 42cm(D)
Double Depth Trays
16cm(H) 31cm(W) 42cm(D)

Vista 60 Clinical Trolleys

Designed to optimise storage where space is limited, Vista 60 Clinical provide additional working space with compact storage using up to 12 trays, removable recessed stainless steel tray and is functional, durable and easy to clean.



Code: Sun-CVT60
Vista 60 Clinical trolley - 12 single depth trays



Code: Sun-CVT61
Vista 60 Clinical trolley
6 single/3 double depth trays



Code: Sun-CVT62
Vista 60 Clinical trolley - 6 double depth trays



Code: Sun-CVT60/LD
2 Locking doors for Vista 60 clinical trolley

product information



Tilt & locking tray runner system



Upto eight trays



Trays integrate with all other Vista products



Single or double depth tray combinations



Adjustable feet for levelling



Optional tamper resistant locking door



Removable top shelf with integrated upstand



Modules can be attached together to form larger systems



Optional castor upgrades available



Durable epoxy coated white frame



Easy clean design and materials



British made



Supplied fully assembled

accessories * These accessories can only be fitted when ordered with the storage module



Sun-SMLC2 *
4 Castors
(front 2 braked)



Sun-SMLC3A
4 Antistatic Castors
(front 2 braked)



* Sun-BUF3 *
4 Protective Buffers



Sun-SMS1/LD *
Tamper Resistant
Locking Door



Sun-DIV1
Single Depth
Tray Dividers (per set)



Sun-DIV4
Double Depth
Tray Dividers (per set)

*** Locking door option includes a secure top and back panel

dimensions



Single Depth Trays
8cm(H) 47cm(W) 42cm(D)
Double Depth Trays
16cm(H) 47cm(W) 42cm(D)



tray colours
please specify when
ordering
(Tray labels and covers
included)

Vista Low Level Storage Modules

Modular, flexible and colourful, Vista Low Level Storage Modules provide a general storage solution for healthcare environments, using premium materials which are functional, durable and easy to clean.



Code: Sun-SM1/8
Vista low level module
8 single depth trays



Code: Sun-SM1/4/2D
Vista low level module
4 single/2 double depth trays



Code: Sun-SM1/4D
Vista low level module
4 double depth trays



Code: Sun-SM1/LD
Locking door for Vista
low level module



Modules can be attached
together to form larger systems

product information



Tilt & locking tray runner system



Up to sixteen trays



Trays integrate with all other Vista products



Single or double depth tray combinations



Adjustable feet for levelling



Optional tamper resistant locking doors



Removable top shelf with integrated upstand



Modules can be attached together to form larger systems



Optional castor upgrades available



Durable epoxy coated white frame



Easy clean design and materials



British made



Supplied fully assembled

accessories * These accessories can only be fitted when ordered with the storage module



Sun-SMLC2 *
4 Castors
(front 2 braked)



Sun-SMLC3A
4 Antistatic Castors
(front 2 braked)



* Sun-BUF3 *
4 Protective Buffers



Sun-SMS2/LD *
2 Tamper Resistant
Locking Doors




Sun-DIV1
Single Depth
Tray Dividers (per set)




Sun-DIV4
Double Depth
Tray Dividers (per set)

*** Locking door option includes a secure top and back panel

dimensions

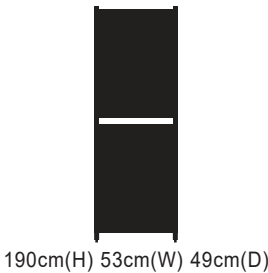

Single Depth Trays
8cm(H) 47cm(W) 42cm(D)
Double Depth Trays
16cm(H) 47cm(W) 42cm(D)

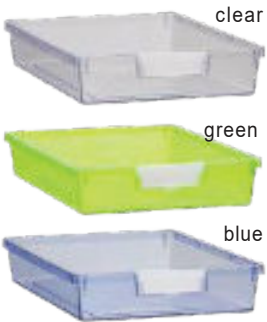

tray colours
please specify when ordering
(Tray labels and covers included)

Vista High Level Storage Modules

Modular, flexible and colourful, Vista High Level Storage Modules provide optimum storage for general healthcare environments, using premium materials which are functional, durable and easy to clean.



Code: Sun-SM2/16
Vista high level module
16 single depth trays



Code: Sun-SM2/8/4D
Vista high level module
8 single/4 double depth trays



Code: Sun-SM2/8D
Vista high level module
8 double depth trays



Code: Sun-SM2/LD
2 locking doors for Vista
high level modules




Modules can be attached
together to form larger systems

product information



Stainless steel trays
(304 grade/DP2 finish)

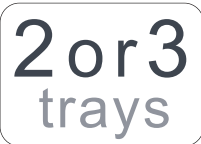

2 or 3 stainless steel
removable trays



Lift-off top tray with
integral handles



Hygienically designed
raised upstands



Trays 'locate' around
frame to prevent
movement



Tray protection
buffer ring



Durable non-marking
75mm castors



Polished stainless
steel frame
(304 grade)



Maximum individual
tray load



Easy clean design
and materials



British made



Supplied fully
assembled

accessories * These accessories can only be fitted when ordered with the trolley



Sun-BUF1*
4 Protective buffers




Sun-LC1*
2 front
braked castors



Sun-LC1A*
4 Antistatic castors
(front 2 braked)

dimensions

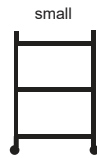

90cm(H) 46cm(W) 52cm(D)



90cm(H) 66cm(W) 52cm(D)

Surgical Trolleys - Stainless Steel Trays

Innovative in design, and in use of materials, our range of surgical trolleys provide hygienic mobile support for clinical environments, whilst being functional, durable and easy to clean.

small - 46cm wide surgical trolleys



Code: Sun-ST4S2
2 stainless steel trays



Code: Sun-ST4S3
3 stainless steel trays

medium - 66cm wide surgical trolleys



Code: Sun-ST6S2
2 stainless steel trays



Code: Sun-ST6S3
3 stainless steel trays

product information



Removable stainless steel shelf (304 grade/DP2 finish))



Lift-off stainless steel trays (304 grade/DP2 finish)



2 or 3 stainless steel shelf/tray combinations



Shelves & trays 'locate' around frame to prevent movement



Shelf/tray protection buffer ring



Durable non-marking 75mm castors



Polished stainless steel frame (304 grade)



Maximum individual tray load



Easy clean design and materials



Supplied fully assembled



British made

accessories * These accessories can only be fitted when ordered with the trolley



Sun-BUF1* 4 Protective buffers



Sun-LC1* 2 front braked castors




Sun-LC1A* 4 Antistatic castors (front 2 braked)

dimensions



86cm(H) 46cm(W) 52cm(D)

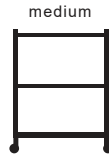

86cm(H) 66cm(W) 52cm(D)

Surgical Trolleys - Stainless Steel Flat Top Shelf & Trays

Innovative in design, and in use of materials, our range of surgical trolleys provide hygienic mobile support for clinical environments, whilst being functional, durable and easy to clean.

small - 46cm wide surgical trolleys



Code: Sun-ST4S2SF
1 stainless steel shelf and 1 tray



Code: Sun-ST4S3SF
1 stainless steel shelf and 2 trays

medium - 66cm wide surgical trolleys



Code: Sun-ST6S2SF
1 stainless steel shelf and 1 tray



Code: Sun-ST6S3SF
1 stainless steel shelf and 2 trays

Surgical Trolleys (Mirror Polished)

product information

premium
quality

Stainless steel shelves,
Trays and frame
(304 grade)

Mirror polished
stainless steel shelves,
Trays and frame

Fully welded
seamless frame

BS 4068 : 1977
Compliant *

2 shelf
options

Flat, fixed or
removable/reversible
shelves

4 Protective buffers

Removable, reversible
folded shelves
(Sun-STFW/RRFS2 only)

Fully welded
fixed shelves
(Sun-STFW/FFS2 only)

Anti-static shelf/tray
protection buffers

75mm antistatic
castors

Supplied fully
assembled

10Kg

Maximum individual
tray load


Easy clean design
and materials

* Designed to comply with the requirements of the hospital market (BS4068: 1977)

accessories * These accessories can only be fitted when ordered with the trolley

Sun-LC1A*
4 Antistatic castors
(front 2 braked)

dimensions

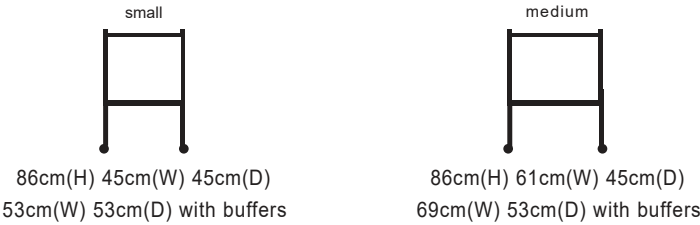

Surgical Trolleys - Stainless Steel (Mirror Polished)

Premium use of materials, and finished to the highest standard, our range of ‘mirror polished’ surgical trolleys provide hygienic mobile support for clinical environments and theatres, whilst being functional, durable and easy to clean.

small - 45cm wide
surgical trolleys (mirror polished)



Code: Sun-STFW4/RRFS2
2 removable reversible folded
shelves (flange up or down)



Code: Sun-STFW4/FFS2
2 fully welded fixed shelves
(flange down)



product information

premium quality

Stainless steel shelves, trays and frame (304 grade)



Mirror polished stainless steel shelves, trays and frame



Fully welded seamless frame



BS 4068 : 1977 Compliant *

2 shelf options

Flat, fixed or removable/reversible shelves



Removable, reversible folded shelves (Sun-STFW/RRFS2 only)



Fully welded fixed shelves (Sun-STFW/FFS2 only)



Anti-static shelf/tray protection buffers



4 Protective buffers



75mm antistatic castors



Supplied fully assembled

10Kg

Maximum individual tray load



Easy clean design and materials

* Designed to comply with the requirements of the hospital market (BS4068: 1977)

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC1A*
4 Antistatic castors (front 2 braked)

Surgical Trolleys - Stainless Steel (Mirror Polished)

Premium use of materials, and finished to the highest standard, our range of ‘mirror polished’ surgical trolleys provide hygienic mobile support for clinical environments and theatres, whilst being functional, durable and easy to clean.

large - 75cm wide
surgical trolleys (mirror polished)

Code: Sun-STFW7/RRFS2
2 removable reversible folded shelves (flange up or down)

extra large - 91cm wide
surgical trolleys (mirror polished)

Code: Sun-STFW9/RRFS2
2 removable reversible folded shelves (flange up or down)

Code: Sun-STFW7/FFS2
2 fully welded fixed shelves (flange down)

Code: Sun-STFW9/FFS2
2 fully welded fixed shelves (flange down)

dimension

large



86cm(H) 75cm(W) 45cm(D)
83cm(W) 53cm(D) with buffers

extra large



86cm(H) 91cm(W) 45cm(D)
99cm(W) 53cm(D) with buffers

51

product information



Unique glass effect safety trays



2 or 3 glass effect removable trays



Lift-off top tray with integral handles



Hygienically designed raised upstands



Trays 'locate' around frame to prevent movement



Tray protection buffer ring



Durable non-marking 75mm castors



Durable epoxy coated white frame



Optional drawers available in a choice of 4 finishes



10Kg
Maximum individual tray load



Easy clean design and materials



British made



Supplied fully assembled

accessories * These accessories can only be fitted when ordered with the trolley



Sun-DT4D*
Single Drawer Unit
(small trolleys)



Sun-DT6D*
Single Drawer Unit
(medium trolleys)



Sun-BUF1*
4 Protective buffers



Sun-LC1*
2 front braked castors



Sun-LC1A*
4 Antistatic castors
(front 2 braked)

dimensions

small



90cm(H) 46cm(W) 52cm(D)

medium



90cm(H) 66cm(W) 52cm(D)

drawer finishes - please specify when ordering



White (high gloss)



Grey

Dressing Trolleys - Glass Effect Safety Trays

Innovative in design, and in use of materials, our range of dressing trolleys provide hygienic mobile support and storage, whilst being functional, durable and easy to clean.

small - 46cm wide dressing trolleys



Code: Sun-DT4G2
2 glass effect safety trays



Code: Sun-DT4G3
3 glass effect safety trays



Sun-DT4G3 with optional
Sun-DT4D single drawer unit

medium - 66cm wide dressing trolleys



Code: Sun-DT6G2
2 glass effect safety trays



Code: Sun-DT6G3
3 glass effect safety trays



Sun-DT6G3 with optional
Sun-DT6D single drawer unit

product information



Stainless steel trays
(304 grade/DP2 finish)



2 or 3
trays



Lift-off top tray with
integral handles



Hygienically designed
raised upstands



Trays 'locate' around
frame to prevent
movement



Tray protection
buffer ring



Durable non-marking
75mm castors



Durable epoxy
coated white frame



Optional drawers
available



10Kg

Maximum individual
tray load



eze

Easy clean design
and materials



British Made
Quality Assured

British made



Supplied Assembled
Supplied Assembled

Supplied fully
assembled

accessories * These accessories can only be fitted when ordered with the trolley



Sun-DT4D*
Single Drawer Unit
(small trolleys)



Sun-DT6D*
Single Drawer Unit
(medium trolleys)



Sun-BUF1*
4 Protective buffers



Sun-LC1*
2 front
braked castors



Sun-LC1A*
4 Antistatic castors
(front 2 braked)

dimensions



drawer finishes - please specify when ordering



White (high gloss)



Grey

Dressing Trolleys - Stainless Steel Trays

Innovative in design, and in use of materials, our range of dressing trolleys provide hygienic mobile support and storage, whilst being functional, durable and easy to clean.

small - 46cm wide dressing trolleys



Code: Sun-DT4S2
2 stainless steel trays



Code: Sun-DT4S3
3 stainless steel trays



Sun-DT4S3 with optional
Sun-DT4D single drawer unit

medium - 66cm wide dressing trolleys



Code: Sun-DT6S2
2 stainless steel trays



Code: Sun-DT6S3
3 stainless steel trays



Sun-DT6S3 with optional
Sun-DT6D single drawer unit

product information

Fully adjustable and configurable laptop/tablet mount allowing comfortable working

Laptop or tablet adjustable securing brackets

Solid and stable laptop plate for error free typing

Dual securing/leashing positions

Optional moulded "Leash/Tether Here" wording on dispensing doors (asstandard) helping agency staff to secure the trolley correctly

Lock supplied with 2 keys, top and bottom boxshare same keys(keyed to differ)

Suited key solutionsincluding pharmacy "master"

Fold down work surface with Hi-Vissafety panel, for protecting hands/fingersfrom being trapped

Hinged lid with slam lock and protective lock escutcheon

Optional integrated lockable lower storage/charging cabinet

Optional integrated recessed storage tray

Pull out storage traysfor waste and consumables

Lift up side table

Gasassisted dispensing door (opening and retention)

Adjustable shelf divider system for drug storage

Infection control friendly

Choice of 9 KYDEX® anti-bacterial door colours

KYDEX® antibacterial moulded construction helping infection control

Free running and easy to manoeuvre with multiple push/pull handles

Robust and secure design with a durable and hygienic silver powder coated frame

Three braked castorsand one anti-static castor with protective buffers

Integrated mouse cradle (mouse not included)

Supplied fully assembled

FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe

*** Keyed to differ as specified by NHSEngland - (keyed to differ means each trolley has its own unique set of keys)

accessories * these accessories can only be fitted when ordered with the trolley ** includes medi-rail

Sun-EPMA/KD/KEY * Additional Keys (Suited)*

Sun-EPMA/STL * Suited Locks

Sun-EPMA/LC3A Four 100mm Antistatic Braked Castors

Sun-EPMA/SBOX1 */** Sharps Box Bracket (Daniels - All Sizes)

Sun-EPMA/SBOX3 */** Sharps Box Bracket (Frontier) 2 & 3 Litre

Sun-EPMA/AGH1 */** Alcohol Gel Holder

Sun-EPMA/GBH1 */** Glove Box Holder

Sun-EPMA/DIV1 Additional Divider Kit (5 Shelf Dividers 100 Base Dividers/10 Labels /5 Clips)

EPMA Ward Drug and Medicine Dispensing Trolleys

The trolleys are designed to provide a complete solution for storing, distributing & dispensing medicine & drugs. Aesthetically pleasing, hygienically designed and combining smart adjustable internal storage with easy access, all trolleys are functional, strong, secure & easy to clean.

Our EPMA Ward Drug Trolleys are available in three sizes, small, medium and large, each size has four options to choose from giving you greater flexibility for your specific needs.

Code: Sun-EPMA1/LAPTOP/BOX/WHITE (Shown)
SMALL EPMA Ward Drug Trolley with Laptop Arm and Bracket, Tray and 1 Secure Bottom Storage Cabinet - White Dispensing Cupboard Doors

Code: Sun-EPMA5/LAPTOP/BOX/WHITE (Shown)
MEDIUM EPMA Ward Drug Trolley with Laptop Arm and Bracket, Trays and 1 Secure Bottom Storage Cabinet - White Dispensing Cupboard Doors

Code: Sun-EPMA11/LAPTOP/TRAYS/WHITE (Shown)
LARGE EPMA Ward Drug Trolley with Laptop Arm and Bracket, Trays and Bottom Shelf - White Dispensing Cupboard Doors

product information



Laminate faced MDF with 2mm edging



Adjustable shelf divider system provides up to 5 Compartments per shelf




Adjustable cabinet base divider system provides up to 16 compartments



Hinged lid with slam lock & protective lock escutcheon



Lock supplied with 2 keys (Keyed to differ)***

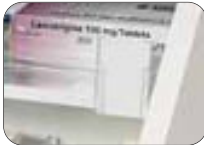

Shelf labels supplied



Gas assisted lid opening & retention



Hygienic stainless steel handle and fittings



Fold down work surface



Lift up side table



4 removable tray positions (2 trays supplied)



Patient notes/file retention strap



Push pull handle (to 3 sides)



Four 100mm braked castors with four protective buffers



Security Leash & Wall Fixing Bracket



Easy clean design and materials



FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe



Supplied fully assembled



British made

*** Keyed to differ as specified by NHS England - (keyed to differ means each trolley has its own unique set of keys)

accessories

* these accessories can only be fitted when ordered with the trolley ** includes medi-rail



Sun-WDT/T2B Additional Double Depth Wide Tray (blue)



Sun-WDT/4Y Additional Double Depth Narrow Tray (yellow)



Sun-WDT/KD/KEY* Additional Keys



Sun-WDT/SBOX1 ** Sharps Box Bracket (Daniels - All Sizes)



Sun-LC3A * Antistatic Braked Castors



Sun-WDT/DIV1 Additional Divider Kit (5 Shelf Dividers, 10 Base Dividers, 10 Labels and 5 Clips)



Sun-WDT/AGH1 *** Alcohol Gel Holder



Sun-WDT/GBH1 *** Glove Box Holder

Ward Drug & Medicine Dispensing Trolleys

Designed to provide a complete solution for storing, distributing and dispensing medicine and drugs. Ward drug & medicine dispensing trolleys are aesthetically pleasing, hygienically designed and combine smart adjustable internal storage with easy access, while being functional, strong, secure and easy to clean.



Code: Sun-WDT20/KD
Ward drug & medicine dispensing trolley (medium capacity) with divider system and 2 storage trays



Sun-WDT20/KD



Wide/Double Depth Trays
16cm(H) 47cm(W) 42cm(D)
Narrow/Double Depth Trays
16cm(H) 31cm(W) 42cm(D)
NB: Maximum tray load 7.5kg
110cm(H) 90cm(W) 58cm(D)
94cm(W) 60cm(D) with buffers



Adjustable storage divider system on shelves and cabinet base



Keyed to Differ as required by NHS England



Code: Sun-WDT50/KD
Ward drug & medicine dispensing trolley (large capacity) with divider system and 2 storage trays



Sun-WDT50/KD



Wide/Double Depth Trays
16cm(H) 47cm(W) 42cm(D)
Narrow/Double Depth Trays
16cm(H) 31cm(W) 42cm(D)
NB: Maximum tray load 7.5kg
128cm(H) 90cm(W) 58cm(D)
94cm(W) 60cm(D) with buffers



Unit Dosage System (UDS) Trolleys

product information - (DT8/UDS Trolley)



9 wide trays and 9 tray dividers



18 narrow trays and 18 tray dividers



Trays provides 54 compartments



6 adjustable height door shelves



Pull-out clear A4 storage tray

product information - (DT7/UDS Trolley)



6 wide trays and 6 tray dividers



12 narrow trays and 12 tray dividers




Trays provides 36 compartments




3 adjustable height door shelves


additional product information



6 pull-out shelves with finger recess



Fold down work surface




Slam lock with 2 keys & lock escutcheon



Integrated cup, pill cup, water jug and pen recess



Integral shelf for patient records file



Push and pull handles, bin & medi-rail




All-round protective rail & front braked castors



Easy clean design and materials



FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe



British made

accessories * these accessories can only be fitted when ordered with the trolley



Sun-DT/KEY* Additional Door Key



Sun-DT/SL1 Security Leash & Wall Fixing Bracket



Sun-DT/B25* Additional 2.5 Litre Bin and Bracket



Sun-DT/AGH1* Alcohol Gel Holder



Sun-DT/WT Additional Wide Tray



Sun-DT/WTD Additional Wide Tray Divider (each)



Sun-DT/WTL Additional Wide Tray Labels (10 pack)



Sun-DT/NT Additional Narrow Tray



Sun-DT/NTD Additional Narrow Tray Divider (each)



Sun-DT/NTL Additional Narrow Tray Labels (10 pack)

Unit Dosage System (UDS) Trolleys

Stylish & contemporary in design, strong, secure and durable in construction, unit dosage drug trolleys are the ideal solution for storing, distributing and dispensing medicines and patient prescriptions. Hygienically designed and finished in attractive beech effect laminate which is easy to clean.



Code: Sun-DT8/UDS
Unit dosage system (UDS) trolley - large



Sun-DT7/UDS
112cm(H)
62cm(W)
60cm(D)



Code: Sun-DT7/UDS
Unit dosage system (UDS) trolley - small



Sun-DT8/UDS
112cm(H)
82cm(W)
60cm(D)



product information



Pull-out shelves with finger recess



Location channel for holding racks



Fold down work surface



3 adjustable height door shelves (DT1/MDS4 DT1/MDS6)



6 adjustable height door shelves (DT2/MDS9)



Pull-out clear A4 storage tray (DT2/MDS9)



Slam lock with 2 keys & lock escutcheon



Integrated cup, pill cup & water jug recesses



Integrated pen recess



Push and pull handles, bin & medi-rail




All-round protective rail



Front push and pull handle



100mm Castors with front brakes

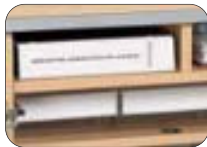

Integral shelf for patients record file



Easy clean design and materials



FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe



Supplied fully assembled



British made

accessories * these accessories can only be fitted when ordered with the trolley



Sun-DT/KEY* Additional Door Key



Sun-DT/SL1 Security Leash & Wall Fixing Bracket



Sun-DT/B25* Additional 2.5 Litre Bin and Bracket



Sun-DT/AGH1* Alcohol Gel Holder



MDS Self Administration Cabinets (see pages 122-123)

• Cards, racks and blister packs are supplied by the pharmacist

dimensions

Sun-DT1/MDS4




112cm(H)
62cm(W)
60cm(D)

Sun-DT1/MDS6



112cm(H)
62cm(W)
60cm(D)

Sun-DT2/MDS9

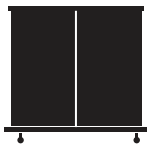

112cm(H)
82cm(W)
60cm(D)

Monitored Dosage System (MDS) Trolleys

Stylish & contemporary in design, strong, secure and durable in construction, monitored dosage drug trolleys are the ideal solution for storing, distributing and dispensing medicines and patient prescriptions stored in blister packs, as supplied by Boots & Venalink. Hygienically designed and finished in attractive beech effect laminate which is easy to clean.



Code: Sun-DT2/MDS9
Monitored dosage system (MDS) trolley - large 9 racks



Code: Sun-DT1/MDS6
Monitored dosage system (MDS) trolley - small 6 racks



Code: Sun-DT1/MDS4
Monitored dosage system (MDS) trolley - small 4 racks



product information



Durable laminate faced MDF with 2mm edging



Secure locking top 2 keys per lock (keyed to differ)***



Digital locking system



Secure solid sides and base



Accepts foolscap or X-Ray size filing pockets (not included)**



Accepts 470mm & 360mm filing pockets (not included)**



Filing pocket identification label for easy re-ordering



Wipe clean filing pocket tabs (supplied with pockets)



Hinged top(s) open to create additional working space



Hygienic stainless steel handles



2 or 4 removable storage trays



Side storage for patients' record file



Push pull handles



Four 100mm braked castors



Four protective buffers



Easy clean design and materials



FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe



Supplied fully assembled



British made

*** Keyed to differ as specified by NHS England - (keyed to differ means each trolley has its own unique set of keys)

accessories * these accessories can only be fitted when ordered with the trolley



Sun-LC3A * Antistatic Braked Castors



Sun-MNT/DCL/KEY* Overriding Key for Digital Locking System



Sun-MNT/FPF/PK10 Filing Pocket Foolscap 360mm(W) x 270mm(D) (Pack of 10)



Sun-MNT/WFPF/PK10 Wide Filing Pocket Foolscap 450mm(W) x 270mm(D) (Pack of 10)



Sun-MNT/FPX/PK10 Filing Pocket X-Ray 360mm(W) x 470mm(D) (Pack of 10)



Sun-MNT/WFPX/PK10 Wide Filing Pocket X-Ray 450mm(W) x 470mm(D) (Pack of 10)



Sun-MNT/FPXF/PK10 Twin Filing Pocket X-Ray & Foolscap 360mm(W) x 270/470mm(D) (Pack of 10)



Sun-MNT/WFPXF/PK10 Twin Wide Filing Pocket X-Ray & Foolscap 450mm(W) x 270/470mm(D) (Pack of 10)

** Sun-MNT10 & Sun-MNT10/DCL Trolley can accept all filing pocket sizes but cannot be combined together.

** Sun-MNT20 & Sun-MNT20/DCL Trolley can accept all filing pocket sizes which can be used together but only in opposite sides of the trolley.

Medical Notes Trolleys

Designed to provide a complete solution for storing, distributing and accessing medical notes and X-Rays. Medical notes trolleys are aesthetically pleasing, hygienically designed and combine smart adjustable filing pockets with easy access, while being functional, strong, secure and easy to clean.



Code: Sun-MNT10
Medical Notes Trolley, enclosed sides with hinged locking top (Shown)

Code: Sun-MNT10/DCL
Medical Notes Trolley, enclosed sides with hinged top and 1 Digital Combination Lock



Optional Digital locking system with optional overriding key (override key(s) must be ordered with trolley)

Sun-MNT10 & 10/DCL



Wide/Single Depth Trays 8cm(H) 47cm(W) 42cm(D)

Wide/Double Depth Trays 16cm(H) 47cm(W) 42cm(D)

NB: Maximum tray load 7.5kg

106cm(H) 64cm(W) 52cm(D)
66cm(W) 57cm(D) with buffers



Keyed to Differ as required by NHS England



Code: Sun-MNT20
Medical Notes Trolley, enclosed sides with hinged locking tops (Shown)

Code: Sun-MNT20/DCL
Medical Notes Trolley, enclosed sides with hinged top and 2 Digital Combination Locks

Sun-MNT20 & 20/DCL



Wide/Single Depth Trays 8cm(H) 47cm(W) 42cm(D)

Wide/Double Depth Trays 16cm(H) 47cm(W) 42cm(D)

NB: Maximum tray load 7.5kg

106cm(H) 121cm(W) 52cm(D)
121cm(W) 57cm(D) with buffers

product information



Durable laminate faced MDF with 2mm edging



Secure locking top (Sun-MRTTL1 only)



Lock is provided with 2 keys (keyed to differ)***



Locking top folds neatly back



Hygienic stainless steel handle



Maximum load 60kg



Strong deep cabinet for optimum file storage



Easy access to files and documents



2 removable storage trays



Secure solid sides and base



Push pull handles



Durable powder coated blue frame



100mm castors (2 braked)



Four protective buffers



Easy clean design and materials




FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe




British made

*** Keyed to differ as specified by NHS England - (keyed to differ means each trolley has its own unique set of keys)

accessories * these accessories can only be fitted when ordered with the trolley




Sun-LC3A * Antistatic Braked Castors




Sun-MRTT/KEY* Additional Key

dimensions



Sun-MRTT1



Sun-MRTTL1

Narrow/Single Depth Trays
8cm(H) 31cm(W) 42cm(D)
NB: Maximum tray load 7.5kg

Storage Box (internal)
25cm(H) 62cm(W) 40cm(D)

98cm(H) 81cm(W) 52cm(D)
83cm(W) 57cm(D) with buffers

Records Transfer Trolleys

Designed to provide a complete solution for storing, distributing and transferring medical notes and files quickly to different locations. Records transfer trolleys are aesthetically pleasing, hygienically designed and easy to access, while being functional, strong, secure and easy to clean.



Code: Sun-MRTT1
Records Transfer Trolley - open



Code: Sun-MRTTL1
Records Transfer Trolley - with folding locking top

*simple
hygienic
disposal*

Clinical Bins & Sack Holders

Our hands-free, flame retardant bins comply with HTM07-01 and HTM83, available in 20L, 50L, and 70L with various lid colours. They are silent closing, rust resistant, easy to clean, and suitable for hospitals, surgeries, care homes, and healthcare environments.



Clinical Bins
Pg 61



Sack Holders
Pg 62

product information

3
sizes

20, 50 & 70
litre capacity

hands
free


Hygienic hands free
operation, foot pedal
must be used to open

HTM07-01
compliant


Safe management
of healthcare
waste compliant

HTM05-03
tested


Designed & tested to
HTM05-03 Part A




Silent, dampened
closing lid as standard




Corrosion resistant
epoxy coated white
galvanised bodies



Silicone bag
retaining cord with
nylon finger loop




Stainless steel foot pedal
with pedal label
(where applicable)




Sealed base
Prevents leakages

wall
mounting
bracket


Integrated wall
mounting bracket




Off the wall design to
help prevent damage



Flame retardant **



Easy clean design
and materials



Delivery normally within
3-5 working days *

** Designed and tested to HTM05-03 Part A * Please contact us regarding availability for quantities of 10 or more.

designated usage - please specify lid colour when ordering

WHITE



General Use

YELLOW



Waste for Incineration

BLACK



Domestic (Municipal) Waste

YELLOW/BLK



Offensive/ Hygiene Waste

ORANGE



Waste Which May Be Treated

PURPLE



Cytotoxic & Cytostatic Waste

BLUE



Medical Waste for Incineration


RED



Anatomical Waste for Incineration


dimensions

20 Litre




51cm(H)
32cm(W)
33cm(D)

50 Litre



59cm(H)
42cm(W)
43cm(D)

70 Litre



78cm(H)
42cm(W)
43cm(D)

Designed for the multi function disposal of waste (HTM07-01), with usage designated by an appropriate coloured lid and pedal label.



Code: Sun-BIN20 - 20 Litre bin



Code: Sun-BIN50 - 50 Litre bin



Code: Sun-BIN70 - 70 Litre bin

product information - sack holders

3
sizes

20, 70 & 90
litre capacity

Hygienic operation,
foot pedal should be
used to open

Elasticated bag
retaining cord with
nylon finger loop

Corrosion resistant
materials

HTM07-01
compliant

Safe management
of healthcare
waste compliant

Easy clean design
and materials

Supplied fully
assembled

British made

* Please contact us regarding availability for quantities of 10 or more.

Sack Holders


Hands-free operation and sturdy design (HTM07-01) is what makes our sack holders hygienic and simple to use. They are ideal for areas where the risk of fire is minimal. Usage is designated by an appropriate coloured lid.



Code: Sun-BINSH-50
50 Litre sack holder

Code: Sun-BINSH-70
70 Litre sack holder

Code: Sun-BINSH-90
90 Litre sack holder

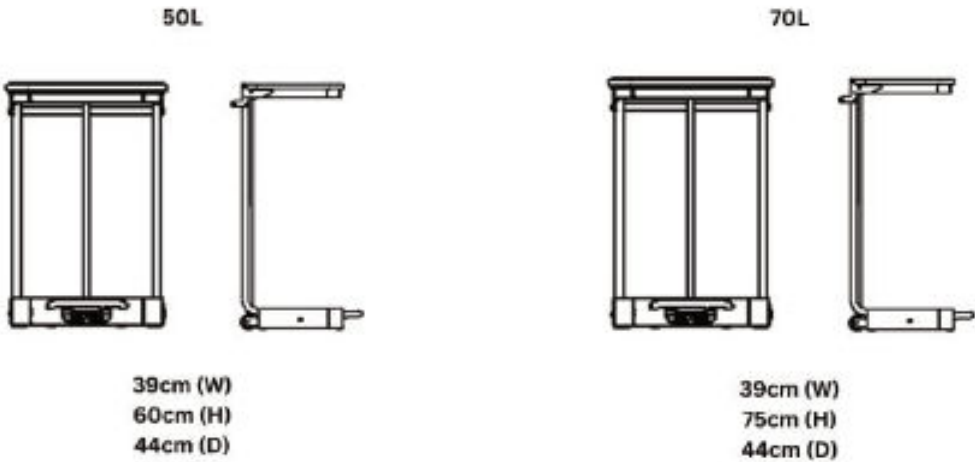

designated usage - please specify lid colour when ordering

<div>WHITE</div> <div></div> <div>General Use</div>	<div>YELLOW</div> <div></div> <div>Waste for Incineration</div>	<div>BLACK</div> <div></div> <div>Domestic (Municipal) Waste</div>	<div>YELLOW/BLK</div> <div></div> <div>Offensive/ Hygiene Waste</div>	<div>ORANGE</div> <div></div> <div>Waste Which May Be Treated</div>
<div>PURPLE</div> <div></div> <div>Cytotoxic & Cytostatic Waste</div>	<div>BLUE</div> <div></div> <div>Medical Waste for Incineration</div>	<div>RED</div> <div></div> <div>Anatomical Waste for Incineration</div>		

*hygienic
mobile
privacy*

Mobile Patient Privacy Screens

Our range of Mobile Folding Screens provide a temporary or permanent privacy solution in hospitals, doctor's surgeries or health care centres. Hygienically designed using our durable and easy to clean panel system, assisting in reducing cross infection.



Mobile Screens - Solid Panels
Pg 64


product information

3/4/5
panel


Available in 3, 4 and 5 panel sizes

8
designs


Unique range of 8 patterns and colours




Hygienic PVC panel system




Flame retardant panels - BS476 part 7




Multiple leg supports for stability




Swivel castors for ease of movement




Durable epoxy coated white frame




Folds flat for storage



Easy clean design and materials


minimal


Minimal assembly required




British made


dimensions

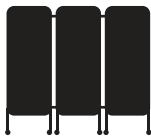

3 Panel
172cm(H) 188cm(W) 3cm(D)
38cm(D) Across Legs

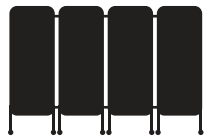

4 Panel
172cm(H) 250cm(W) 3cm(D)
38cm(D) Across Legs

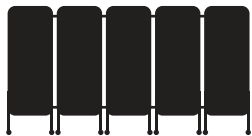

5 Panel
172cm(H) 312cm(W) 3cm(D)
38cm(D) Across Legs


Customisation and branding available!




As a UK manufacturer, we offer many options at little or no extra cost (order size dependent)

Contact us to discuss your requirements.


screen panel patterns and colours - please specify when ordering




white



grey



sky blue



jungle design

Mobile Screens - Solid Panels

Designed to provide a temporary or mobile privacy solution, our range of folding screens incorporate our unique hygienic panel system in 8 patterns and colours and are functional, durable and easy to clean.




3 Panel
172cm(H) 188cm(W) 3cm(D)
38cm(D) Across Legs

4 Panel
172cm(H) 250cm(W) 3cm(D)
38cm(D) Across Legs



5 Panel
172cm(H) 312cm(W) 3cm(D)
38cm(D) Across Legs

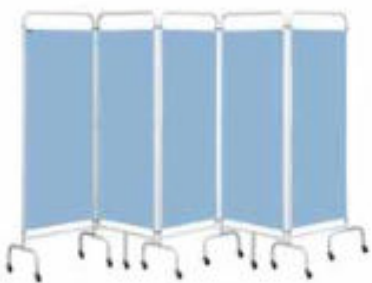

*strong
secure
protection*

Medicine Cabinets

Our British made medicine cabinets provide secure and compliant storage for all pharmaceuticals including controlled drugs. Durable and robust, they feature secure locks and accessible design, ideal for hospitals, clinics, care homes, and healthcare facilities.



Medicine Cabinets
Pg 66-67



Controlled Drug Cabinets
Pg 68-69



Medicine Cabinets with
Controlled Drugs Inner Cabinets
Pg 70-72

product information



Independently tested



BS2881:1989 Security Level 1 compliant



High security 8 lever lock designed to EN1300 and Vds 2396



Anti-pick levers & deadlock mechanism



3 keys supplied (per lock)
Extra keys available



Protective key escutcheon



Adjustable shelves and door trays



Internal door retaining flange



Durable and hygienic white epoxy coating



Single door cabinets are right hand hinged (left hand optional)



Optional warning light and buzzer alarm



Easy clean design and materials



British made



Supplied fully assembled



Delivery normally within 3-5 working days *
(Sun-MC1, MC2 & MC3 only)

* Please contact us regarding availability for quantities of 10 or more.

accessories * Factory fit - only available at time of ordering cabinet



Sun-MC/LIGHT1 *
(For Sun-MC1/MC2/MC3 Cabinets)
Sun-MC/LIGHT2 *
(For Sun-MC6/MC7 Cabinets)



Sun-MC/BUZZ *
Buzzer alarm (has to be combined with Sun-MC/Light)



Sun-MC/LHH *
Left hand hinge option



Sun-MC/KEY *
Additional key (when ordered with cabinet)

5V DC external warning light shows green when closed, red when open and flashes red when left open for 2.5 minutes. Nurse call output, (volt free contact) for connection to nurse call system. Internal illumination when door is opened. Must be hard wired to a 220/240V AC mains supply. Buzzer alarm activates when cabinet is left open for 5 minutes (only available when combined with optional warning light)

dimensions - NB: Optional warning light transformer box adds 55mm to overall height



MC1: 60cm(H), 30cm(W), 30cm(D)
MC2: 60cm(H), 50cm(W), 30cm(D)
MC3: 60cm(H), 60cm(W), 30cm(D)
MC7: 85cm(H), 60cm(W), 30cm(D)

Medicine Cabinets

Independently tested to conform to or exceed all the required standards and legislation, all our medicine cabinets provide secure storage, which can be adapted internally to your exact requirements.



Code: Sun-MC1
Medicine cabinet with 3 shelves, 2 door trays, shown with optional warning light



Code: Sun-MC2
Medicine cabinet with 3 shelves, 2 door trays, shown with optional warning light



Code: Sun-MC3
Medicine cabinet with 3 shelves, 2 door trays, shown with optional warning light



Code: Sun-MC7
Medicine cabinet with 4 shelves, 4 door trays, shown with optional warning light



product information



Independently tested



BS2881:1989 Security Level 1 compliant



High security 8 lever lock designed to EN1300 and Vds 2396



Anti-pick levers & deadlock mechanism



3 keys supplied (per lock) Extra keys available



Protective key escutcheon



Adjustable shelves and door trays



Internal door retaining flange



Durable and hygienic white epoxy coating



Optional warning light and buzzer alarm



Easy clean design and materials



British made



Supplied fully assembled

accessories * Factory fit - only available at time of ordering cabinet



Sun-MC/LIGHT3 *
(For Sun-MC4/MC5/MC8/MC9 Cabinets)



Sun-MC/BUZZ *
Buzzer alarm (has to be combined with Sun-MC/Light)



Sun-MC/KEY *
Additional key (when ordered with cabinet)

5V DC external warning light shows green when closed, red when open and flashes red when left open for 2.5 minutes. Nurse call output, (volt free contact) for connection to nurse call system. Internal illumination when door is opened. Must be hard wired to a 220/240V AC mains supply. Buzzer alarm activates when cabinet is left open for 5 minutes (only available when combined with optional warning light)


dimensions - NB: Optional warning light transformer box adds 55mm to overall height




60cm(H)
80cm(W)
30cm(D)



60cm(H)
100cm(W)
30cm(D)

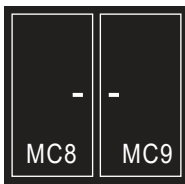

85cm(H)
80cm(W)
30cm(D)

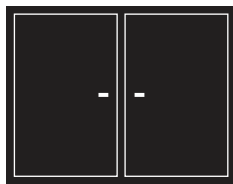

85cm(H)
100cm(W)
30cm(D)

Medicine Cabinets

Independently tested to conform to or exceed all the required standards and legislation, all our medicine cabinets provide secure storage, which can be adapted internally to your exact requirements.



Code: Sun-MC4
Medicine cabinet with 6 shelves, 4 door trays, shown with optional warning light



Code: Sun-MC5
Medicine cabinet with 6 shelves, 4 door trays, shown with optional warning light



Code: Sun-MC8
Medicine cabinet with 8 shelves, 8 door trays, shown with optional warning light



Code: Sun-MC9
Medicine cabinet with 8 shelves, 8 door trays, shown with optional warning light

product information



Independently tested

misuse
of drugs

Meets the Misuse of
Drugs (Safe Custody
Regulation 1973)



BS 2881:1989
Security Level 1
compliant



High security 8 lever
lock designed to
EN1300 and Vds 2396



Anti-pick levers &
deadlock mechanism



3 keys supplied
(per lock)
Extra keys available



Protective key
escutcheon



Adjustable shelves
and door trays



Internal door
retaining flange



Durable and hygienic
white epoxy coating



Single door cabinets
are right hand hinged
(left hand optional)



Optional warning light
and buzzer alarm



Easy clean design
and materials



Supplied fully
assembled




British made



Delivery normally within
3-5 working days *
(Sun-CDC21, 22 & 23 only)

* Please contact us regarding availability for quantities of 10 or more.

accessories * Factory fit - only available at time of ordering cabinet

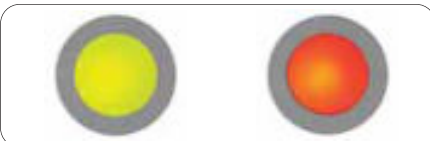

Sun-CDC/LIGHT1 *
(For Sun-CDC21/CDC22/CDC23 Cabinets)
Sun-CDC/LIGHT2 *
(For Sun-CDC24/CDC25 Cabinets)



Sun-CDC/BUZZ *
Buzzer alarm (has to
be combined with
Sun-CDC/Light)



Sun-CDC/LHH *
Left hand hinge option



Sun-CDC/KEY *
Additional key
(when ordered
with cabinet)

5V DC external warning light shows green when closed, red when open and flashes red when left open for 2.5 minutes. Nurse call output, (volt free contact) for connection to nurse call system. Internal illumination when door is opened. Must be hard wired to a 220/240V AC mains supply. Buzzer alarm activates when cabinet is left open for 5 minutes (only available when combined with optional warning light)

dimensions - NB: Optional warning light transformer box adds 55mm to overall height

CDC21	CDC22	CDC23	CDC24	CDC25
30cm(H) 21cm(W) 27cm(D)	30cm(H) 34cm(W) 27cm(D)	55cm(H) 34cm(W) 27cm(D)	55cm(H) 50cm(W) 30cm(D)	60cm(H) 60cm(W) 30cm(D)

Controlled Drug Cabinets

Independently tested to conform to or exceed all the required standards and legislation, all our controlled drug cabinets provide secure storage, which can be adapted internally to your exact requirements.



Code: Sun-CDC21
Controlled drug cabinet with 1 shelf



Code: Sun-CDC22
Controlled drug cabinet with 1 shelf and 1 door tray,
shown with optional warning light



Code: Sun-CDC23
Controlled drug cabinet with
2 shelves and 2 door trays,
shown with optional warning light



Code: Sun-CDC24
Controlled drug cabinet with
2 shelves and 2 door trays,
shown with optional warning light



Code: Sun-CDC25
Controlled drug cabinet with 2 shelves, 2 door trays and 2 locks,
shown with optional warning light

product information



Independently tested



Meets the Misuse of Drugs (Safe Custody Regulation 1973)



BS 2881:1989 Security Level 1 compliant



High security 8 lever lock designed to EN1300 and Vds 2396



Anti-pick levers & deadlock mechanism



3 keys supplied (per lock)
Extra keys available



Protective key escutcheon



Adjustable shelves and door trays



Internal door retaining flange



Durable and hygienic white epoxy coating



Single door cabinets are right hand hinged (left hand optional)



Optional warning light and buzzer alarm



Easy clean design and materials



British made



Supplied fully assembled

accessories * Factory fit - only available at time of ordering cabinet



Sun-CDC/LIGHT2 *
(CDC28 Cabinets)
Sun-CDC/LIGHT3 *
(For Sun-CDC27 Cabinet)



Sun-CDC/BUZZ *
Buzzer alarm (has to be combined with Sun-CDC/Light)




Sun-CDC/LHH *
Left hand hinge option



Sun-CDC/KEY *
Additional key (when ordered with cabinet)

5V DC external warning light shows green when closed, red when open and flashes red when left open for 2.5 minutes. Nurse call output, (volt free contact) for connection to nurse call system. Internal illumination when door is opened. Must be hard wired to a 220/240V AC mains supply. Buzzer alarm activates when cabinet is left open for 5 minutes (only available when combined with optional warning light)

dimensions - NB: Optional warning light transformer box adds 55mm to overall height

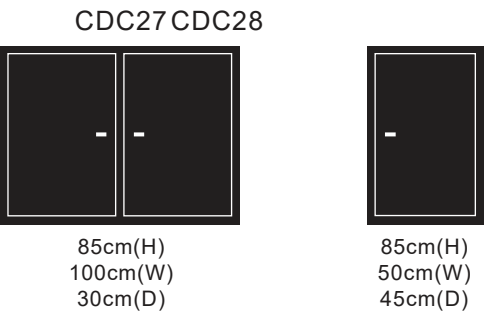

Controlled Drug Cabinets

Independently tested to conform to or exceed all the required standards and legislation, all our controlled drug cabinets provide secure storage, which can be adapted internally to your exact requirements.



Code: Sun-CDC28
Controlled drug cabinet with 4 shelves and 4 door trays, shown with optional warning light



Code: Sun-CDC27
Controlled drug cabinet with 8 shelves and 8 door trays, shown with optional warning light



product information



Independently tested



Meets the Misuse of Drugs (Safe Custody Regulation 1973)



BS 2881:1989 Security Level 1 compliant



High security 8 lever lock designed to EN1300 and Vds 2396



Anti-pick levers & deadlock mechanism



3 keys supplied (per lock)
Extra keys available



Protective key escutcheon on outer cabinet



5V External warning light is green when closed, red when open - fitted as standard on outer medicine cabinet



Optional buzzer alarm



Adjustable shelves and door trays



Identification coloured key fobs



Identification key escutcheon on Inner cabinet



Internal door retaining flange



Durable and hygienic white epoxy coating



Single door cabinets are right hand hinged (left hand optional)



Easy clean design and materials




British made




Supplied fully assembled

accessories * Factory fit - only available at time of ordering cabinet



Sun-MC/BUZZ *
Buzzer alarm



Sun-MC/LHH *
Left hand hinge option




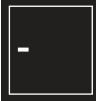




Sun-MC/KEY *
Additional key for outer cabinet (when ordered with cabinet)



Sun-CDC/KEY *
Additional key for controlled drug inner cabinet (when ordered with cabinet)

Warning Light (fitted as standard)
5V DC external warning light shows green when closed, red when open and flashes red when left open for 2.5 minutes. Nurse call output, (volt free contact) for connection to nurse call system. Internal illumination when door is opened. Must be hard wired to a 220/240V AC mains supply. Buzzer alarm activates when cabinet is left open for 5 minutes (only available when combined with optional warning light)

dimensions - NB: Height measurement for outer cabinets includes transformer box

Controlled Drug Inner Cabinets			Medicine Outer Cabinets		
CDC21	CDC22	CDC23	MC2	MC3	MC7
					
30cm(H) 21cm(W) 27cm(D)	30cm(H) 34cm(W) 27cm(D)	55cm(H) 34cm(W) 27cm(D)	66cm(H) 50cm(W) 30cm(D)	66cm(H) 60cm(W) 30cm(D)	91cm(H) 60cm(W) 30cm(D)

Independently tested to conform to or exceed all the required standards and legislation. Our range of combined medicine and drug cabinets provide secure storage, which can be adapted internally to your exact requirements. A warning light and buzzer provide additional security.



CODE: Sun-MCDC221
MC2 Medicine outer cabinet combined with CDC21 controlled drugs inner cabinet



CODE: Sun-MCDC321
MC3 Medicine outer cabinet combined with CDC21 controlled drugs inner cabinet



CODE: Sun-MCDC322
MC3 Medicine outer cabinet combined with CDC22 controlled drugs inner cabinet



CODE: Sun-MCDC722
MC7 Medicine outer cabinet combined with CDC22 controlled drugs inner cabinet



CODE: Sun-MCDC723
MC7 Medicine outer cabinet combined with CDC23 controlled drugs inner cabinet

product information



Independently tested



Meets the Misuse of Drugs (Safe Custody Regulation 1973)



BS 2881:1989 Security Level 1 compliant



High security 8 lever lock designed to EN1300 and Vds 2396



Anti-pick levers & deadlock mechanism



3 keys supplied (per lock)
Extra keys available



Protective key escutcheon on outer cabinet



5V External warning light is green when closed, red when open - fitted as standard on outer medicine cabinet



Optional buzzer alarm



Adjustable shelves and door trays



Identification coloured key fobs



Identification key escutcheon on Inner cabinet



Internal door retaining flange



Durable and hygienic white epoxy coating



Easy clean design and materials




British made




Supplied fully assembled


accessories * Factory fit - only available at time of ordering cabinet



Sun-MC/BUZZ *
Buzzer alarm






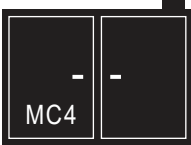

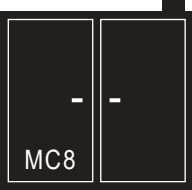
Sun-MC/KEY *
Additional key for outer cabinet
(when ordered with cabinet)



Sun-CDC/KEY *
Additional key for controlled drug inner cabinet
(when ordered with cabinet)

Warning Light (fitted as standard)
5V DC external warning light shows green when closed, red when open and flashes red when left open for 2.5 minutes. Nurse call output, (volt free contact) for connection to nurse call system. Internal illumination when door is opened. Must be hard wired to a 220/240V AC mains supply. Buzzer alarm activates when cabinet is left open for 5 minutes (only available when combined with optional warning light)

dimensions - NB: Height measurement for outer cabinets includes transformer box

Controlled Drug Inner Cabinets			Medicine Outer Cabinets		
CDC21	CDC22	CDC23			
					
30cm(H) 21cm(W) 27cm(D)	30cm(H) 34cm(W) 27cm(D)	55cm(H) 34cm(W) 27cm(D)	66cm(H) 80cm(W) 30cm(D)	66cm(H) 100cm(W) 30cm(D)	91cm(H) 80cm(W) 30cm(D)

Independently tested to conform to or exceed all the required standards and legislation. Our range of combined medicine and drug cabinets provide secure storage, which can be adapted internally to your exact requirements. A warning light and buzzer provide additional security.



CODE: Sun-MCDC421
MC4 Medicine outer cabinet combined with CDC21 controlled drugs inner cabinet



CODE: Sun-MCDC521
MC5 Medicine outer cabinet combined with CDC21 controlled drugs inner cabinet



CODE: Sun-MCDC522
MC5 Medicine outer cabinet combined with CDC22 controlled drugs inner cabinet



CODE: Sun-MCDC821
MC8 Medicine outer cabinet combined with CDC21 controlled drugs inner cabinet



product information



Independently tested



Meets the Misuse of Drugs (Safe Custody Regulation 1973)



BS 2881:1989 Security Level 1 compliant



High security 8 lever lock designed to EN1300 and Vds 2396



Anti-pick levers & deadlock mechanism



3 keys supplied (per lock)
Extra keys available



Protective key escutcheon on outer cabinet



5V External warning light is green when closed, red when open - fitted as standard on outer medicine cabinet



Optional buzzer alarm



Adjustable shelves and door trays



Identification coloured key fobs



Identification key escutcheon on Inner cabinet



Internal door retaining flange



Durable and hygienic white epoxy coating



Easy clean design and materials




British made




Supplied fully assembled


accessories * Factory fit - only available at time of ordering cabinet



Sun-MC/BUZZ *
Buzzer alarm



Sun-MC/KEY *
Additional key for outer cabinet
(when ordered with cabinet)




Sun-CDC/KEY *
Additional key for controlled drug inner cabinet
(when ordered with cabinet)

Warning Light (fitted as standard)
5V DC external warning light shows green when closed, red when open and flashes red when left open for 2.5 minutes. Nurse call output, (volt free contact) for connection to nurse call system. Internal illumination when door is opened. Must be hard wired to a 220/240V AC mains supply. Buzzer alarm activates when cabinet is left open for 5 minutes (only available when combined with optional warning light)

dimensions - NB: Height measurement for outer cabinets includes transformer box


Controlled Drug Inner Cabinets

CDC21




30cm(H)
21cm(W)
27cm(D)

CDC22




30cm(H)
34cm(W)
27cm(D)

CDC23



55cm(H)
34cm(W)
27cm(D)

MC9



91cm(H)
100cm(W)
30cm(D)

Independently tested to conform to or exceed all the required standards and legislation. Our range of combined medicine and drug cabinets provide secure storage, which can be adapted internally to your exact requirements. A warning light and buzzer provide additional security.



CODE: Sun-MCDC921
MC9 Medicine outer cabinet combined with CDC21 controlled drugs inner cabinet



CODE: Sun-MCDC922
MC9 Medicine outer cabinet combined with CDC22 controlled drugs inner cabinet



CODE: Sun-MCDC923
MC9 Medicine outer cabinet combined with CDC23 controlled drugs inner cabinet



coordinating
personal
storage

Ward Furniture

Our mobile bedside cabinets, wardrobes, and overbed tables meet NHS standards, offering durable, easy to clean, and flexible storage solutions with drawers, cupboards, side doors, and attached lockers, ideal for healthcare settings.



Axis Bedside Cabinets
Pg 74



Bedside Combination Units
Pg 75



Wardrobes
Pg 76



Overbed/Overchair Tables
Pg 77-79

product information - (Sun-BSL1)



Large locking top drawer



Locking small middle drawer for extra storage



Easy access cupboard with twin saloon doors



Built in shelf inside cupboard for storage

product information - (Sun-BSL2)



Large locking fold down shelf



Open access shelf for extra storage



Easy access cupboard with twin saloon doors



Built in shelf inside cupboard for storage

product information - (Sun-BSL3)



Large locking top drawer



Large open access shelf for extra storage



Locking bottom cupboard (left or right hand hinged)



Built in shelf inside cupboard for storage

product information - All Models



Wrapped matching rounded vinyl tops for added safety (Axis)




60mm Ø castors (front 2 braked)




"Bolted through" castors & door fixings for strength, hygiene and longevity



Stylish satin finish easy clean handles

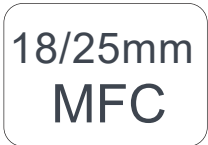

18/25mm MFC

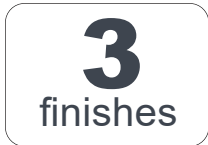

3 finishes



BS EN 16121: 2013 Level 2 Safety, Strength, Durability & Stability



FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe



Supplied fully assembled




Easy clean design and materials




British made


dimensions

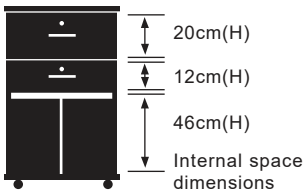

Sun-BLS1
99cm(H) 45cm(W) 49cm(D)

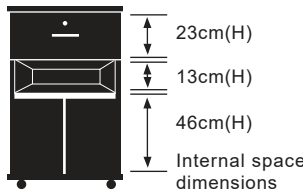

Sun-BLS2
99cm(H) 45cm(W) 49cm(D)

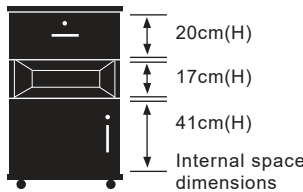

Sun-BLS3
99cm(H) 45cm(W) 49cm(D)

Axis Bedside Cabinets

Stylish and contemporary in design, strong and durable in construction our Axis Hospital Bedside Locker range provide the ideal solution for storing and securing medicines and personal belongings on the ward. Manufactured from robust MFC material and designed to BS1765-1:1990. Our Axis bedside lockers come with rounded vinyl wrapped matching tops. They have been designed to be cost effective, strong, easy to clean and offers good infection control performance.

Axis Bedside Lockers



Code: Sun-BSL1/COLOUR/KA
Axis Bedside Locker, keyed alike, large locking top drawer, locking small middle drawer, saloon bottom doors



Code: Sun-BSL2/COLOUR/KA
Axis Bedside Locker, keyed alike, large locking fold down shelf, open access shelf, saloon bottom doors



Code: Sun-BSL3/COLOUR/KA
Axis Bedside Locker, keyed alike, large locking top drawer, large open access shelf, right or left hand hinged locking bottom cupboard

Lock Options



Keyed Alike Locks (KA)



Keyed to Differ Locks (KD)



Digital Combination Locks (DCL)

Finishes



White (WHT) Oak (OAK)



Beech (BCH)



Right or left hand hinged cupboard doors (Sun-BSL3 only)

Bedside Combination Units (Mobile)

product information



Available with or without locks
(2 keys per lock)

18/25mm
MFC


Manufactured
chip board
construction

18/25mm
LFMDF

Manufactured from
laminate faced MDF to
BS1765-1 : 1990



Conforms to
British Standard
BS1765-1 : 1990
(LFMDF models only)

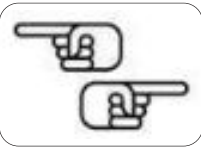

Available in left
or right hand



Durable 2mm panel
edging all round



Chrome hanging
rail in wardrobe



60mm castors
(front 2 braked)



Stylish flush
fitting handles

3
finishes

Maple, beech, oak
and white



FIRA Standard 6250
2005 Specification:
Adhesion Performance
General/Severe



Supplied fully
assembled




Easy clean design
and materials




British made

dimensions

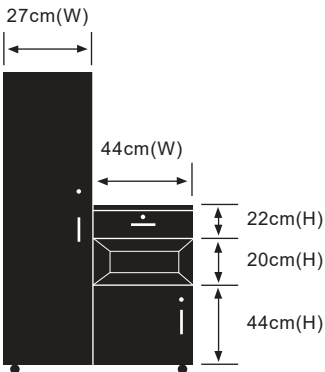

Bedside Combination Unit (Mobile)
144cm(H) 71cm(W) 49cm(D)

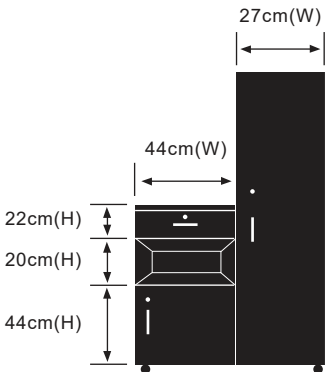

Bedside Combination Unit (Mobile)
144cm(H) 71cm(W) 49cm(D)

finishes



beech oak white

Bedside Combination Units (Mobile)

Stylish & contemporary in design, strong and durable in construction, our bedside combination units provide the ideal solution for storing and securing personal items. Choose from standard MFC material or LFMDf for compliance with BS1765-1:1990. Available in 3 finishes, all of which are designed to be easy to clean.



Left or right hand bedside combination unit with single locking wardrobe & hanging rail, locking top drawer, middle open storage recess & large locking cupboard with built in shelf.

Code: Sun-CBHBC4/MFC/LOCKS

Left hand bedside combination unit with locks.
Manufactured from MFC material.

Code: Sun-CBHBC4/LFMDF/LOCKS

Left hand bedside combination unit with locks.
Manufactured from laminate faced MDFmaterial.



Code: Sun-CBHBC5/MFC/LOCKS

Right hand bedside combination unit with locks.
Manufactured from MFC material.

Code: Sun-CBHBC5/LFMDF/LOCKS

Right hand bedside combination unit with locks.
Manufactured from laminate faced MDFmaterial.



Left or right hand bedside combination unit with single wardrobe & hanging rail, top drawer, middle open storage recess, & large cupboard with built in shelf.

Code: Sun-CBHBC4/MFC

Left hand bedside combination unit.
Manufactured from MFC material.

Code: Sun-CBHBC4/LFMDF

Left hand bedside combination unit.
Manufactured from laminate faced MDFmaterial.



Code: Sun-CBHBC5/MFC

Right hand bedside combination unit.
Manufactured from MFC material.

Code: Sun-CBHBC5/LFMDF

Right hand bedside combination unit.
Manufactured from laminate faced MDFmaterial.



product information

18mm
MFC

Manufactured
chip board
construction



Two shelves
(One above hanging rail
and one at base)



Chrome
hanging rail



Durable 2mm panel
edging all round



Stylish flush
fitting handles

3
finishes

Maple, beech, oak
and white



Easy clean design
and materials



British made



Supplied fully
assembled

dimensions


HEAVY
ITEM



Single
wardrobe
185cm(H)
43cm(W)
60cm(D)




Gents single
wardrobe
185cm(H)
43cm(W)
60cm(D)

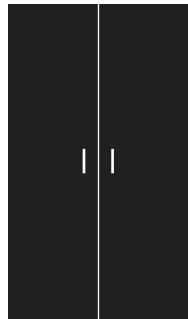

Double
wardrobe
185cm(H)
80cm(W)
60cm(D)

finishes

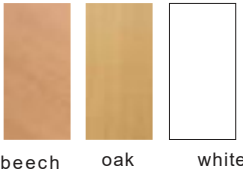

Wardrobes

Stylish & contemporary in design, strong and durable in construction, our range of wardrobes provide the ideal solution for storing personal items. Available in three finishes, all of which are designed to be easy to clean.



Code: Sun-CBSW1/MFC
Single wardrobe with hanging
rail and 2 shelves



Code: Sun-CBSW2/MFC
Gents single wardrobe with hanging
rail, 2 shelves and 2 drawers



Code: Sun-CBDW1/MFC
Double wardrobe with
hanging rail and 2 shelves



product information



MFC is strong, durable, and resistant to stains and heat



Compact grade laminate is extremely strong & durable, resistant to stains, moisture, heat & scratching



Compact grade laminate has a smooth finish, rounded corners & infection control friendly



C-Shaped Base for overbed or chair use




Parallel leg base for overbed or chair use



1 Piece flat top



Quick and easy 35cm height adjustment (between 74cm-109cm)

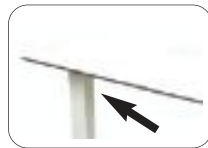

Anti-trapping safety feature preventstrapping of patient between the bed & table



Sealed tubular steel base - infection control friendly



4 High quality swivel castors - for ease of movement




Durable and hygienic white epoxy coating base & frame



Safe working load (Spread evenly)



Product Weight

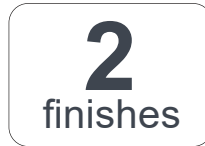

2 finishes
Beech and White



FIRA Standard 6250 2005
Specification:
Adhesion Performance
General/Severe



Easy clean design and materials



British made



Supplied fully assembled

dimensions



109cm
Height Adjustment 35cm
74cm

Overbed/Overchair Table
109cm(H)
95cm(W)
45cm(D)



Table Top
95cm(W) 45cm(D)

finishes



beech white



Code: Sun-OBT2C
Overbed/Overchair Table with
C-Shaped Base & Flat Top



Code: Sun-OBT2P
Overbed/Overchair Table with
Parallel Base & Flat Top

product information



Compact grade laminate is extremely strong & durable, resistant to stains, moisture, heat & scratching



Compact grade laminate has a smooth finish, rounded corners & infection control friendly



C-Shaped Base for overbed or chair use



Parallel leg base for overbed or chair use




1 Piece flat top, with raised lip (Sun-OBT3)



Sealed tubular steel base - infection control friendly




4 High quality swivel castors - for ease of movement

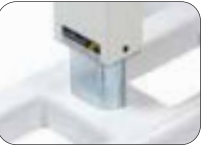

Durable and hygienic white epoxy coating base & frame



Quick and easy 35cm height adjustment (between 74cm-109cm)

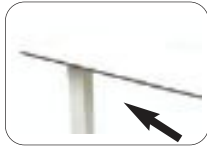

Anti-trapping safety feature prevents trapping of patient between the bed & table




25kg

Safe working load



20kg

Product Weight

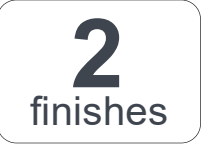

2 finishes

beech, blue and white



Supplied fully assembled



British made




FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe



Easy clean design and materials

dimensions

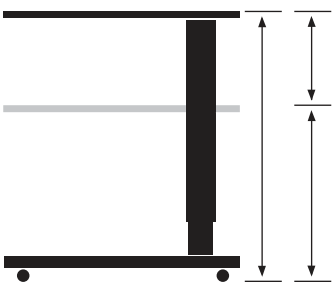


109cm
Height Adjustment 35cm
74cm

Overbed/Overchair Table
109cm(H)
95cm(W)
45cm(D)



Table Top
95cm(W) 45cm(D)

finishes

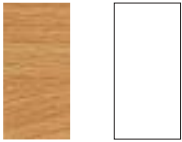

beech

white

Astro Overbed/Overchair Tables
(Compact Grade Laminate Top)



Stylish & contemporary in design, strong and durable in construction, our overbed tables provide the ideal daily solution for multiple uses. Available in 5 specifically designed table tops in a choice of 2 materials, (Compact Grade Laminate and High Impact PVC/Kydex), 2 finishes and parallel or C-shaped base.



Code: Sun-OBT3C
Overbed/Overchair Table with C-Shaped Base & Flat Top with 1 Raised Lip



Code: Sun-OBT3P
Overbed/Overchair Table with Parallel Base & Flat Top with 1 Raised Lip

product information



High impact PVC is strong & durable, resistant to stains & heat



High impact PVC has a smooth finish, rounded corners & infection control friendly



C-Shaped Base for overbed or chair use



Parallel leg base for overbed or chair use




1 Piece recessed top (Sun-OBT5)



Recessed table top helps to retain liquids and other items



Quick and easy 35cm height adjustment (between 74cm-109cm)

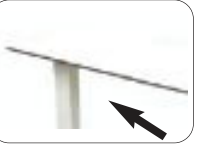

Anti-trapping safety feature preventstrapping of patient between the bed & table



Safe working load




Product Weight



Sealed tubular steel base - infection control friendly



4 High quality swivel castors - for ease of movement

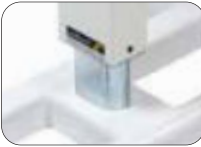


Durable and hygienic white epoxy coating base & frame



Supplied fully assembled



British made

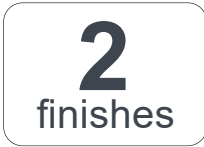

2 finishes
beech, blue and white



FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe



Easy clean design and materials

dimensions



109cm
Height Adjustment 35cm
74cm
Overbed/Overchair Table
109cm(H)
95cm(W)
45cm(D)



Table Top
95cm(W) 45cm(D)

finishes



beech white

Overbed/Overchair Tables
(High Impact PVC Top)

Stylish & contemporary in design, strong and durable in construction, our overbed tables provide the ideal daily solution for multiple uses. Available in 5 specifically designed table tops in a choice of 2 materials, (Compact Grade Laminate and High Impact PVC/Kydex), 3 finishes and parallel or C-shaped base.



Code: Sun-OBT5C
Overbed/Overchair Table with C-Shaped Base & Flat Recessed Top



Code: Sun-OBT5P
Overbed/Overchair Table with Parallel Base & Flat Recessed Top

High Impact PVC is a heavy duty material which can be formed into recesses, providing a sealed area which can retain liquids and other items. It is strong and durable, resistant to stains and heat. Its smooth finish and rounded edges make it easy to clean and is infection control friendly.

hygienic
safe
procedures

Stainless Steel Healthcare Furniture

Simple, practical, easy to manoeuvre, our trolleys and drip stands are suitable for healthcare environments. We provide a range of designs and sizes, available in a choice of materials to suit your budget and hygiene requirements.

Vista Trolleys & Storage Modules



Surgical Trolley
Pg 81-82



Surgical Trolleys
(Mirror Polished)
Pg 83-84

Drip Stands



Mobile Drip Stands
Pg 85

product information



Stainless steel trays
(304 grade/DP2 finish)



2 or 3 stainless steel
removable trays



Lift-off top tray with
integral handles



Hygienically designed
raised upstands



Trays 'locate' around
frame to prevent
movement



Tray protection
buffer ring



Durable non-marking
75mm castors



Polished stainless
steel frame
(304 grade)



Maximum individual
tray load



Easy clean design
and materials



British made



Supplied fully
assembled

accessories * These accessories can only be fitted when ordered with the trolley



Sun-BUF1*
4 Protective buffers




Sun-LC1*
2 front
braked castors




Sun-LC1A*
4 Antistatic castors
(front 2 braked)

dimensions

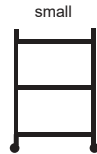

90cm(H) 46cm(W) 52cm(D)

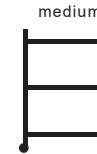

90cm(H) 66cm(W) 52cm(D)

Surgical Trolleys - Stainless Steel Trays

Innovative in design, and in use of materials, our range of surgical trolleys provide hygienic mobile support for clinical environments, whilst being functional, durable and easy to clean.

small - 46cm wide surgical trolleys



Code: Sun-ST4S2
2 stainless steel trays



Code: Sun-ST4S3
3 stainless steel trays

medium - 66cm wide surgical trolleys



Code: Sun-ST6S2
2 stainless steel trays



Code: Sun-ST6S3
3 stainless steel trays

product information



Removable stainless steel shelf (304 grade/DP2 finish))



Lift-off stainless steel trays (304 grade/DP2 finish)



2 or 3 stainless steel shelf/tray combinations



Shelves & trays 'locate' around frame to prevent movement



Shelf/tray protection buffer ring



Durable non-marking 75mm castors



Polished stainless steel frame (304 grade)



Maximum individual tray load



Easy clean design and materials



Supplied fully assembled



British made

accessories * These accessories can only be fitted when ordered with the trolley



Sun-BUF1* 4 Protective buffers




Sun-LC1* 2 front braked castors



Sun-LC1A* 4 Antistatic castors (front 2 braked)

dimensions

small

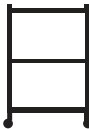

86cm(H) 46cm(W) 52cm(D)

medium



86cm(H) 66cm(W) 52cm(D)

Surgical Trolleys - Stainless Steel Flat Top Shelf & Trays

Innovative in design, and in use of materials, our range of surgical trolleys provide hygienic mobile support for clinical environments, whilst being functional, durable and easy to clean.

small - 46cm wide surgical trolleys



Code: Sun-ST4S2SF
1 stainless steel shelf and 1 tray



Code: Sun-ST4S3SF
1 stainless steel shelf and 2 trays

medium - 66cm wide surgical trolleys



Code: Sun-ST6S2SF
1 stainless steel shelf and 1 tray



Code: Sun-ST6S3SF
1 stainless steel shelf and 2 trays

Surgical Trolleys (Mirror Polished)

product information

premium quality

Stainless steel shelves, trays and frame (304 grade)



Mirror polished stainless steel shelves, trays and frame



Fully welded seamless frame



BS 4068 : 1977 Compliant *

2 shelf options

Flat, fixed or removable/reversible shelves



4 Protective buffers



Removable, reversible folded shelves (Sun-STFW/RRFS2 only)



Fully welded fixed shelves (Sun-STFW/FFS2 only)



Anti-static shelf/tray protection buffers



75mm antistatic castors



Supplied fully assembled

10Kg

Maximum individual tray load



Easy clean design and materials


* Designed to comply with the requirements of the hospital market (BS4068: 1977)

accessories * These accessories can only be fitted when ordered with the trolley



Sun-LC1A*
4 Antistatic castors (front 2 braked)

dimensions

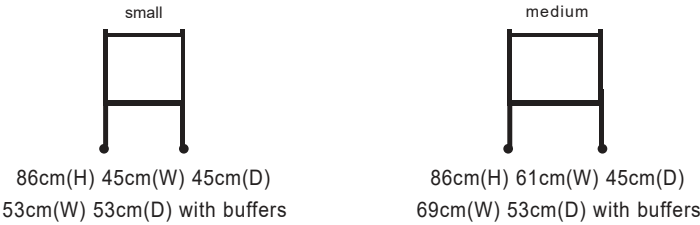

Surgical Trolleys - Stainless Steel (Mirror Polished)

Premium use of materials, and finished to the highest standard, our range of ‘mirror polished’ surgical trolleys provide hygienic mobile support for clinical environments and theatres, whilst being functional, durable and easy to clean.

small - 45cm wide
surgical trolleys (mirror polished)



Code: Sun-STFW4/RRFS2
2 removable reversible folded shelves (flange up or down)



Code: Sun-STFW4/FFS2
2 fully welded fixed shelves (flange down)



product information

premium
quality

Stainless steel shelves, trays and frame (304 grade)

Mirror polished stainless steel shelves, trays and frame

Fully welded seamless frame

BS 4068 : 1977 Compliant *

2 shelf
options

Flat, fixed or removable/reversible shelves

Removable, reversible folded shelves (Sun-STFW/RRFS2 only)

Fully welded fixed shelves (Sun-STFW/FFS2 only)

Anti-static shelf/tray protection buffers

4 Protective buffers

75mm antistatic castors

Supplied fully assembled

10Kg

Maximum individual tray load


Easy clean design and materials

* Designed to comply with the requirements of the hospital market (BS4068: 1977)

accessories * These accessories can only be fitted when ordered with the trolley

Sun-LC1A*
4 Antistatic castors
(front 2 braked)

dimension

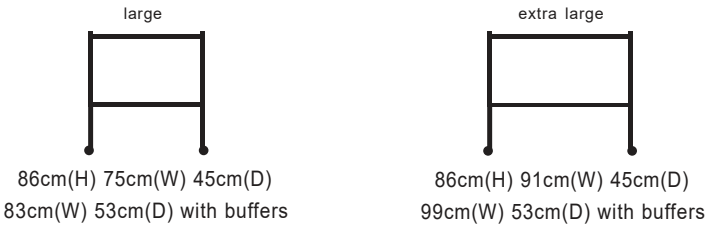

Surgical Trolleys - Stainless Steel (Mirror Polished)

Premium use of materials, and finished to the highest standard, our range of ‘mirror polished’ surgical trolleys provide hygienic mobile support for clinical environments and theatres, whilst being functional, durable and easy to clean.

large - 75cm wide
surgical trolleys (mirror polished)

extra large - 91cm wide
surgical trolleys (mirror polished)



Code: Sun-STFW7/RRFS2
2 removable reversible folded shelves (flange up or down)



Code: Sun-STFW9/RRFS2
2 removable reversible folded shelves (flange up or down)



Code: Sun-STFW7/FFS2
2 fully welded fixed shelves (flange down)



Code: Sun-STFW9/FFS2
2 fully welded fixed shelves (flange down)

product information




Sun-IV01, IV02



Simple locking height adjustment



Complies with BS3619:1976

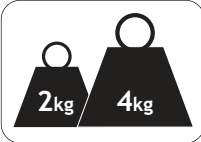

2kg hook load or 4kg total



Lightweight design approx 3kg (Sun-IV01, IV02)



Stable 5 star Plastic castor base (Sun-IV01 & IV02)



Simple self assembly required




Easy clean design and materials



British made

* Please contact us regarding availability for quantities of 10 or more.

dimensions

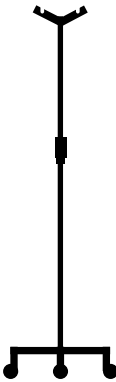

Sun-IV01
160cm - 245cm(H)
51cm(W)



Sun-IV02
160cm - 245cm(H)
51cm(W)

Mobile Drip Stands

Simple, practical, easy to adjust and manoeuvre, our drip stands are suitable for most healthcare environments. We can provide a lightweight design through to ones with heavy duty weighted bases, available in a choice of materials to suit your budget and hygiene requirements.



Code: Sun-IV01
Chrome steel drip stand,
plastic base, 2 plastic hooks



Code: Sun-IV02
Chrome steel drip stand,
plastic base, 4 chrome hooks

*stylish
modular
storage*

HTM Compliant Cabinets & Worksurfaces

Our HTM compliant cabinets, worksurfaces, and mobile units provide durable, hygienic, and flexible storage for clinical environments. Practitioner drawers, pull-out units, and project cabinets ensure safe, organised, and efficient solutions for healthcare settings.

Cabinetry



Beech Cabinet Range
Pg 87



White (High Gloss) Cabinet Range
Pg 88

Plumbing Solutions



Integrated Plumbing Systems
Pg 89



Sinks (HTM64 Compliant)
Pg 90



Taps (HTM64 Compliant)
Pg 91

Fridges



Pharmacy Fridges
Pg 92

Beech Cabinet Range

product information



Quality, 18mm MFC carcass with 2mm edging



Locking doors supplied with 2 keys



Contemporary door and drawer profile



Stylish flush fitting handles



38mm deep dark grey laminate work top



Base cabinets have one adjustable shelf, wall cabinets have two



Easy clean design and materials



FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe

- Base cabinets can be combined and supplied with a one piece work top in lengths of up to 4.1m
- Cut to length worktop can be supplied upon request

accessories



Sun-DT1 Drawer Insert for 50cm Wide Drawers



Integrated Plumbing System (IPS) (see pages 150-151)



Sinks (HTM64) (see pages 152-155)




Taps (HTM64) (see pages 156-159)




Sinks & Taps (General) (see pages 160-161)

dimensions All single door base & wall cabinets are supplied right handed (unless otherwise specified)


Sun-BU8B 30cm Base Cabinet

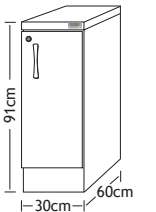


Sun-BU1B 50cm Base Cabinet

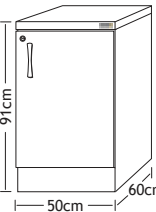


Sun-BU2B 50cm Drawerline Cabinet

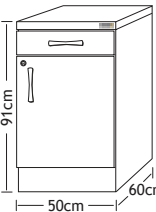


Sun-BU3B 50cm 5 Drawer Pack

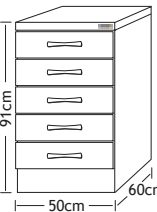


Sun-BU4B 100cm Base Cabinet

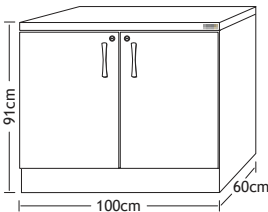

Sun-BU5B 100cm Drawerline Cabinet

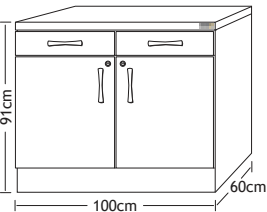


Sun-BU6B 100cm Sink Cabinet (excludes sink)

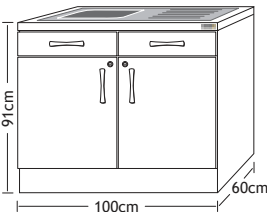


Also available in 40cm, 60cm and 80cm widths


Sun-WU1B 30cm Wall Cabinet

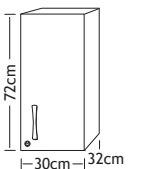


Sun-WU2B 50cm Wall Cabinet

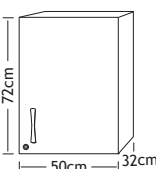

Sun-WU3B 100cm Wall Cabinet

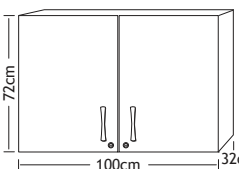


Sun-WU4B Corner Wall Cabinet

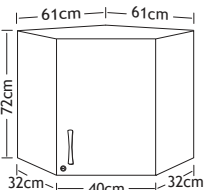

Also available in 40cm, 60cm and 80cm widths

Beech Cabinet Range

Aesthetically designed, contemporary range of modular cabinetry incorporating wall, base, corner and sink units finished in beech effect. Easy to clean, adjustable shelves and locks as standard.

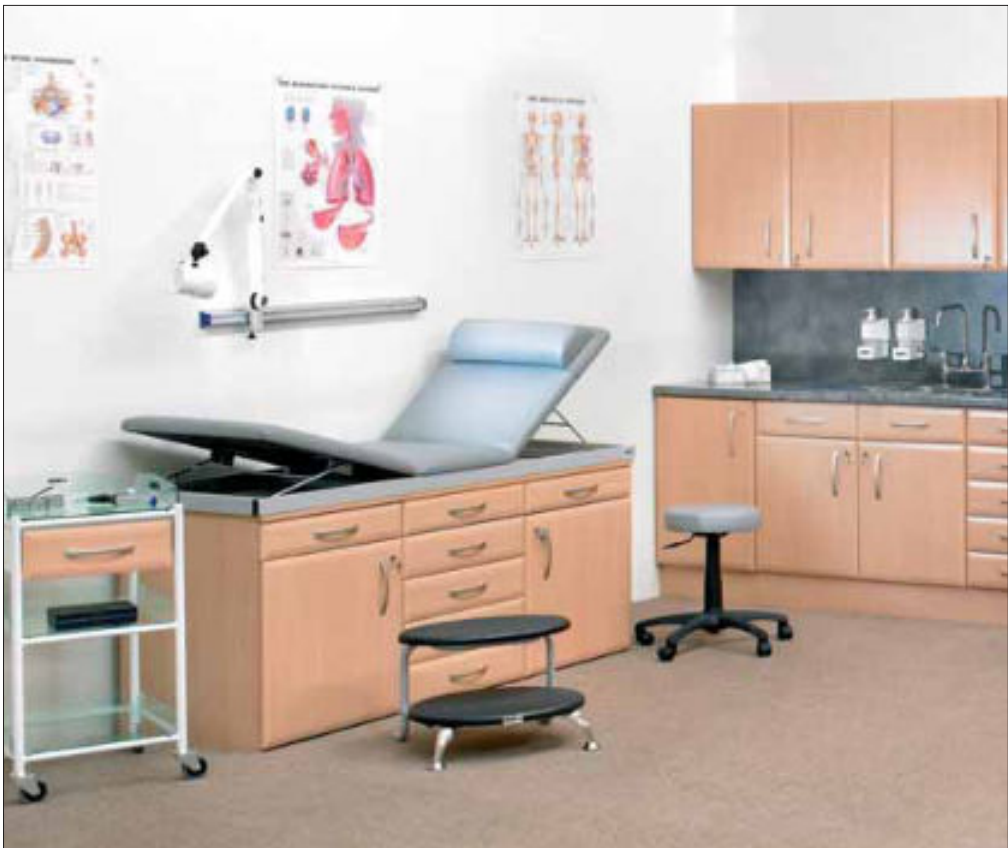

White (High Gloss) Cabinet Range


product information



Quality, 18mm MFC carcass with 2mm edging




Locking doors supplied with 2 keys

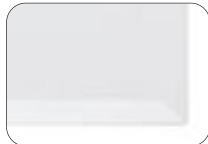

Contemporary door and drawer profile



Stylish flush fitting handles

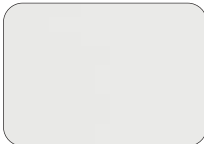

38mm deep crystal white laminate work top



Base cabinets have one adjustable shelf, wall Cabinets have two



Easy clean design and materials



FIRA Standard 6250 2005 Specification: Adhesion Performance General/Severe

- Base cabinets can be combined and supplied with a one piece work top in lengths of up to 4.1m
- Cut to length worktop can be supplied upon request

accessories



Sun-DT1 Drawer Insert for 50cm Wide Drawers



Integrated Plumbing System (IPS) (see pages 150-151)



Sinks (HTM64) (see pages 152-155)




Taps (HTM64) (see pages 156-159)




Sinks & Taps (General) (see pages 160-161)

dimensions All single door base & wall cabinets are supplied right handed (unless otherwise specified)


Sun-BU8W 30cm Base Cabinet

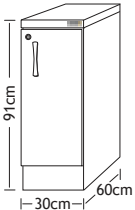


Sun-BU1W 50cm Base Cabinet

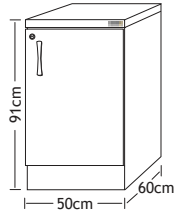


Sun-BU2W 50cm Drawerline Cabinet

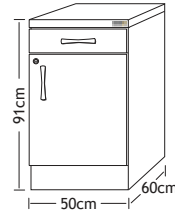


Sun-BU3W 50cm 5 Drawer Pack

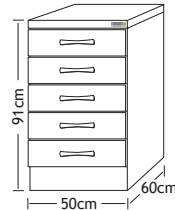


Sun-BU4W 100cm Base Cabinet

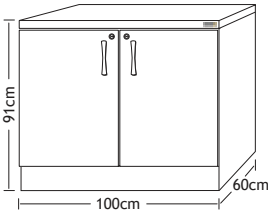

Sun-BU5W 100cm Drawerline Cabinet

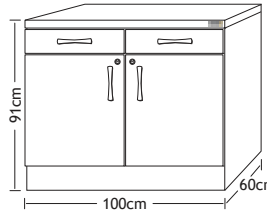


Sun-BU6W 100cm Sink Cabinet (excludes sink)

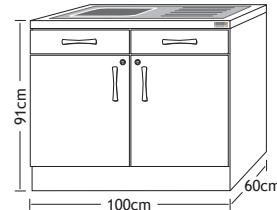


Also available in 40cm, 60cm and 80cm widths


Sun-WU1W 30cm Wall Cabinet

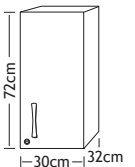


Sun-WU2W 50cm Wall Cabinet

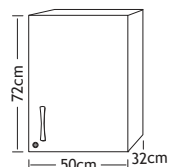

Sun-WU3W 100cm Wall Cabinet

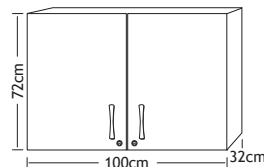

Sun-WU4W Corner Wall Cabinet

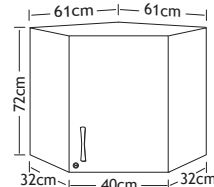

Also available in 40cm, 60cm and 80cm widths

White (High Gloss) Cabinet Range

Aesthetically designed, contemporary range of modular cabinetry incorporating wall, base, corner and sink units finished in high gloss white. Easy to clean, adjustable shelves and locks as standard.



Supplied fully assembled



product information

HTM64
compliant

For use in
clinical areas



Hygienic and durable
design



Removable
postformed panels

18mm
LFMDF

Laminate faced
MDF panels

4 Widths
3 Heights

4 standard widths
in 3 different heights



Easy clean design
and materials



FIRA
Independently
tested by FIRA

installation

Professional
installation available

FIRA
Independently
tested by FIRA

FIRA Standard 6250 2005
Specification:
Adhesion Performance
General/Severe

accessories * Required for sink installation and for HTM64 compliance



Sun-SNK14
Washbasin - medium
(see pages 154-155)

FIXING
kit

Sun-SNK14/KIT*
Panel fixing kit for
Sun-SNK14

SINK
adapter

Sun-SNK14/KIT1*
Sink outlet adapter
for Sun-SNK14



Sun-SNK15
Washbasin - large
(see pages 154-155)

FIXING
kit

Sun-SNK15/KIT*
Panel fixing kit for
Sun-SNK15

SINK
adapter

Sun-SNK15/KIT1*
Sink outlet adapter
for Sun-SNK15



Sun-DISPENSER4
Paper towel
dispenser




Sun-TAP15
Electronic thermostatic
mixer tap with
time flow sensor



Sun-TAP17
Thermostatic
sequential mixer tap

dimensions & codes

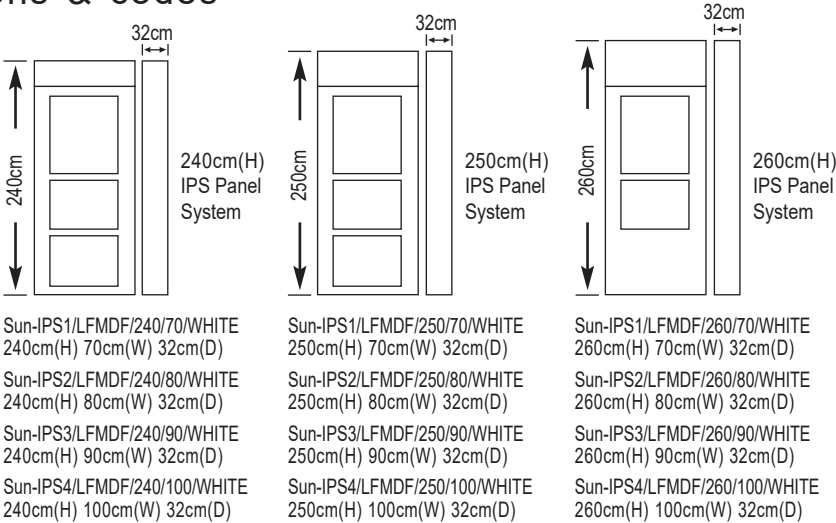

Our range of HTM64 compliant integrated plumbing systems (IPS), combine a range of healthcare sanitary ware fitted into a specially designed pre-assembled panel arrangement. They are suitable for use in both clinical and personal washing areas and provide a cost-effective, functional and easy to maintain solution.



Sinks (HTM64 Compliant)

product information - Sun-SNK24 & 25

HTM64
compliant

For use in clinical areas

No plug, chain or overflow

Available in left or right hand drainer

Premium quality stainless steel

200mm tap centres

Fixing kit supplied

Bottle Trap & Waste Kit

Bottle trap and waste kit needed (see accessories)*

Easy clean design and materials

British made

* Required for sink installation and for HTM64 compliance

accessories - (please see pages 167-172 for our full range & specification of taps)

Sun-SNK24 & 25/WK*
Waste kit


Sun-SNK24 & 25/BT*
Bottle trap

Sun-TAP13
Deck mounted sequential
Thermostatic mixer tap

* Required for sink installation and for HTM64 compliance


dimensions

Sun-SNK24 & Sun-SNK25

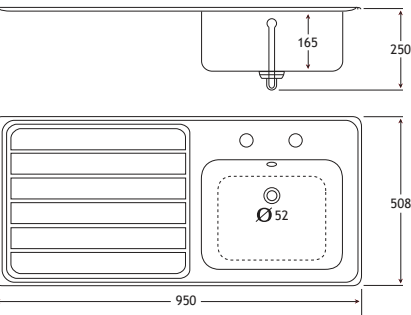

dimensions

Sun-SNK14 & Sun-SNK15

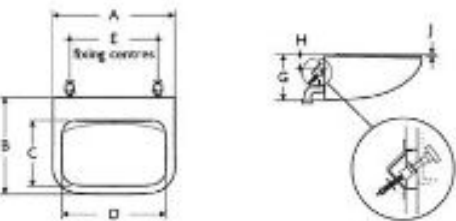

	A	B	C	D	E	G	H
Sun-SNK14 (M)	500	400	270	425	360	187	53
Sun-SNK15 (L)	600	480	342	500	450	196	55

J

product information - Sun-SNK14 & 15

HTM64
compliant

For use in clinical areas

No plug, chain or overflow

Suitable for IPS systems

Back Outlet Adapter

Back outlet adapter needed (see accessories)*

Sinks (HTM64 Compliant)

Our range of HTM64 compliant inset stainless steel sinks and bowls are suitable for clean and dirty areas as well as clinical healthcare applications.

Inset Sinks

HTM64 compliant inset stainless steel sink, available with left or right hand drainer.



Code: Sun-SNK24 - Inset stainless steel sink with right hand drainer



Code: Sun-SNK25 - Inset stainless steel sink with left hand drainer

Washbasins

Our range of HTM64 compliant washbasins are suitable for use in clinical areas and for duct panel on integrated plumbing systems (IPS). Made from vitreous china to BS3402.



Code: Sun-SNK14
Washbasin (rear draining) - medium



Code: Sun-SNK15
Washbasin (rear draining) - large

Taps (HTM64 Compliant)

product information - Sun-TAP13



Suitable for clinical use



200mm tap centres



Integrated thermostat regulates temperature



180° Span single lever operation



Built in thermostat cleansing feature (requires bridging hose)



Durable attractive chrome finish



Easy clean design and materials



1 year warranty



British made

product information - Sun-TAP15



Suitable for clinical use



Pre-determined water running length time (set up on installation)



Suitable for scrub up applications



200mm tap centres



Integral thermostat



Built in thermostat cleansing feature (requires bridging hose)



Durable attractive chrome finish



Easy clean design and materials



British made

product information - Sun-TAP17



Suitable for clinical use



For use in general areas



200mm tap centres



Integral thermostat



180° Span single lever operation



Built in thermostat cleansing feature (requires bridging hose)



Durable attractive chrome finish




Easy clean design and materials




British made

Taps (HTM64 Compliant)


Our range of HTM64 compliant taps are suitable for use in clinical areas and in duct panels on integrated plumbing systems (IPS).

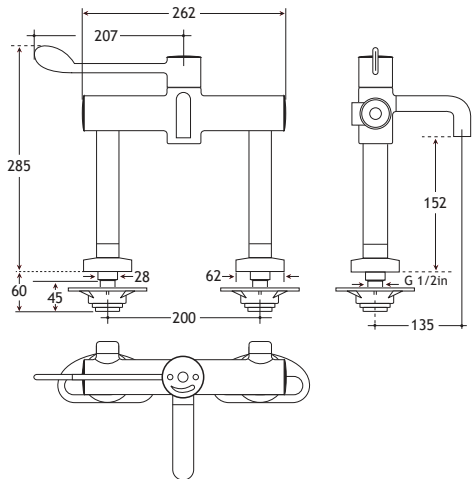

Code: Sun-TAP13
Deck mounted sequential thermostatic mixer tap

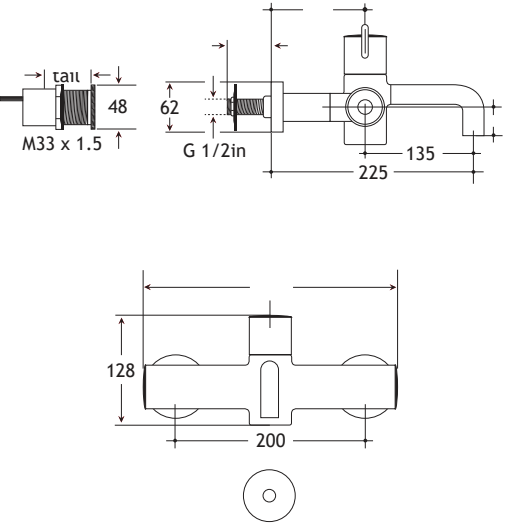

Code: Sun-TAP15
Electronic thermostatic mixer tap with time flow sensor

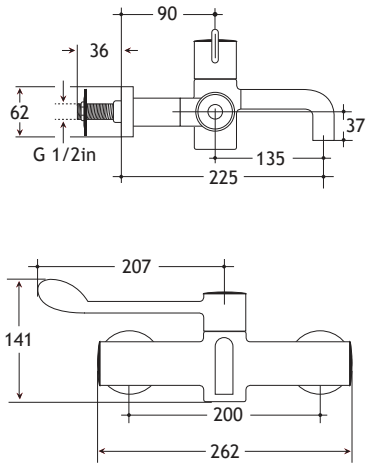

Code: Sun-TAP17
Thermostatic sequential mixer tap

product information



5 sizes

5 sizes/capacities available



Wall mounting, worktop mounted, under worktop or freestanding



Lockable doors (2 keys supplied)



Hi-Low temperature & door open alarm



mains failure alarm

Via battery backup



Off cycle automatic defrost



ECO friendly refrigerant



DH Green Book & RPSGB/GPhC compliant



Fan assisted air cooling for rapid temperature recovery



Internal light fitted as standard on all models



Temperature range 2°C to 8°C



Glass doors for visible stock control



Real time temperature display



Min/max temperature recording



Accurate temperature control



Solid doors provide extra security



Temperature data storage



Adjustable shelves



White painted corrosion proof aluminium interior



Low energy consumption



Easy clean design and materials




British made



Supplied fully assembled

* 5 year warranty upon registration by the end user (UK only) 1 year warranty for export orders

dimensions

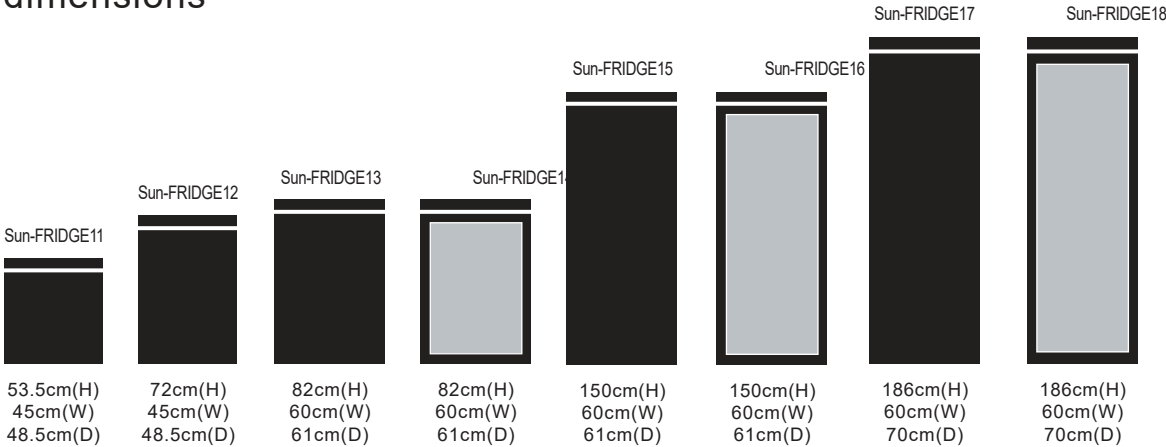

Pharmacy Fridges

Designed specifically to meet the demands of health care professionals and for the storage of temperature sensitive vaccines and pharmaceuticals in accordance with manufacturers licensed temperature ranges and guidance. Low energy use, fan assisted cooling, flexible shelving and a 5 year warranty all come as standard.



Code: Sun-FRIDGE11
36 Litre wall/worktop mounted fridge - 2 shelves



Code: Sun-FRIDGE12
66 Litre wall/worktop mounted fridge - 3 shelves



Code: Sun-FRIDGE13
150 Litre under worktop fridge - 3 shelves



Code: Sun-FRIDGE14
150 Litre under worktop, glass door fridge - 3 shelves



Code: Sun-FRIDGE15
300 Litre free standing fridge - 6 shelves
Code: Sun-FRIDGE17 (not shown)
430 Litre free standing fridge - 6 shelves



Code: Sun-FRIDGE16
300 Litre free standing, glass door fridge - 6 shelves
Code: Sun-FRIDGE18 (not shown)
430 Litre free standing, glass door fridge - 6 shelves

upholstery colours

Anti-bacterial Vinyl Upholstery

Quality

- BS 5790
- ISO 7617 Anex D

Flame Retardancy

- BS 7176:2007 + A1:2011 Medium Hazard
- BS EN 1021-1:2006 (smouldering cigarette)
- BS EN 1021-2:2006 (match flame equivalent)
- BS 5852:2006 Clause 11, Ignition Source 5
- FAR/JAR Part 25:App. F part 1 (Vertical Test) / EASA CS.25.853
- IMO Resolution MSC 307(88) (2010 FTP Code : Annex 1:Part 8

Martindale rubs

- up to 1.000.000

Cleaning & Maintenance Instructions

Our vinyl-coated fabrics are designed for durability and easy maintenance. To keep them looking fresh and in top condition, follow these simple cleaning steps:

- Regular Cleaning – The best way to maintain the fabric's appearance is to clean it regularly.
- Use a Mild Soapy Water Solution – Avoid harsh chemicals; a gentle soap and water mix is sufficient.
- Soft Brush or Microfibre Cloth – For best results, use a soft brush or a microfibre cloth to clean the surface without causing damage.
- Quick Wipe-Downs – A damp cloth with soapy water is all it takes to remove dust and minor stains.
- Seek Assistance if Needed – If unsure about specific stains or cleaning methods, professional advice may help.

Perfect for busy environments and everyday use, these materials are designed to be easy to clean and maintain, ensuring long-lasting durability.

This antimicrobial vinyl provides protection against C. difficile, MRSA, E. coli, and Covid-19.

To support infection control efforts, an antimicrobial agent is incorporated into the vinyl during the production process. Unlike surface treatments, this integration ensures the agent remains effective even after regular cleaning. In addition to its antimicrobial properties, the material resists the growth of mould and mildew, helping to prevent musty or stale odours.

Everyday Range



APPLE



BLACK



BLUEBELL



BLUEBERRY




FEATHER



LAPIS



SNOW-WHITE

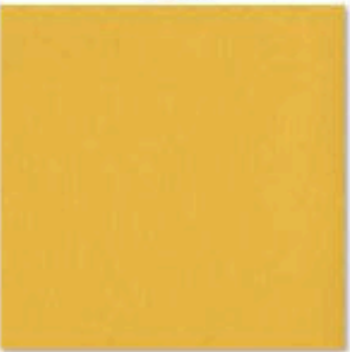

SUNFLOWER

Kaleidoscope Range

Lifestyle Range



BLACKCURRANT



BLOSSOM



CHOCOLATE



CORNFLOWER



DAMSON



FIRE-CRACKER



FLAKE



LEAF-GREEN



LIME



MANGO



MUSH-ROOM



NAVY



ORCHID




PEACH




ROSE-IP




TRUFFLE

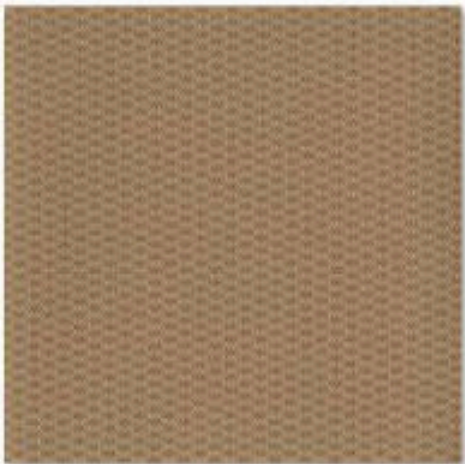


ALDWARK

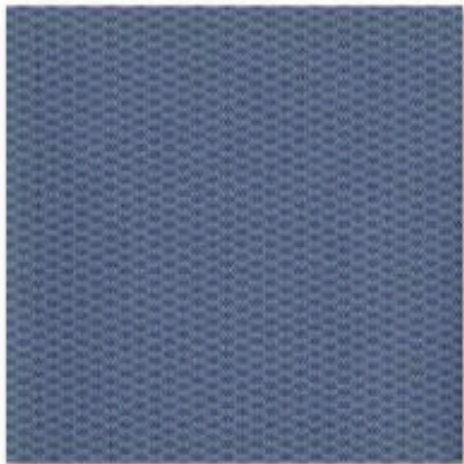


DEANGATE

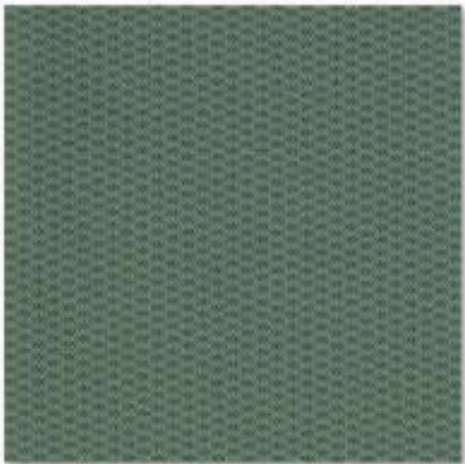

FOSS-BANK

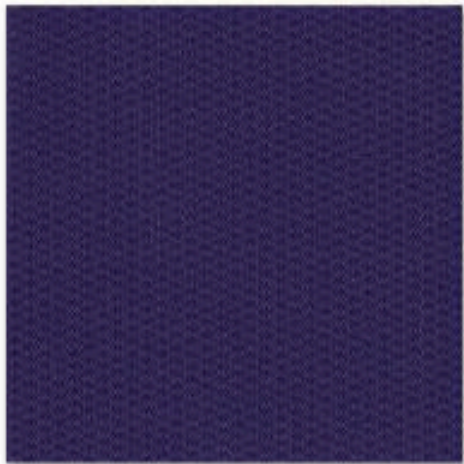

GRAPE

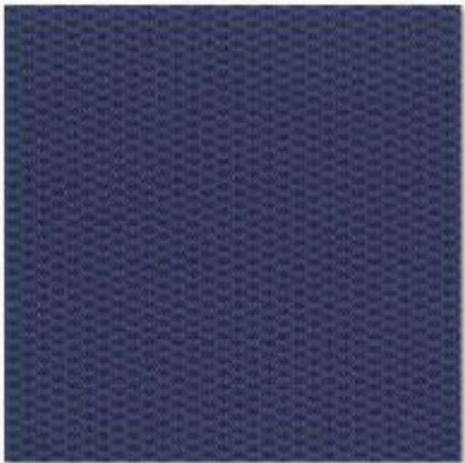

MINSTER



PRIORY



TOFT-GREEN





UK's leading **manufacturer**
of healthcare furniture



Based in Bradford, West Yorkshire we are proud of our UK manufacturing facility which undertakes its design, manufacturing, assembly and distribution all from one central site. Having the ability to design new innovative products, manufacture and supply them from start to finish allows us to provide bespoke solutions and outstanding levels of service.



our promise

To create smart, contemporary, functional furniture, we start with intuitive design and solid engineering. Then, we use the best materials and subject our products to rigorous testing. That's what makes a Sunflower product and why our warranty and service is a **real promise**.




We **manufacture**
on site and in the UK



Our advanced manufacturing facilities enable us to produce high-quality items with precision and efficiency, whether for single pieces or large orders in a variety of colours and finishes.

With a combination of cutting-edge technology and skilled craftsmanship, we ensure consistency, durability, and a flawless finish. Our streamlined processes also allow for quick turnaround times while maintaining the flexibility to meet specific customer delivery requirements.

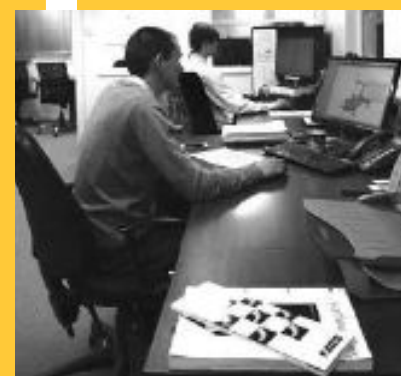

In-house design with
3D CAD technology



Our specialised product designers have many years of experience in the healthcare market. Designing from 'concept to reality' using the latest solid modelling CAD software our designers can create anything from a new single product, something bespoke or a whole new room or environment to your exact requirements.



our products & service

We design medical furniture to help healthcare become more efficient, combined with high service levels, extensive warranties and consideration for the environment.

We believe passionately in innovative design and continuous development, utilise modern materials and technology to produce stylish and hygienic new products and solutions.

Ease of cleaning and hygienic design is carefully incorporated into everything we do.

We are committed to reducing our impact on the environment and utilise sustainable materials and sources wherever possible, achieving ISO and FSC standards. We are committed to recognised industry standards, giving you confidence in our company, our products and our people.



Suppliers to the NHS

NHS Supply Chain Contract:
Medical Furniture & Sack Holders (FAG15887–OJEU ref: 2014/S 198-349358)

NHS Supply Chain Contract:
Ward & Residential Furniture (FAG000016544–OJEU ref: 2017/S 056-104552)

NHS Shared Business Services Contract:
Medical Furniture (SBS/16/SB/TBP/9050/12 – OJEU 2016/S 174-312172)

For more details visit my.supplychain.nhs.uk/catalogue



Warranty

We provide a minimum 1 year warranty on all products, with the majority having between 2 to 5 years and some components for life, giving you complete peace of mind.



Express Delivery

Need an item quickly? – We offer a range of products through our express system with delivery normally within 3-5 working days, just look for the Express symbol.



Delivery & Installation

We can deliver and install your new furniture anywhere in the UK using our own vehicles and trained installers, please enquire for further details.



Bespoke Projects & Product Personalisation

We offer a bespoke design, manufacturing, fitting, and project management service for new or existing healthcare premises. Many of our products can be customised with bespoke graphics, such as hospital logos or ward names, to meet your specific needs.



Export

We can export our products worldwide and have established distributors around the globe. Contact our export department at export@sunflowermedical.co.uk for further details, or if you would like to become part of our global distributor network.

Contact Us & How To Order

T: +44(0)1274684004

E: sales@sunflowermedical.co.uk

Contact us with your requirements or for details of our approved national distributors.

Alternatively buy online at: www.sunflowermedical.co.uk



Recommended Material Cleaning Guidelines:



We make every effort during our design process to ensure our products are easy to clean and utilise modern materials and processes where possible to assist in the prevention of cross infection. We recommend that all surfaces should be thoroughly cleaned and disinfected regularly, and have compiled the following guidelines to assist you with your cleaning process:

- Vinyl & Intervene Upholstery:**

 - Alcohol based wipes are acceptable as a cleaning method such as Clinell Universal Wipes. Vigorous wiping can be used on the vinyl upholstery but the use of abrasive materials is not recommended as it will lead to scratching and damage to the surface.
 - Hypochlorite based cleaning can be undertaken up to 10000ppm. Particular care must be taken to avoid pooling of cleaning substance as concentrations increase with evaporation and cause damage.
 - Hypochlorite cleaning can, over time, cause discolouration and eventually hardening and cracking. In real world use in hospitals Sunflower has not seen this to be a problem but it is a potential issue regardless of the manufacturer of the furniture.
 - All wetness must be dried after cleaning with care been taken to dry into joints.

KYDEX and Plastic Components:

 - Alcohol based wipes are acceptable as a cleaning method such as Clinell Universal Wipes. Vigorous wiping can be used on the Kydex and Plastic but the use of abrasive materials is not recommended as it will lead to scratching and damage to the surface.
 - Hypochlorite based cleaning can be undertaken up to 10000ppm. Particular care must be taken to avoid pooling of cleaning substance as concentrations increase with evaporation and cause damage.
 - All wetness must be dried after cleaning with care been taken to dry into joints.

Powder Coating (Paint):

 - Powder coated surfaces(paint) should be cleaned using a soft damp non-abrasive cloth & wiped dry. Increased or continual exposure to wet/damp conditions or environments of high humidity may result in rusting/or cause the coating to peel, and therefore is not recommended. Some dulling of the finish may occur over time if cleaning products are consistently used.
 - Hypochlorite based cleaning can be undertaken up to 10000ppm. Particular care must be taken to avoid pooling of cleaning substance as concentrations increase with evaporation and cause damage.
 - Bleach, bleach based sprays, solvents, abrasives, synthetic detergents, wax polishes or aerosols are not recommended.

Melamine, Laminates & Vinyl Wrapping:

 - Melamine & Laminate faced surfaces(used on MFC & MDF) should be cleaned using a soft damp non-abrasive cloth and wiped dry.
 - Vinyl wrapped surfaces (e.g. Doors & Drawers) should be cleaned using a soft damp non-abrasive cloth and wiped dry.
 - Hypochlorite based cleaning can be undertaken up to 10000ppm. Particular care must be taken to avoid pooling of cleaning substance As concentrations increase with evaporation and cause damage.
 - Bleach, bleach based sprays, solvents, abrasives, synthetic detergents, wax polishes or aerosols are not recommended.

Stainless Steel:

 - Stainless steel surfaces should be cleaned using a soft damp non-abrasive cloth and wiped dry. Discolouration may occur if allowed to dry naturally. Stainless steel can still rust and certain chemicals can still damage the material. Iodine is one particular chemical which is known to corrode/stain and should be avoided for cleaning. Increased or continual exposure to wet/damp conditions or environments of high humidity may result in some rusting and therefore is not recommended.
 - Hypochlorite based cleaning can be undertaken up to 10000ppm. Particular care must be taken to avoid pooling of cleaning substance as concentrations increase with evaporation and cause damage.
- Aluminium:**

 - Aluminium surfaces should be cleaned using a soft damp non-abrasive cloth and wiped dry. Discolouration may occur if allowed to dry naturally. Increased or continual exposure to wet/damp Conditions or environments of high humidity may result in corrosion and therefore is not recommended.
 - Bleach, bleach based sprays, solvents, abrasives, synthetic detergents, wax polishes or aerosols are not recommended.

Plastics, PVC & Rubber Based Materials:

 - Plastics, PVC (e.g. screen panels) & rubber should be cleaned using a soft damp non-abrasive cloth and wiped dry.
 - Bleach, bleach based sprays, solvents, abrasives, synthetic detergents, wax polishes or aerosols are not recommended.

Electrical Components, Wiring & Actuators:

 - Any electrical components should not be exposed to water, liquids or areas of high humidity. Plastic casings, wires and handsets can be cleaned using a soft damp non-abrasive cloth and wiped dry.
 - Bleach, bleach based sprays, solvents, abrasives, synthetic detergents, wax polishes or aerosols are not recommended.

■ Always switch off the mains supply before cleaning.

Metal or MFC Components:

 - Alcohol based wipes are acceptable as a cleaning method such as Clinell Universal Wipes. Vigorous wiping can be used on the body but the use of abrasive materials is not recommended as it will lead to scratching and damage to the surface.
 - Hypochlorite based cleaning can be undertaken up to 10000ppm. Particular care must be taken to avoid pooling of cleaning substance as concentrations increase with evaporation and cause damage.
 - All wetness must be dried after cleaning with care been taken to dry into joints. Failure to dry full can lead to paint work damage and rusting over time.
 - Particular care must be taken no ensure that no liquids or cleaning materials enter the mechanism of locks. This will cause premature failure.

Castors:

 - Castors due to their nature pick up “fluff and bits”. Initially before wet cleaning all “fluff and bits” should be removed by hand. Once these are removed then further cleaning can be undertaken as below.
 - Alcohol based wipes are acceptable as a cleaning method such as Clinell Universal Wipes. Vigorous wiping can be used on the castor but the use of abrasive materials is not recommended as it will lead to scratching and damage to the surface.
 - Hypochlorite based cleaning can be undertaken up to and 10000ppm. Particular care must be taken to avoid pooling of cleaning substance as concentrations increase with evaporation and cause damage.
 - All wetness must be dried after cleaning with care been taken to dry into joints.

Shelves and Trays:

 - The shelves and trays should be cleaned individually. They should be fully wiped on all surfaces and repositioned in the trolley as cleaning is undertaken.
 - Alcohol based wipes are acceptable as a cleaning method such as Clinell Universal Wipes. Vigorous wiping can be used on the shelves and trays but the use of abrasive materials is not recommended as it will lead to scratching and damage to the surface.
 - Hypochlorite based cleaning can be undertaken less up to and including 10000ppm. Particular care must be taken to avoid pooling of cleaning substance as concentrations increase with evaporation and cause damage.

We design medical furniture to help healthcare become more efficient, combined with high service levels, extensive warranties and consideration for the environment.

We believe passionately in innovative design and continuous development, utilise modern materials and technology to produce stylish and hygienic new products and solutions.

Ease of cleaning and hygienic design is carefully incorporated into everything we do.



Infection Control

The requirement for furniture that provides functional, ergonomic hygienic storage which assists in minimising infection control issues has always been a significant factor in our design and manufacturing processes and one of continuous improvement. Smart design such as sloping top wall units help facilitate cleaning and the build-up of dust as well as preventing the storage of items which should ideally be stored elsewhere. Materials which are durable and hard wearing, moisture resistant and anti-bacterial can all be utilised as part of our furniture and storage solutions.



Environment

We are committed to reducing our impact on the environment and utilise sustainable materials and sources wherever possible. Our cabinet and worktop manufacturing processes are optimised to be as efficient as possible, reducing energy costs and minimising waste. Reducing our carbon footprint through an on-going reduction programme is part of our company process and our ISO 14001 Environmental Quality Standard Certification and FSC Accreditation (the Mark of Responsible Forestry) shows our commitment to reducing our environmental impact even further.



FIRA Standard 6250
2005 Specification:
Adhesion Performance
General/Severe



BS EN 16139: 2013
Furniture Strength,
Durability & Safety



Anti-bacterial Vinyl Strength Test
BS EN ISO 13934-1: 2013
BS EN ISO 13973-3: 2000
BS EN ISO 12947-2: 1998



Fire (Flammability) Testing
Flammability BS7176 Medium Hazard
BS5852 : 2006 Clause 11
(composites)

